



INTOXILYZER® 8000 INSTALLATION AND REPAIR CHECKOUT
 OFFICE OF ATTORNEY GENERAL
 CRIME LABORATORY DIVISION
 Toxicology Section/Breath Alcohol Program
 SFN 59281 (03-2014)

Serial Number <i>80-006509</i>	Instrument Location <i>Lanore</i>
Reason for Install/Repair <input checked="" type="checkbox"/> Install After Receiving From Crime Laboratory <input type="checkbox"/> Install After Location Change <input type="checkbox"/> Other (Specify) _____	

Check When Done:

- 1. Surge Protector Installed/Properly Grounded.
- 2. Telephone Line Connected to Intoxilyzer® 8000.
- 3. Breath Tube Heated.
- 4. Enter Preliminary Data (i.e. Date, Time, DST (Y), and Location; Level 2, Function E).
- 5. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 6. Run Tests:
 - A. Print Test (Level 1, Function P).
 - B. ACA Test (Level 1, Function C).
 - C. RFI Test (CMS Mode or Level 1, Function B or C; Key Radio During Test).
- 7. Repair and/or Maintenance Performed (if Any): *N/A*
- 8. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 9. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 10. Send the Following to the Crime Laboratory: Completed Intoxilyzer®8000 Installation and Repair Checkout (SFN59281, Form 104-G), Print Test, ACA Test, and RFI Test.

Field Inspector Signature <i>[Signature]</i>	Date <i>5-20-16</i> ✓
Reviewed By (Crime Laboratory Use Only) <i>Deb Shanover</i>	Date <i>25 May 16</i>

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006509
Location = LAMO 8164.13.00 06/09
05/20/2016 14:21

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNopqrstuvwxyz!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNopqrstuvwxyz!@#\$\$%^&*()_+?

Current Instrument Setup


Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 01316080A2
Standard Cyl #? 32
Standard Expiration? 02/05/2018
Oper No? 039706

Flow Cal. Date: 08/20/2015
Slope 695
Intercept -685441

IR Calibration Date: 04/01/2016
 3um 9um

0th Coef(*100): -42983 -22063
1st Coef(*100): 281397 138045
2nd Coef(*100): 2616 1525
H2O adj(mg/l*10k): 643 374

***** Printer Test End *****

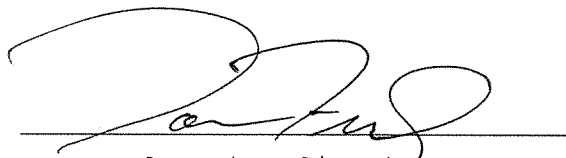

Operator Signature
THOMAS HERZIG

Remarks:


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006509
Location = LAMO 8164.13.00 06/09
05/20/2016 14:13

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:15
02 Std. Gas	0.075	14:15
03 Room Air	0.000	14:16
04 Std. Gas	0.076	14:16
05 Room Air	0.000	14:17
06 Std. Gas	0.076	14:17
07 Room Air	0.000	14:18

Lot No = 01316080A2
Cyl No = 32
Exp Date = 02/05/2018
County = 23 Oper No. = 039706

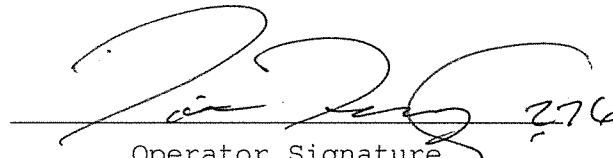
Operator Signature
THOMAS HERZIG

Remarks:

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006509
Location = LAMO 8164.13.00 06/09
05/20/2016 14:19

DRY CAL CHECK

Test	AC	Time
01 Room Air	RFI*	14:19
02 Room Air	0.000	14:20

*Invalid Test
Inhibited - RFI

Lot No = 01316080A2
Cyl No = 32
Exp Date = 02/05/2018
County = 23 Oper No. = 039706



Operator Signature
THOMAS HERZIG

Remarks: