



INTOXILYZER® 8000 INSTALLATION AND REPAIR CHECKOUT
NORTH DAKOTA OFFICE OF ATTORNEY GENERAL
CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM
SFN 59281 (06/2018)

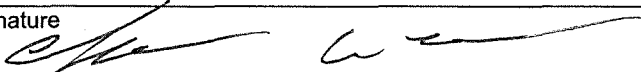
Serial Number 80-006678	Instrument Location 5FPD
Reason for Install/Repair ☐ Install After Receiving From Crime Laboratory ☐ Install After Location Change ☒ Other (Specify) Taken to UND for download and then back to 5FPD	

Check When Done:

- ☒ 1. Surge Protector Installed/Property Grounded.
- ☒ 2. Telephone Line Connected to Intoxilyzer® 8000.
- ☒ 3. Breath Tube Heated.
- ☒ 4. Enter Preliminary Data (i.e. Date, Time, DST (Y), and Location; Level 2, Function E).
- ☒ 5. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- ☒ 6. Run Tests.
 - ☐ A. Print Test (Level 1, Function P).
 - ☒ B. ACA Test (Level 1, Function C).
 - ☒ C. Radio Frequency Interference (RFI) Test (CMS Mode or Level 1, Function B or C; Key Radio During Test).
- ☒ 7. Repair and/or Maintenance Performed (if any): **Download @ UND**
- ☒ 8. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- ☒ 9. File Previous Intoxilyzer® Record (SFN504096, Form 120-G) at the Intoxilyzer® Location at the Agency.
- ☒ 10. Send the Following to the Crime Laboratory: Completed Intoxilyzer® 8000 Installation and Repair Checkout (SFN59281, Form 104-G), Print Test, ACA Test, and RFI Test.

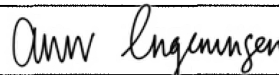
AEI

10/8/26

Field Inspector Signature 	Date 10-8-26
--	------------------------

Crime Laboratory Use Only

This installation has been reviewed and the instrument is approved to be used for the analysis of breath to determine alcohol concentration from the date the Field Inspector performed the installation. This record on file at the Office of Attorney General, Crime Laboratory Division, in the County of Burleigh, North Dakota, is certified to be a true and correct copy of the documents received.

Reviewed/Certified By Anna Ingemansen 	Certified Date 14Oct2026
---	------------------------------------

C. was checked by Anna Ingemansen. This printout is attached to the Installation Paperwork.
This does not affect the install of this instrument. AEI 10/08/2026

CMI, Inc. Intoxilyzer	Alcohol Analyzer
North Dakota Model 8000	SN 80-006678
Location = GFPD	8164.14.00 09/16
10/08/2026	14:34

```
abcdefghijklmnopqrstuvwxyz1234567890-=|
ABCDEFGHIJKLMNOPQRSTUVWXYZ!@#$%^&*()_+?
```

0th Coef(*100):	-15445	-21650
1st Coef(*100):	251854	131766
2nd Coef(*100):	3094	1378
H2O adj(mg/l*10k):	549	513

Ch L #642

AEI

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006678
Location = GFPD 8164.14.00 09/16
10/08/2026 14:35

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:36
02 Std. Gas	0.081	14:36
03 Room Air	0.000	14:37
04 Std. Gas	0.081	14:37
05 Room Air	0.000	14:37
06 Std. Gas	0.081	14:38
07 Room Air	0.000	14:38

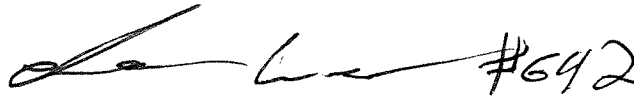
Lot No = 264437

Cyl No = 123

Exp Date = 01/29/2028

County = 18

Oper No. = 133510



Operator Signature
ADAM COLE

Remarks:

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006678
Location = GFPD 8164.14.00 09/16
10/08/2026 14:40

Test	AC	Time
01 Room Air	RFI*	14:41
02 Room Air	0.000	14:42

*Invalid Test
Inhibited - RFI

Sub Name = TEST, TEST NA
Sub DOB = 01/01/2000
Sub Sex = Male Weight = 100
Test = DUI Cit = 1
Dr. Lic. = ND/NA
Lot No = 264437
Cyl No = 123
Expiration Date = 01/29/2028
County = 18 Oper No. = 133510



Operator Signature
ADAM COLE

Remarks:

Form 106-I8000