



ETHANOL GAS STANDARD CYLINDER REPORT
 OFFICE OF ATTORNEY GENERAL
 CRIME LABORATORY DIVISION
 Toxicology Section/Breath Alcohol Program
 SFN 59282 (10-2013)

Chemical Test Operator Name (Print) EBERLE, FRANK THEODORE		
Location BOWMAN COUNTY 104 1ST. NW SUITE #9 BOWMAN, ND	Intoxilyzer® Serial Number 80-004188	
Gas Lot Number 14214080A6	Gas Cylinder Number 067	Gas Expiration Date 07/05/2016

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C).
 Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
 - A. 0.083 AC
 - B. 0.082 AC
 - C. 0.082 AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 5. Send the Following to the Crime Laboratory:
 - A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).
 - B. ACA Test Record.

Chemical Test Operator Signature <i>Frank Eberle</i>	Date 1-12-2015 ✓
Reviewed By (Crime Laboratory Use Only) <i>Deb Kasher</i>	Date 21 Jan 2015



ETHANOL GAS STANDARD CYLINDER REPORT
 OFFICE OF ATTORNEY GENERAL
 CRIME LABORATORY DIVISION
 Toxicology Section/Breath Alcohol Program
 SFN 59282 (10-2013)

Chemical Test Operator Name (Print) EBERLE, FRANK THEODORE		
Location BOWMAN COUNTY 104 1ST. NW SUITE #9 BOWMAN, ND		Intoxilyzer® Serial Number 80-004188
Gas Lot Number 14214080A6	Gas Cylinder Number 067	Gas Expiration Date 07/05/2016

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C).
 Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
 - A. 0.083AC
 - B. 0.082AC
 - C. 0.082AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 5. Send the Following to the Crime Laboratory:
 - A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).
 - B. ACA Test Record.

Chemical Test Operator Signature <i>Frank Eberle</i>	Date 1-12-2016 <i>wrong year DK</i>
Reviewed By (Crime Laboratory Use Only)	Date

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = BOWM 8164.13.00 06/09
01/12/2015 10:30

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:31
02 Std. Gas	0.083	10:32
03 Room Air	0.000	10:32
04 Std. Gas	0.082	10:33
05 Room Air	0.000	10:33
06 Std. Gas	0.082	10:33
07 Room Air	0.000	10:34


Lot No = 14214080A6

Cyl No = 067

Exp Date = 07/05/2016

County = 06

Oper No. = 130015

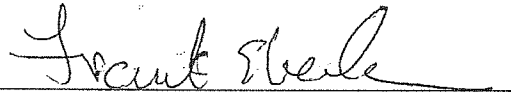

Operator Signature
FRANK EBERLE

Remarks:

Form 106-I8000