

**INTOXILYZER® 8000 INITIAL INSPECTION**

Intoxilyzer® 8000 Serial Number: 80-006495

**A. GENERAL Setup and Checks:**

1.  Power On
  - 5 min. countdown during initial warm-up
  - 1 min. countdown to begin test after standby mode
2.  Display Intensity (Ctrl + ↑ or ↓)
3. Password Security:
  - Level 1 (Operator)     Level 2 (Field Inspector)     Level 3 (Lab Personnel)
4.  Printer Setup (Level 3,S,P)
5.  Forms Load (CMI, Inc., CoBrA Software)

**B. LEVEL 1 Setup and Checks:**

1.  Internal Printer Test and Paper Feed (Level 1,P). Attach record.
2.  External Printer Test (Level 1,P). Attach record.
3. Barcode Reader (Level 1,B):
  - Operator ID     ND License     R key: Refusal. Attach record.

**C. LEVEL 2 Setup and Checks:**

1.  Continuous Air Blank (Level 2,A)
2.  Diagnostic (Level 2,D). Attach record.
3.  Prelim Data Entry (Level 2,E)
4. Ser Number/Version Info (Level 2,V):
  - Inst Ser#: 80-006495    DSP Ver: 1362.14
  - Core Ver: 8164.13    CPLD Ver: 7

**D. LEVEL 3 Setup and Checks:**

1. DVM Monitor (Level 3,D,D): 3: 12718 s: 1 C: 47.00  
9: 13502 s: 1 B: 45 0  
(RFI indicator with keyed radio): 1
2. Function (F) Keys :  F5 – Tone                       F9 – Fan  
 F6 – Green LED                       F10 – Pump  
 F7 – Red LED                           F11 – Solenoid
3. Atmospheric Monitor (Level 3,D,A): Atmo Sensor, P: 965 CF: 1.04
4. Flow Monitor (Level 3,D,F): Flow Sensor, R: 230 F: 0.000 L/s

5. Temperature Monitor (Level 3,D,T): C: 47.0 A: 49.5 I: 0.0  
B: 45.0 V: 0.0
6. Tank Monitor (Level 3,D,G): P (display): 1000 psi P (gauge): ~1000 psi
7.  General Setup (Level 3,S,G)
8.  Configure Start Test (Level 3,S,T)
9.  Network (Level 3,S,N)

**E. TESTS (attach all test records):**

**Configure wet bath simulator standards before running tests (Level 1,S):**

1. Zero AC (Level 1,C):  
Sim ser #: MP3067 Lot #: MILLI Q H<sub>2</sub>O Exp date: NA  
Known value: 0.000 AC  
Results: 0.000 AC 0.000 AC 0.000 AC
2. Low AC (Level 1,C):  
Sim ser #: MP3066 Lot #: 15160 Exp date: 7/14/17  
Known value: 0.020 AC  
Results: 0.019 AC 0.019 AC 0.019 AC
3. High AC (Level 1,C):  
Sim ser #: MP3068 Lot #: 14102 Exp date: 4/23/16  
Known value: 0.400 AC  
Results: 0.395 AC 0.394 AC 0.396 AC

**Configure dry gas standard before running tests (Level 1,S):**

4. Ethanol + 0.05% Acetone (Level 1,B):  
Sim ser #: DR5135 Lot #: ICS 4 Exp date: 9/9/2017  
Known value: 0.10 AC + 0.05% ACETONE  
 Display reads "Interferent Detect"
5. Acetone (Level 1,B):  
Sim ser #: DR5143 Lot #: ACB 2 Exp date: 9/29/17  
Known value: 0.05% ACETONE  
 Display reads "Interferent Detect"
6. Methanol (Level 1,B):  
Sim ser #: DR3377 Lot #: MET 2 Exp date: 9/29/17  
Known value: 0.10% METHANOL  
 Display reads "Interferent Detect"
7. Isopropanol (Level 1,B):  
Sim ser #: DR3453 Lot #: IPA 2 Exp date: 9/29/17  
Known value: 0.10% ISOPROPRANOL  
 Display reads "Interferent Detect"

- 8.  Ambient Fail (Level 1,B): Display and test record read correct.
- 9.  Invalid Sample (Level 1,B): Display and test record read correct.
- 10.  Sequence Aborted (Level 1,B): Display and test record read correct.
- 11.  ABA (Level 1,B): Test record reads correct.
- 12.  Calibration Check Out of Tolerance (CMS mode): Display and test record read correct.
- 13.  RFI Detected (B,C, or CMS mode): Display and test record read correct.
- 14.  Tank Pressure Below Minimum (CMS mode): Display reads correct.
- 15.  Deficient Samples (CMS mode): Display and test record read correct.
- 16.  Improper Sample (CMS mode): Display and test record read correct.
- 17.  Difference Too Great (CMS mode): Display and test record read correct.
- 18.  CMS (CMS mode): Test record reads correct.
- 19. Stability Test (Level 3,T,S):

Ave: 0.0814 Std Dev: 0.0005 Rel Std Dev: 0.62%

Maintenance/Repairs/Notes: \_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

Instrument is acceptable to be used in the field.

Charles E. Eden  
Field Inspector's Signature

21 JAN 16  
Date

Deb Shanaver  
Reviewed by

21 Jan 16  
Date

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer  
North Dakota Model 8000 SN 80-006495  
Location = TOXL 8164.13.00 06/09  
01/21/2016 13:34

\*\*\*\*\* Printer Test \*\*\*\*\*

12 char/line  
abcdefghijkl  
ABCDEFGHIJKL

16 char per line  
abcdefghijklmnop  
ABCDEFGHIJKLMN

21 chars per line  
abcdefghijklmnopqrstu  
ABCDEFGHIJKLMNQRSTU

24 characters per line  
abcdefghijklmnopqrstuwx  
ABCDEFGHIJKLMNQRSTUWX

32 characters per line  
abcdefghijklmnopqrstuwx123456  
ABCDEFGHIJKLMNQRSTUWXYZ123456

42 characters per line  
abcdefghijklmnopqrstuwx1234567890abcdef  
ABCDEFGHIJKLMNQRSTUWXYZ1234567890ABCDEF

Current Instrument Setup

Data Entry Mode: Enabled  
Start Test Sequence: DABACABA  
Display Prelim Rslt? Yes  
Display Third Digit? Yes  
Inhib Printer(Y/N)? No  
Display Volume? No  
Disable On Menfull? No  
# of Print Copies? 1  
Select Std (D/W/I)? Dry  
Standard Value? 0.000  
Standard Lot #? 29614080A2  
Standard Cyl #? 6  
Standard Expiration? 12/05/2016  
Oper No? 666666

Flow Cal. Date: 08/18/2015  
Slope 678  
Intercept -657244

IR Calibration Date: 08/18/2015  
3um 9um

0th Coef(*100):	-22405	-18664
1st Coef(*100):	270469	136654
2nd Coef(*100):	3040	1360
H2O adj(ng/l*10K):	699	443

\*\*\*\*\* Test End \*\*\*\*\*



Operator Signature  
CHARLES EDER

Remarks:

PRINT TEST

Form 106-18000




Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006495  
Location = TOXL      8164.13.00 06/09  
01/21/2016      13:42

Test	AC	Time
01 Room Air	0.000	13:43
02 *Subject Test	REF*	13:43
03 Room Air	0.000	13:44

\*Subject Test Refused

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female      Weight = 150  
Test = DUI      Cit = REFUSAL  
Dr. Lic. = ND/DIS821456  
Lot No = 29614080A2  
Cyl No = 6  
Expiration Date = 12/05/2016  
County = 08      Oper No. = 666666

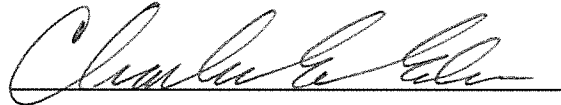

Operator Signature  
CHARLES EDER

Remarks:

REFUSAL

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006495  
Location = TOXL      8164.13.00 06/09  
01/21/2016      13:48

DIAGNOSTIC

Voltage/Current Test	Pass
RAM Test	Pass
EEPROM Checksum Test	Pass
Real Time Clock Test	Pass
DSP Test	Pass
Analytical Stability Test	Pass
Modem Test	Pass
Temp Regulation Test	Pass

Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: DIAGNOSTIC

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006495  
Location = TOXL      8164.13.00 06/09  
01/21/2016      14:07

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:08
02 Std. Sol.	0.000	14:09
03 Room Air	0.000	14:09
04 Std. Sol.	0.000	14:10
05 Room Air	0.000	14:10
06 Std. Sol.	0.000	14:11
07 Room Air	0.000	14:12

08 Sim Temp = 34.0°C

Simul Ser No = MP3067  
Std Sol No = MILLIQ H2O  
County = 08      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks:      ZERO AC  
                    0.000 AC



Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006495  
Location = TOXL      8164.13.00 06/09  
01/21/2016      14:21

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:22
02 Std. Sol.	0.019	14:23
03 Room Air	0.000	14:24
04 Std. Sol.	0.019	14:24
05 Room Air	0.000	14:25
06 Std. Sol.	0.019	14:26
07 Room Air	0.000	14:26

08 Sim Temp = 34.0°C

Simul Ser No = MP3066

Std Sol No = 15160

County = 08

Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks:

*Low AC*

*0.020 AC*

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006495  
Location = TOXL      8164.13.00 06/09  
01/21/2016      14:28

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:29
02 Std. Sol.	0.395	14:30
03 Room Air	0.000	14:31
04 Std. Sol.	0.394	14:31
05 Room Air	0.000	14:32
06 Std. Sol.	0.396	14:33
07 Room Air	0.000	14:34

08 Sim Temp = 34.0°C

Simul Ser No = MP3068  
Std Sol No = 14102  
County = 08      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: *HIGH AC*  
*0.400 AC*


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006495  
Location = TOXL      8164.13.00 06/09  
01/21/2016      14:36

Test	AC	Time
01 Room Air	0.000	14:37
02 *Subject Test	INT*	14:38
03 Room Air	0.000	14:38

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female      Weight = 150  
Test = DUI      Cit = 0.10AC+0.05ACE  
Dr. Lic. = ND/DIS821456  
Lot No = 29614080A2  
Cyl No = 6  
Expiration Date = 12/05/2016  
County = 08      Oper No. = 666666

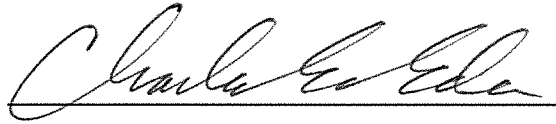

Operator Signature  
CHARLES EDER

Remarks: *INTERFERENT DETECT*

*0.10 AC + 0.05% ACETONE*


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-006495  
Location = TOXL                    8164.13.00 06/09  
01/21/2016                            14:40

Test	AC	Time
01 Room Air	0.000	14:41
02 *Subject Test	INT*	14:41
03 Room Air	0.000	14:42

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = 0.05 ACETONE  
Dr. Lic. = ND/DIS821456  
Lot No = 29614080A2  
Cyl No = 6  
Expiration Date = 12/05/2016  
County = 08                            Oper No. = 666666

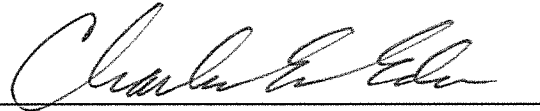

Operator Signature  
CHARLES EDER

Remarks: *INTERFERENT DETECT*  
*0.05% ACETONE*

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-006495  
Location = TOXL                8164.13.00 06/09  
01/21/2016                        14:43

Test	AC	Time
01 Room Air	0.000	14:44
02 *Subject Test	INT*	14:44
03 Room Air	0.000	14:45

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                Weight = 150  
Test = DUI                        Cit = 0.10 METHANOL  
Dr. Lic. = ND/DIS821456  
Lot No = 29614080A2  
Cyl No = 6  
Expiration Date = 12/05/2016  
County = 08                        Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: 0.10% METHANOL

*INTERFERENT DETECT*

Form 106-I8000


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-006495  
Location = TOXL                    8164.13.00 06/09  
01/21/2016                            14:45

Test	AC	Time
01 Room Air	0.000	14:46
02 *Subject Test	INT*	14:47
03 Room Air	0.000	14:47

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = 0.10 IPA  
Dr. Lic. = ND/DIS821456  
Lot No = 29614080A2  
Cyl No = 6  
Expiration Date = 12/05/2016  
County = 08                            Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: 0.10% ISOPROPANOL  
INTERFERENT DETECT

Form 106-I8000




Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-006495  
Location = TOXL                 8164.13.00 06/09  
01/21/2016                        14:56

Test	AC	Time
01 Room Air	0.000	14:57
02 Invalid Sample	X.XXX	14:57
03 Room Air	0.000	14:58

\*Invalid Test - Mouth Alcohol

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                 Weight = 150  
Test = DUI                        Cit = INVALID SAMPLE  
Dr. Lic. = ND/DIS821456  
Lot No = 29614080A2  
Cyl No = 6  
Expiration Date = 12/05/2016  
County = 08                        Oper No. = 666666

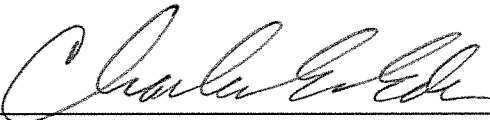

Operator Signature  
CHARLES EDER

Remarks:

*INVALID  
SAMPLE*

Form 106-I8000




Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-006495  
Location = TOXL                    8164.13.00 06/09  
01/21/2016                            15:09

Test	AC	Time
01 Room Air	0.000	15:10
02 *Subject Test	ABT*	15:10
03 Room Air	0.000	15:11

\*Sequence Aborted

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = SEQUENCE ABORT  
Dr. Lic. = ND/DIS821456  
Lot No = 29614080A2  
Cyl No = 6  
Expiration Date = 12/05/2016  
County = 08                            Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: SEQUENCE ABORTED

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-006495  
Location = TOXL                    8164.13.00 06/09  
01/21/2016                            15:13

Test	AC	Time
01 Room Air	0.000	15:13
02 Subject Test 1	0.053	15:14
03 Room Air	0.000	15:15
04 Reported AC	0.053	15:14

No RFI Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = ABA TEST  
Dr. Lic. = ND/DIS821456  
Lot No = 29614080A2  
Cyl No = 6  
Expiration Date = 12/05/2016  
County = 08                            Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks:

*ABA TEST*

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006495  
Location = TOXL      8164.13.00 06/09  
01/21/2016      15:17

Test	AC	Time
01 Diagnostic	OK	15:18
02 Room Air	0.000	15:19
03 Subject Test 1	0.040	15:19
04 Room Air	0.000	15:21
05 Std. Gas	0.087*	15:23
06 Room Air	0.000	15:23

\*Cal Check Out of Tolerance

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female      Weight = 150  
Test = DUI      Cit = CAL OUT OF TOL  
Dr. Lic. = ND/DIS821456  
Lot No = 29614080A2  
Cyl No = 6  
Expiration Date = 12/05/2016  
County = 08      Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

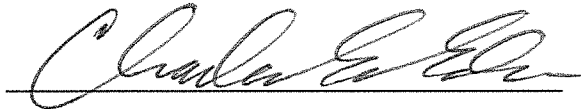

Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
OUT OF TOLERANCE

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006495  
Location = TOXL      8164.13.00 06/09  
01/21/2016      15:26

Test	AC	Time
01 Diagnostic	OK	15:27
02 Room Air	0.000	15:27
03 *Subject Test	RFI*	15:27
04 Room Air	0.000	15:28

\*Invalid Test  
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female      Weight = 150  
Test = DUI      Cit = RFI  
Dr. Lic. = ND/DIS821456  
Lot No = 29614080A2  
Cyl No = 6  
Expiration Date = 12/05/2016  
County = 08      Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

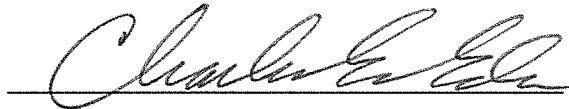

Operator Signature  
CHARLES EDER

Remarks:

RFI CHECK

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-006495  
Location = TOXL                    8164.13.00 06/09  
01/21/2016                            15:29

Test	AC	Time
01 Diagnostic	OK	15:30
02 Room Air	0.000	15:30
03 *Subject Test	0.000*	15:33
04 Room Air	0.000	15:36
05 Std. Gas	0.081	15:37
06 Room Air	0.000	15:38
07 *Subject Test	0.000*	15:41
08 Room Air	0.000	15:42


No RFI Detected

\*Deficient Sample - Value Printed was  
Highest Obtained

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = DEFICIENT SAMP  
Dr. Lic. = ND/DIS821456  
Lot No = 29614080A2  
Cyl No = 6  
Expiration Date = 12/05/2016  
County = 08                            Oper No. = 666666

20 minute waiting period ascertained? Y  
Clean Mouthpiece used and disposed? Y  
Clean Mouthpiece used and disposed? Y

I followed the Approved Method and the  
instructions displayed by the Intoxilyzer  
in conducting this test.

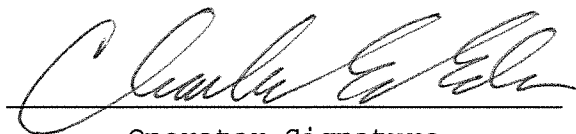

Operator Signature  
CHARLES EDER

Remarks: DEFICIENT  
SAMPLES

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006495  
Location = TOXL      8164.13.00 06/09  
01/21/2016      15:47

Test	AC	Time
01 Diagnostic	OK	15:47
02 Room Air	IPS*	15:48
03 Room Air	0.000	15:48

\*Improper Sample

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female      Weight = 150  
Test = DUI      Cit = IMPROPER SAMPLE  
Dr. Lic. = ND/DIS821456  
Lot No = 29614080A2  
Cyl No = 6  
Expiration Date = 12/05/2016  
County = 08      Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature  
CHARLES EDER

Remarks: *IMPROPER  
SAMPLE*

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006495  
Location = TOXL      8164.13.00 06/09  
01/21/2016      16:05

Test	AC	Time
01 Diagnostic	OK	16:06
02 Room Air	0.000	16:06
03 Subject Test 1	0.078*	16:06
04 Room Air	0.000	16:09
05 Std. Gas	0.080	16:10
06 Room Air	0.000	16:11
07 Subject Test 2	0.138*	16:12
08 Room Air	0.000	16:13

\*Difference Too Great

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female      Weight = 150  
Test = DUI      Cit = DIFF TOO GREAT  
Dr. Lic. = ND/DIS821456  
Lot No = 29614080A2  
Cyl No = 6  
Expiration Date = 12/05/2016  
County = 08      Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature  
CHARLES EDER

Remarks:

DIFFERENCE TOO  
CEL GREAT GREAT

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-006495  
Location = TOXL                    8164.13.00 06/09  
01/21/2016                            16:16


Test	AC	Time
01 Diagnostic	OK	16:17
02 Room Air	0.000	16:18
03 Subject Test 1	0.096	16:18
04 Room Air	0.000	16:20
05 Std. Gas	0.078	16:22
06 Room Air	0.000	16:23
07 Subject Test 2	0.096	16:24
08 Room Air	0.000	16:25
09 Reported AC	0.096	16:18

Difference OK  
No RFI Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = CMS TEST  
Dr. Lic. = ND/DIS821456  
Lot No = 29614080A2  
Cyl No = 6  
Expiration Date = 12/05/2016  
County = 08                            Oper No. = 666666


20 minute waiting period ascertained? Y  
Clean Mouthpiece used and disposed? Y  
Clean Mouthpiece used and disposed? Y

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

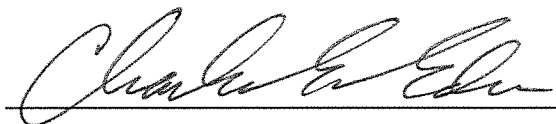

Operator Signature  
CHARLES EDER

Remarks:

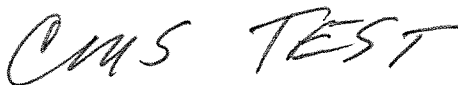



Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-006495  
Location = TOXL                8164.13.00 06/09  
01/21/2016                        16:33

Test	AC	Time
01 Room Air	0.000	16:34
02 Std. Gas	0.082	16:34
03 Room Air	0.000	16:35
04 Std. Gas	0.081	16:35
05 Room Air	0.000	16:36
06 Std. Gas	0.081	16:36
07 Room Air	0.000	16:37
08 Std. Gas	0.081	16:37
09 Room Air	0.000	16:38
10 Std. Gas	0.081	16:38
11 Room Air	0.000	16:38
12 Std. Gas	0.081	16:39
13 Room Air	0.000	16:39
14 Std. Gas	0.081	16:40
15 Room Air	0.000	16:40
16 Std. Gas	0.081	16:41
17 Room Air	0.000	16:41
18 Std. Gas	0.082	16:42
19 Room Air	0.000	16:42
20 Std. Gas	0.082	16:43
21 Room Air	0.000	16:43
22 Std. Gas	0.081	16:43
23 Room Air	0.000	16:44
24 Std. Gas	0.081	16:44
25 Room Air	0.000	16:45
26 Std. Gas	0.081	16:45
27 Room Air	0.000	16:46
28 Std. Gas	0.082	16:46
29 Room Air	0.000	16:47
30 Std. Gas	0.082	16:47
31 Room Air	0.000	16:48
32 Std. Gas	0.082	16:48
33 Room Air	0.000	16:49
34 Std. Gas	0.082	16:49
35 Room Air	0.000	16:50
36 Std. Gas	0.081	16:50
37 Room Air	0.000	16:50
38 Std. Gas	0.081	16:51
39 Room Air	0.000	16:51
40 Std. Gas	0.082	16:52
41 Room Air	0.000	16:52
42 Cal Check Stats:		
43 Average	0.0814	
44 Std Dev	0.0005	
45 Rel Std Dev	0.62%	

STABILITY TEST  
*Charles E. Edr*

Operator Signature  
CHARLES EDER

Remarks:

LOT: 29614 080A2  
EXP. DATE: 12/5/2016  
CYL: 006  
Form 106-I8000

STABILITY TEST  
Charles E. Eder