

80-004184

28 September 2011

I called the Grant Co SO's office to talk to Steven Bay about the Intoxilyzer 8000 Installation paperwork that he hasn't sent in yet. He is the Field Inspector for that location and I reminded him about this a month or so back and he told me he would look for the paperwork. He wasn't in today so I left a message for him to call me back.

I also called Terrick Chase, HP, to inform him not to run anymore tests on the 8000 until I hear from Steve. If Steve can't find the install paperwork, I would like him to do another installation and send that paperwork in to me as soon as he completes it.

*Deb Kashur*

Deb Kashur

18 Oct 11  
called Grant Co SO @ 8:40 AM  
left msg for SB. to call me back  
w/someone  
regarding the intox.

## Kashur, Deborah

---

**From:** Kashur, Deborah  
**Sent:** Wednesday, October 19, 2011 12:33 PM  
**To:** Bay, Steve E.  
**Subject:** RE: Intox 8000 Install  
**Attachments:** image001.gif

Thank you, Steve! Deb

---

**From:** Steve Bay [mailto:sbay@nd.gov]  
**Sent:** Wednesday, October 19, 2011 11:59 AM  
**To:** Kashur, Deborah  
**Subject:** Re: Intox 8000 Install

Deb, sorry it is my fault for not getting this done. I have printed out the documents. I am leaving the office at noon but will complete this on the 20th. If I have trouble I will call. Again I apologize.. Steve

-----Original Message-----

**From:** Kashur, Deborah  
**Date:** 10/18/2011 9:41:11 AM  
**To:** Bay, Steve E.  
**Cc:** Olson, Hope R.; Eder, Charles E.  
**Subject:** Intox 8000 Install

Steven,

I left a message this morning for you to call me at the Crime Laboratory. On 28 September 2011 I also called your office and left a message for you to call me back. This was in regards to Intoxilyzer 8000 installation paperwork that I have never received. I have spoken with you about this matter, but I am unsure of the date we talked. You told me you would look into it. I have also called Officer Chase telling him not to use that Intoxilyzer until a proper installation has been completed and sent to me.

If an installation was *not* performed since you received this intox on 03 May 2011 then any and every test that was run on this intox after that date will be invalid. If an installation was performed please send in the paperwork. If you cannot find the paperwork please perform another installation and send it in immediately. I have attached all necessary paperwork to complete an installation.

I'm afraid that if an installation will not or cannot be completed for any reason, then we have no choice but to

pull that instrument from your agency.

Please contact me regarding this matter.

Deb Kashur, Forensic Scientist

Office of Atty General, Crime Lab Division

2641 E Main Ave

Bismarck, ND 58501-5044

701-328-6159

Fax: 701-328-6185

email: [dkashur@nd.gov](mailto:dkashur@nd.gov)

**FREE Animations for Your Email - by IncrediMail!**

***Click Here!***



328-6185

# **GRANT COUNTY SHERIFF'S DEPT.**

## **FAX COVER SHEET**

**FAX NUMBER 701-622-3343**

**BUSINESS NUMBER 701-622-3331**

**106. MAIN STREET**

**CARSON, ND 58529**

**FROM:** Steve Grant Co

**TO:** ATTN- Deb

**RE:** Installation CK

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

**COVER + 4 PAGES**



### INTOXILYZER® 8000 INSTALLATION AND REPAIR CHECKOUT

OFFICE OF ATTORNEY GENERAL  
CRIME LABORATORY DIVISION  
Toxicology Section/Breath Alcohol Program  
SFN 59281 (5-2011)

Serial Number <i>80-004184</i>	Instrument Location <i>Grant Co</i>
Reason for Install/Repair <input type="checkbox"/> Install After Receiving From Crime Laboratory <input type="checkbox"/> Install After Location Change <input checked="" type="checkbox"/> Other (Specify) <u><i>Installation ck</i></u>	

Check When Done:

- 1. Surge Protector Installed/Properly Grounded.
- 2. Telephone Line Connected to Intoxilyzer® 8000.
- 3. Breath Tube Heated.
- 4. Enter Preliminary Data (i.e. Date, Time, DST (Y), and Location; Level 2, Function E).
- 5. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 6. Run Tests:
  - A. Print Test (Level 1, Function P).
  - B. ACA Test (Level 1, Function C).
  - C. RFI Test (Use CMS Mode; Key Radio During First Room Air).
- 7. Repair and/or Maintenance Performed (if Any): *None Needed*
- 8. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 9. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 10. Send the Following to the Crime Laboratory: Completed Intoxilyzer®8000 Installation and Repair Checkout (SFN59281, Form 104-G), Print Test, ACA Test, and RFI Test.

Field Inspector Signature <i>[Signature]</i>	Date <i>10-20-11</i>
Reviewed By (Crime Laboratory Use Only)	Date

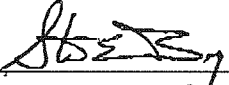
Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004184  
Location = TOXL      8164.13.00 06/09  
10/20/2011      10:50

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:50
02 Std. Gas	0.081	10:51
03 Room Air	0.000	10:51
04 Std. Gas	0.080	10:51
05 Room Air	0.000	10:52
06 Std. Gas	0.081	10:52
07 Room Air	0.000	10:53

Lot No = 15410080A1  
Cyl No = 12  
Exp Date = 07/01/2012  
County = 19      Oper No. = 039219

  
\_\_\_\_\_  
Operator Signature  
STEVEN E. BAY

Remarks:

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer  
North Dakota Model 8000 SN 80-004184  
Location = TOXL 8164.13.00 06/09  
10/20/2011 10:34

\*\*\*\*\* Printer Test \*\*\*\*\*

abcdefghijklmnopqrstuvwxy1234567890-=  
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&\*()\_+?

abcdefghijklmnopqrstuvwxy1234567890-=  
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&\*()\_+?

Current Instrument Setup


Data Entry Mode: Enabled  
Start Test Sequence: DABACABA  
Display Prelim Rslt? Yes  
Display Third Digit? Yes  
Inhib Printer(Y/N)? No  
Display Volume? No  
Disable On Memfull? Yes  
# of Print Copies? 1  
Select Std (D/W/I)? Dry  
Standard Value? 0.080  
Standard Lot #? 15410080A1  
Standard Cyl #? 12  
Standard Expiration? 07/01/2012  
Oper No? 039219

Flow Cal. Date: 07/30/2010  
Slope 804  
Intercept -925008

IR Calibration Date: 07/30/2010  
3um 9um

	3um	9um
0th Coef(*100):	-19612	-18975
1st Coef(*100):	270159	137378
2nd Coef(*100):	1174	963
H2O adj(mg/l*10k):	427	372

\*\*\*\*\* Printer Test End \*\*\*\*\*

  
Operator Signature  
STEVEN E. BAY

Remarks:


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004184  
Location = TOXL      8164.13.00 06/09  
10/20/2011      10:44

Test	AC	Time
01 Room Air	0.000	10:47
02 Subject Test 1	0.000	10:47
03 Room Air	0.000	10:48
04 Reported AC	0.000	10:47

No RFI Detected

Sub Name = BAY, STEVEN E  
Sub DOB = 03/22/1952  
Sub Sex = Male      Weight = 175  
Test = OTH      Cit = NA  
Dr. Lic. = ND/BAY527877  
Lot No = 15410080A1  
Cyl No = 12  
Expiration Date = 07/01/2012  
County = 19      Oper No. = 039219

  
\_\_\_\_\_  
Operator Signature  
STEVEN E. BAY

Remarks:





### INTOXILYZER® 8000 INSTALLATION AND REPAIR CHECKOUT

OFFICE OF ATTORNEY GENERAL  
CRIME LABORATORY DIVISION  
Toxicology Section/Breath Alcohol Program  
SFN 59281 (5-2011)

Serial Number <i>80-004184</i>	Instrument Location <i>Grant Co</i>
Reason for Install/Repair <input type="checkbox"/> Install After Receiving From Crime Laboratory <input type="checkbox"/> Install After Location Change <input checked="" type="checkbox"/> Other (Specify) <u><i>Installation ck</i></u>	

Check When Done:

- 1. Surge Protector Installed/Properly Grounded.
- 2. Telephone Line Connected to Intoxilyzer® 8000.
- 3. Breath Tube Heated.
- 4. Enter Preliminary Data (i.e. Date, Time, DST (Y), and Location; Level 2, Function E).
- 5. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 6. Run Tests:
  - A. Print Test (Level 1, Function P).
  - B. ACA Test (Level 1, Function C).
  - C. RFI Test (Use CMS Mode; Key Radio During First Room Air).
- 7. Repair and/or Maintenance Performed (if Any): *None Needed*
- 8. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 9. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 10. Send the Following to the Crime Laboratory: Completed Intoxilyzer®8000 Installation and Repair Checkout (SFN59281, Form 104-G), Print Test, ACA Test, and RFI Test.

Field Inspector Signature <i>[Signature]</i>	Date <i>10-20-11</i>
Reviewed By (Crime Laboratory Use Only)	Date

↑  
 write new date RFI is run please. This is the date the testing is complete.



Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

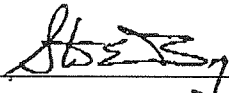
CMI, Inc. Intoxilyzer Alcohol Analyzer  
North Dakota Model 8000 SN 80-004184  
Location = TOXL ~~XL~~ CARS 8164.13.00 06/09  
10/20/2011 10:50



DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:50
02 Std. Gas	0.081	10:51
03 Room Air	0.000	10:51
04 Std. Gas	0.080	10:51
05 Room Air	0.000	10:52
06 Std. Gas	0.081	10:52
07 Room Air	0.000	10:53

Lot No = 15410080A1  
Cyl No = 12  
Exp Date = 07/01/2012  
County = 19 Oper No. = 039219

  
Operator Signature  
STEVEN E. BAY

Remarks:

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer  
North Dakota Model 8000 SN 80-004184  
Location = ~~TOX~~ CARS 8164.13.00 06/09  
10/20/2011 10:44

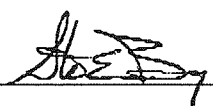


Test	AC	Time
01 Room Air	0.000	10:47
02 Subject Test 1	0.000	10:47
03 Room Air	0.000	10:48
04 Reported AC	0.000	10:47

*RFI test needs to be run in CMS mode - Esc Esc Start Test. For subj. name just type in "test" or simply "enter" through that information. Key radio during the "room air".*

No RFI Detected

Sub Name = BAY, STEVEN E  
Sub DOB = 03/22/1952  
Sub Sex = Male Weight = 175  
Test = OTH Cit = NA  
Dr. Lic. = ND/BAY527877  
Lot No = 15410080A1  
Cyl No = 12  
Expiration Date = 07/01/2012  
County = 19 Oper No. = 039219

  
\_\_\_\_\_  
Operator Signature  
STEVEN E. BAY

Remarks:

SN 80-004184 DK

Kashur, Deborah

To: Bay, Steve E.  
Subject: FW:  
Attachments: doc02425420111021153142.pdf

Steve,

Thanks for faxing me the installation paperwork. I would have responded sooner except that I was in Jamestown court all day yesterday.

There are a few things that need to be fixed. I've attached your install checkout so you can see what I'm talking about.

1. You had checked box 4 for entering the preliminary data. However, you didn't change the location code. Every test record still has - Location = TOXL, when it should be - Location = CARS. That is a big court issue and needs to be changed immediately. Remember to get into that level the password is - *DK 26 Apr 12*  
*- crossed out password*

2. Under #6, C you ran an RFI test, but it looks like you didn't key your radio near the breath hose or you didn't use a radio at all? And it was run in the wrong mode. You ran it by entering - Esc Esc Enter - when it should have been run in CMS mode by entering - Esc Esc Start Test button.

To be honest, it would probably be less confusing and easier if you completely started over with your agency location code and a new installation checkout. However, if you don't want to start over you can do this:

Change the location code in the 8000, on every test record cross out TOXL and write CARS and initial, rerun the RFI test in the correct mode and make sure the radio interferes with the test so it is aborted, then cross off the date of 10-20-11 and rewrite the date with the date you ran the RFI test.

Please call me if you have questions. I wanted to email you this info instead of calling you so you can look things over first.

Thank you, Steve!

Deb Kashur, Forensic Scientist  
Office of Atty General, Crime Lab Division  
2641 E Main Ave  
Bismarck, ND 58501-5044  
701-328-6159  
Fax: 701-328-6185  
email: [dkashur@nd.gov](mailto:dkashur@nd.gov)

-----Original Message-----

From: ndagweb@nd.gov [mailto:ndagweb@nd.gov]  
Sent: Friday, October 21, 2011 3:32 PM  
To: Kashur, Deborah  
Subject:

-----  
TASKalfa 420i  
[00:c0:ee:7a:2b:ff]  
-----