



ETHANOL GAS STANDARD CYLINDER REPORT

NORTH DAKOTA OFFICE OF ATTORNEY GENERAL

CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM

SFN 59282 (08/2018)

Chemical Test Operator Name (Print) <i>Kohler, Jedediah Scott</i>		
Location <i>Mott</i>	Intoxilyzer® Serial Number <i>80-007091</i>	
Gas Lot Number <i>05620080A1</i>	Gas Cylinder Number <i>41</i>	Gas Expiration Date <i>04/05/2022</i>

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C).
Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
 - A. *0.079* AC
 - B. *0.078* AC
 - C. *0.079* AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 5. Send the Following to the Crime Laboratory:
 - A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).
 - B. ACA Test Record.

Chemical Test Operator Signature <i>[Signature]</i>	Date <i>11/8/20</i>
Reviewed By (Crime Laboratory Use Only) <i>[Signature]</i>	Date <i>11/17/20</i>


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-007091
Location = MOTT 8164.16.00 09/18
11/08/2020 13:12

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:14
02 Std. Gas	0.079	13:14
03 Room Air	0.000	13:15
04 Std. Gas	0.078	13:15
05 Room Air	0.000	13:16
06 Std. Gas	0.079	13:16
07 Room Air	0.000	13:16

Lot No = 05620080A1
Cyl No = 41
Exp Date = 04/05/2022
County = 21 Oper No. = 132666

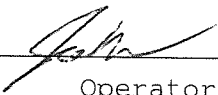

Operator Signature
JEDEDIAH KOHLER

Remarks: *ACA Function Test*

Form 106-I8000