

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-006669

Location: TOXL

- Check When Done
- A. Is the breath tube heated? ✓
 - B. Is the diagnostic check complete? ✓
 - C. Is the time, date, and year correct (re-set if necessary)? ✓
 - D. Print test. (Attach test record.) ✓
 - E. Low AC. Use ≤ 0.03 AC in ACA mode. (Attach test record.)
Sim SN MP3003 Lot No. 201610A Rep AC 0.020 ✓
 - F. High AC. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN MP3063 Lot No. 16102 Rep AC 0.400 ✓
 - G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.
(Attach test record.)
Sim SN DR5135 Lot No. ICS4 Solution 0.10AC + 0.05% ACETONE ✓
 - H. RFI check. Run CMS mode. Key radio during CMS test.
(Attach test record.) ✓
 - I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) 0.080 AC

Lot No. 1731608042 Cyl No. 026 Exp. Date 8/5/2018

Test 1 <u>0.081</u>	Test 4 <u>0.082</u>	Test 7 <u>0.082</u>
Test 2 <u>0.081</u>	Test 5 <u>0.082</u>	Test 8 <u>0.082</u>
Test 3 <u>0.081</u>	Test 6 <u>0.082</u>	Test 9 <u>0.082</u>

Average 0.082

J. Tank pressure, Level 3 Func D Sub G P: 325 psi Reg ~350 psi

K. Remarks/Maintenance record: ANNUAL INSPECTION

4.6.2017
Date

11 Apr 17
Date

Charles E. Ehr
Field Inspector's Signature

Dub Shanaver
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006669
Location = TOXL 8164.14.00 09/16
04/06/2017 12:59

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:59
02 Std. Sol.	0.017	13:00
03 Room Air	0.000	13:01
04 Std. Sol.	0.017	13:02
05 Room Air	0.000	13:02
06 Std. Sol.	0.017	13:03
07 Room Air	0.000	13:03

08 Sim Temp = 34.0°C

Simul Ser No = MP3003
Std Sol No = 201610A
County = 08 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: *Low AC*
0.017 0.020 AC
CEE
Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006669
Location = TOXL 8164.14.00 09/16
04/06/2017 13:14

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:14
02 Std. Sol.	0.396	13:15
03 Room Air	0.000	13:16
04 Std. Sol.	0.397	13:16
05 Room Air	0.000	13:17
06 Std. Sol.	0.396	13:18
07 Room Air	0.000	13:18

08 Sim Temp = 34.0°C

Simul Ser No = MP3063
Std Sol No = 16102
County = 08 Oper No. = 666666

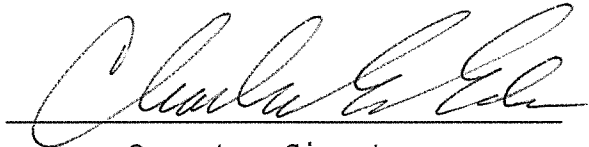

Operator Signature
CHARLES EDER

Remarks: HIGH AC
0.400 AC


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006669
Location = TOXL 8164.14.00 09/16
04/06/2017 13:25

Test	AC	Time
01 Room Air	0.000	13:26
02 *Subject Test	INT*	13:26
03 Room Air	0.000	13:27

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = INTERFERENT
Dr. Lic. = ND/DIS821456
Lot No = 17316080A2
Cyl No = 26
Expiration Date = 08/05/2018
County = 08 Oper No. = 666666

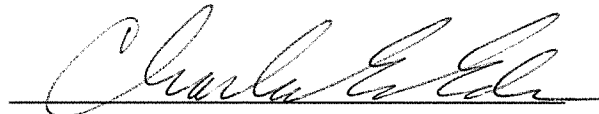

Operator Signature
CHARLES EDER

Remarks:

Interferent
0.10 AC + 0.05% ACETONE

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006669
Location = TOXL 8164.14.00 09/16
04/06/2017 13:27

Test	AC	Time
01 Diagnostic	OK	13:28
02 Room Air	0.000	13:28
03 *Subject Test	RFI*	13:29
04 Room Air	0.000	13:29

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = RFI CHECK
Dr. Lic. = ND/DIS821456
Lot No = 17316080A2
Cyl No = 26
Expiration Date = 08/05/2018
County = 08 Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

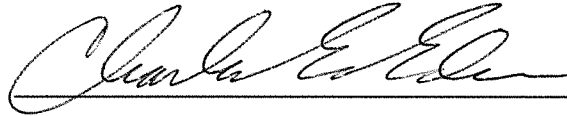

Operator Signature
CHARLES EDER

Remarks:

RFI CHECK


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006669
Location = TOXL 8164.14.00 09/16
04/06/2017 13:30

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:30
02 Std. Gas	0.081	13:31
03 Room Air	0.000	13:31
04 Std. Gas	0.081	13:32
05 Room Air	0.000	13:32
06 Std. Gas	0.081	13:32
07 Room Air	0.000	13:33

Lot No = 17316080A2
Cyl No = 26
Exp Date = 08/05/2018
County = 08 Oper No. = 666666

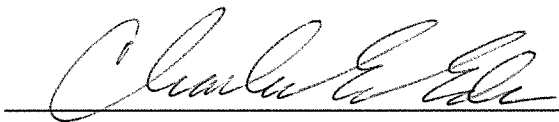

Operator Signature
CHARLES EDER

Remarks:

0.080 AC
CALIBRATION CHECK

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006669
Location = TOXL 8164.14.00 09/16
04/06/2017 13:37

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:37
02 Std. Gas	0.082	13:38
03 Room Air	0.000	13:38
04 Std. Gas	0.082	13:38
05 Room Air	0.000	13:39
06 Std. Gas	0.082	13:39
07 Room Air	0.000	13:40

Lot No = 17316080A2
Cyl No = 26
Exp Date = 08/05/2018
County = 08 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
0.080 AC