

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-006677

Location: TOXL

- Check When Done
- A. Is the breath tube heated? ✓
- B. Is the diagnostic check complete? ✓
- C. Is the time, date, and year correct (re-set if necessary)? ✓
- D. Print test. (Attach test record.) ✓
- E. Low AC. Use  $\leq 0.03$  AC in ACA mode. (Attach test record.)  
Sim SN MP3003 Lot No. 201610A Rep AC 0.020 ✓
- F. High AC. Use  $\geq 0.25$  AC in ACA mode. (Attach test record.)  
Sim SN MP3063 Lot No. 16102 Rep AC 0.400 ✓
- G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.  
(Attach test record.) ✓  
Sim SN DR5135 Lot No. ICS4 Solution 0.10 AC + 0.05% ACETONE
- H. RFI check. Run CMS mode. Key radio during CMS test.  
(Attach test record.) ✓
- I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.  
(Attach test record.) 0.080 AC

Lot No. 17316080A2 Cyl No. 26 Exp. Date 8/5/2018

Test 1 <u>0.082</u>	Test 4 <u>0.082</u>	Test 7 <u>0.082</u>
Test 2 <u>0.082</u>	Test 5 <u>0.082</u>	Test 8 <u>0.082</u>
Test 3 <u>0.082</u>	Test 6 <u>0.083</u>	Test 9 <u>0.082</u>

Average 0.082

J. Tank pressure, Level 3 Func D Sub G

P: 178 psi Reg 225 psi

K. Remarks/Maintenance record: ANNUAL INSPECTION

4-12-2017  
Date

12 Apr 17  
Date

Charles E. Ed  
Field Inspector's Signature

Dib Shanaver  
Reviewed by

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-006677  
Location = TOXL                    8164.14.00 09/16  
04/12/2017                            14:13

\*\*\*\*\* Printer Test \*\*\*\*\*

abcdefghijklmnopqrstuvwxyz1234567890-=\_|  
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&^\*()\_+?

abcdefghijklmnopqrstuvwxyz1234567890-=\_|  
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&^\*()\_+?

Current Instrument Setup

Data Entry Mode:            Enabled  
Start Test Sequence:        DABACABA  
Display Prelim Rslt?        Yes  
Display Third Digit?        Yes  
Inhib Printer(Y/N)?        No  
Display Volume?            No  
Disable On Memfull?        Yes  
# of Print Copies?          1  
Select Std (D/W/I)?        Dry  
Standard Value?            0.080  
Standard Lot #?            17316080A2  
Standard Cyl #?            26  
Standard Expiration?        08/05/2018  
Oper No?                    666666

Flow Cal. Date:            08/24/2016  
Slope                        696  
Intercept                   -655873

IR Calibration Date:        08/24/2016  
                                  3um                        , 9um

-----  
0th Coef(\*100):            -12424                    -22862  
1st Coef(\*100):            263434                    131353  
2nd Coef(\*100):            5070                      1477  
H2O adj(mg/l\*10k):        569                        590

\*\*\*\*\* Printer Test End \*\*\*\*\*



Operator Signature  
CHARLES EDER

Remarks:                    PRINT TEST

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006677  
Location = TOXL      8164.14.00 09/16  
04/12/2017      14:16

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:16
02 Std. Sol.	0.017	14:17
03 Room Air	0.000	14:18
04 Std. Sol.	0.017	14:18
05 Room Air	0.000	14:19
06 Std. Sol.	0.018	14:20
07 Room Air	0.000	14:20

08 Sim Temp = 34.0°C

Simul Ser No = MP3003  
Std Sol No = 201610A  
County = 08      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: *Low AC 0.0204c*

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006677  
Location = TOXL      8164.14.00 09/16  
04/12/2017      14:22

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:22
02 Std. Sol.	0.393	14:23
03 Room Air	0.000	14:24
04 Std. Sol.	0.395	14:24
05 Room Air	0.000	14:25
06 Std. Sol.	0.395	14:26
07 Room Air	0.000	14:26

08 Sim Temp = 34.0°C

Simul Ser No = MP3063  
Std Sol No = 16102  
County = 08      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks:      HIGH AC  
                    0.400AC


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-006677  
Location = TOXL                    8164.14.00 09/16  
04/12/2017                            14:29

Test	AC	Time
01 Room Air	0.000	14:29
02 *Subject Test	INT*	14:30
03 Room Air	0.000	14:30

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = INTERFERENT  
Dr. Lic. = ND/DIS821456  
Lot No = 17316080A2  
Cyl No = 26  
Expiration Date = 08/05/2018  
County = 08                            Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: *Interferent*  
*0.10AC + 0.05% ACETONE*

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-006677  
Location = TOXL                    8164.14.00 09/16  
04/12/2017                            14:31

Test	AC	Time
01 Diagnostic	OK	14:31
02 Room Air	0.000	14:32
03 *Subject Test	RFI*	14:32
04 Room Air	0.000	14:32

\*Invalid Test  
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = RFI CHECK  
Dr. Lic. = ND/DIS821456  
Lot No = 17316080A2  
Cyl No = 26  
Expiration Date = 08/05/2018  
County = 08                            Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

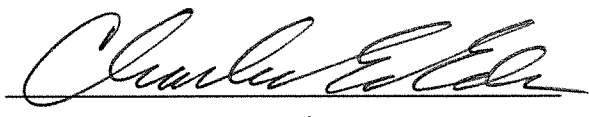

Operator Signature  
CHARLES EDER

Remarks: *RFI CHECK*


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006677  
Location = TOXL      8164.14.00 09/16  
04/12/2017      14:33

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:34
02 Std. Gas	0.082	14:34
03 Room Air	0.000	14:34
04 Std. Gas	0.082	14:35
05 Room Air	0.000	14:35
06 Std. Gas	0.082	14:36
07 Room Air	0.000	14:36

Lot No = 17316080A2  
Cyl No = 26  
Exp Date = 08/05/2018  
County = 08      Oper No. = 666666

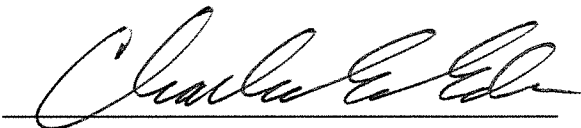

Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
CCE  
0.082 AC  
0

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006677  
Location = TOXL      8164.14.00 09/16  
04/12/2017      14:36

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:37
02 Std. Gas	0.082	14:37
03 Room Air	0.000	14:38
04 Std. Gas	0.082	14:38
05 Room Air	0.000	14:39
06 Std. Gas	0.083	14:39
07 Room Air	0.000	14:39

Lot No = 17316080A2  
Cyl No = 26  
Exp Date = 08/05/2018  
County = 08      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
0.080 AC




Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-006677  
Location = TOXL                    8164.14.00 09/16  
04/12/2017                            14:40

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:40
02 Std. Gas	0.082	14:41
03 Room Air	0.000	14:41
04 Std. Gas	0.082	14:41
05 Room Air	0.000	14:42
06 Std. Gas	0.082	14:42
07 Room Air	0.000	14:43

Lot No = 17316080A2  
Cyl No = 26  
Exp Date = 08/05/2018  
County = 08                            Oper No. = 666666

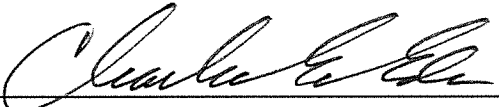  
Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
0.080 AC