

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-006685

Location: TOXL

- |  | Check When Done |
|--|-----------------|
| A. Is the breath tube heated?  | <u>X</u>        |
| B. Is the diagnostic check complete?   | <u>X</u>        |
| C. Is the time, date, and year correct (re-set if necessary)?  | <u>X</u>        |
| D. Print test. (Attach test record.)   | <u>X</u>        |
| E. Low AC. Use $\leq 0.03$ AC in ACA mode. (Attach test record.)<br>Sim SN <u>DR3846</u> Lot No. <u>16200</u> Rep AC <u>0.020</u>  | <u>X</u>        |
| F. High AC. Use $\geq 0.25$ AC in ACA mode. (Attach test record.)<br>Sim SN <u>DR7111</u> Lot No. <u>16060</u> Rep AC <u>0.300</u>   | <u>X</u>        |
| G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.<br>(Attach test record.)<br>Sim SN <u>G10329</u> Lot No. <u>1CS4</u> Solution <u>0.10AC + 0.05%</u><br><span style="margin-left: 150px;">Acetone</span> | <u>X</u>        |
| H. RFI check. Run CMS mode. Key radio during CMS test.<br>(Attach test record.)  | <u>X</u>        |
| I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.<br>(Attach test record.) <u>0.080</u> AC  |                 |

Lot No. 17316080A2 Cyl No. 27 Exp. Date 5 Aug 18

Test 1 <u>0.080</u>	Test 4 <u>0.081</u>	Test 7 <u>0.081</u>
Test 2 <u>0.081</u>	Test 5 <u>0.081</u>	Test 8 <u>0.081</u>
Test 3 <u>0.081</u>	Test 6 <u>0.081</u>	Test 9 <u>0.081</u>

Average 0.081

J. Tank pressure, Level 3 Func D Sub G P: 1005 psi Reg 1000 psi

K. Remarks/Maintenance record: Annual Inspection.

13 Apr 17  
Date

13 APR 17  
Date

Roberta Yeger-Timmis  
Field Inspector's Signature

Charles E. Elie  
Reviewed by

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer           Alcohol Analyzer  
North Dakota Model 8000       SN 80-006685  
Location = TOXL                8164.14.00 09/16  
04/13/2017                      14:08

\*\*\*\*\* Printer Test \*\*\*\*\*

abcdefghijklmnopqrstuvwxy1234567890-=  
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#%^&\*()\_+?

abcdefghijklmnopqrstuvwxy1234567890-=  
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#%^&\*()\_+?

Current Instrument Setup

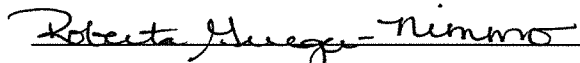
Data Entry Mode:            Enabled  
Start Test Sequence:        DABACABA  
Display Prelim Rslt?        Yes  
Display Third Digit?         Yes  
Inhib Printer(Y/N)?         No  
Display Volume?             No  
Disable On Memfull?         Yes  
# of Print Copies?          1  
Select Std (D/W/I)?         Dry  
Standard Value?             0.080  
Standard Lot #?             17316080A2  
Standard Cyl #?             27  
Standard Expiration?        08/05/2018  
Oper No?                    888888

Flow Cal. Date:             09/06/2016  
Slope                       705  
Intercept                  -687370

IR Calibration Date:        09/06/2016

	3um	9um
0th Coef(*100):	-14155	-22915
1st Coef(*100):	271436	135296
2nd Coef(*100):	2016	1251
H2O adj(mg/l*10k):	372	505

\*\*\*\*\* Printer Test End \*\*\*\*\*

  
Operator Signature  
ROBERTA GRIEGER-NIMMO

Remarks:

Print test

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

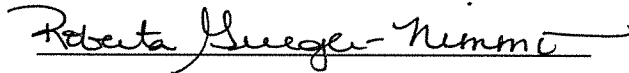
CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006685  
Location = TOXL      8164.14.00 09/16  
04/13/2017      14:10

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:10
02 Std. Sol.	0.020	14:11
03 Room Air	0.000	14:12
04 Std. Sol.	0.020	14:13
05 Room Air	0.000	14:13
06 Std. Sol.	0.020	14:14
07 Room Air	0.000	14:15

08 Sim Temp = 34.0°C

Simul Ser No = DR3846  
Std Sol No = 16200  
County = 08      Oper No. = 888888

  
Operator Signature  
ROBERTA GRIEGER-NIMMO

Remarks: Low AC  
0.020 AC

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006685  
Location = TOXL      8164.14.00 09/16  
04/13/2017      14:17

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:17
02 Std. Sol.	0.306	14:19
03 Room Air	0.000	14:19
04 Std. Sol.	0.304	14:20
05 Room Air	0.000	14:21
06 Std. Sol.	0.303	14:22
07 Room Air	0.000	14:22

08 Sim Temp = 34.0°C

Simul Ser No = DR7111

Std Sol No = 16060

County = 08

Oper No. = 888888

Roberta Grieger-Nimmo

Operator Signature  
ROBERTA GRIEGER-NIMMO

Remarks: High AC  
0.300 AC

Form 106-I8000

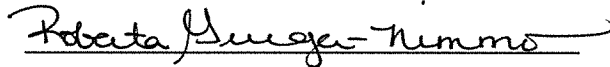
Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-006685  
Location = TOXL                    8164.14.00 09/16  
04/13/2017                            14:24

Test	AC	Time
01 Room Air	0.000	14:25
02 *Subject Test	INT*	14:27
03 Room Air	0.000	14:28

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = 08  
Dr. Lic. = ND/DIS821456  
Lot No = 17316080A2  
Cyl No = 27  
Expiration Date = 08/05/2018  
County = 08                            Oper No. = 888888



Operator Signature  
ROBERTA GRIEGER-NIMMO

Remarks: Interferent Check  
0.10 AC + 0.05% Acetone



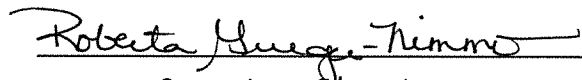
Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006685  
Location = TOXL      8164.14.00 09/16  
04/13/2017      14:32

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:33
02 Std. Gas	0.080	14:33
03 Room Air	0.000	14:34
04 Std. Gas	0.081	14:34
05 Room Air	0.000	14:35
06 Std. Gas	0.081	14:35
07 Room Air	0.000	14:36

Lot No = 17316080A2  
Cyl No = 27  
Exp Date = 08/05/2018  
County = 08      Oper No. = 888888

  
Operator Signature  
ROBERTA GRIEGER-NIMMO

Remarks: Calibration Check  
0.080 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006685  
Location = TOXL      8164.14.00 09/16  
04/13/2017      14:36

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:37
02 Std. Gas	0.081	14:37
03 Room Air	0.000	14:38
04 Std. Gas	0.081	14:38
05 Room Air	0.000	14:39
06 Std. Gas	0.081	14:39
07 Room Air	0.000	14:40

Lot No = 17316080A2  
Cyl No = 27  
Exp Date = 08/05/2018  
County = 08      Oper No. = 888888

Roberta Grieger-Nimmo  
Operator Signature  
ROBERTA GRIEGER-NIMMO

Remarks: Calibration Check  
0.080AC



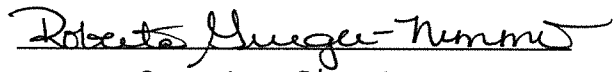
Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006685  
Location = TOXL      8164.14.00 09/16  
04/13/2017      14:40

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:40
02 Std. Gas	0.081	14:41
03 Room Air	0.000	14:42
04 Std. Gas	0.081	14:42
05 Room Air	0.000	14:43
06 Std. Gas	0.081	14:43
07 Room Air	0.000	14:44

Lot No = 17316080A2  
Cyl No = 27  
Exp Date = 08/05/2018  
County = 08      Oper No. = 888888

  
Operator Signature  
ROBERTA GRIEGER-NIMMO

Remarks: Calibration Check  
0.080 AC