

INTOXILYZER® 8000 INSPECTION

Intoxilyzer® 8000 Serial Number: 80-006690

Location: RUGB
HEART OF AMERICA

A. General Setup and Checks:

1. Diagnostic Tests Pass and Instrument in Ready Mode CORRECTIONS-RUGBY
2. Breath tube heated
3. Date, time and location code (Level 2,E). Re-set if necessary.
4. Print test (Level 1,P). Sign and attach test record.
5. Tank monitor (Level 3,D,G). Display and Regulator within 50 psi.
 - a. Display: 352 psi Regulator: 350 psi

B. Tests (Sign and attach test records):

- Configure simulator for the following tests (Level 1,S).
1. Wet Calibration Check
 - a. Low AC (Level 1,C): Known Value ≤ 0.03 AC: 0.020 AC
Sim. Ser #: DR3377 Lot #: 201610A Exp. Date: 10-4-18
 - b. High AC (Level 1,C): Known Value ≥ 0.25 AC: 0.300 AC
Sim. Ser #: DR5191 Lot #: 17350 Exp. Date: 10-11-19

Configure dry gas standard for the following tests (Level 1,S).

2. Interferent Check (Level 1,B): Known Value: 0.10 AC + 0.05% Acetone
Sim. Ser #: DR3451 Lot #: ICS#5 Exp. Date: 5-15-19
 Display reads "Interferent Detect"
3. RFI Check (CMS Mode) Display reads "RFI Detect"
4. Dry Calibration Check (Level 1,C): Known Value 0.080 AC
Gas Cyl Lot #: 34917080A3 Cyl #: 5 Exp. Date: 2/5/20

Test 1 <u>0.081</u>	Test 4 <u>0.080</u>	Test 7 <u>0.081</u>
Test 2 <u>0.081</u>	Test 5 <u>0.081</u>	Test 8 <u>0.081</u>
Test 3 <u>0.080</u>	Test 6 <u>0.081</u>	Test 9 <u>0.080</u>
Average <u>0.081</u>		

C. Remarks/Maintenance: ANNUAL INSPECTION

Instrument is acceptable to be used in the field.

Charles E. Edr
Breath Analyst Signature

5-17-18
Date

Roberta Mueggner-Timmer
Reviewed by

24 May 18
Date

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006690
Location = RUGB 8164.14.00 09/16
05/17/2018 10:39

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:39
02 Std. Sol.	0.017	10:40
03 Room Air	0.000	10:41
04 Std. Sol.	0.017	10:41
05 Room Air	0.000	10:42
06 Std. Sol.	0.017	10:43
07 Room Air	0.000	10:43

08 Sim Temp = 34.0°C

Simul Ser No = DR3377
Std Sol No = 201610A
County = 35 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: *Low AC*
 0.020 AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006690
Location = RUGB 8164.14.00 09/16
05/17/2018 10:44

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:45
02 Std. Sol.	0.291	10:46
03 Room Air	0.000	10:46
04 Std. Sol.	0.295	10:47
05 Room Air	0.000	10:47
06 Std. Sol.	0.295	10:48
07 Room Air	0.000	10:49

08 Sim Temp = 34.0°C

Simul Ser No = DR5191
Std Sol No = 17350
County = 35 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: *HIGH AC*
0.300 AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006690
Location = RUGB 8164.14.00 09/16
05/17/2018 10:21

Test	AC	Time
01 Room Air	0.000	10:22
02 *Subject Test	INT*	10:22
03 Room Air	0.000	10:23

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = INTERFERENT
Dr. Lic. = ND/DIS821456
Lot No = 34917080A3
Cyl No = 5
Expiration Date = 02/05/2020
County = 35 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

Interferent Check
0.10 AC + 0.05% ACETONE

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006690
Location = RUGB 8164.14.00 09/16
05/17/2018 10:18

Test	AC	Time
01 Diagnostic	OK	10:19
02 Room Air	0.000	10:19
03 *Subject Test	RFI*	10:19
04 Room Air	0.000	10:20

*Invalid Test
 Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = RFI CHECK
Dr. Lic. = ND/DIS821456
Lot No = 34917080A3
Cyl No = 5
Expiration Date = 02/05/2020
County = 35 Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

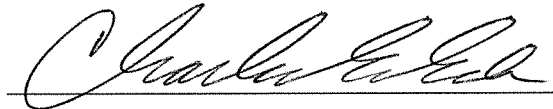

Operator Signature
CHARLES EDER

Remarks:

RFI CHECK

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006690
Location = RUGB 8164.14.00 09/16
05/17/2018 10:23

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:23
02 Std. Gas	0.081	10:24
03 Room Air	0.000	10:24
04 Std. Gas	0.081	10:25
05 Room Air	0.000	10:25
06 Std. Gas	0.080	10:26
07 Room Air	0.000	10:26

Lot No = 34917080A3
Cyl No = 5
Exp Date = 02/05/2020
County = 35 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
0.080 AC


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006690
Location = RUGB 8164.14.00 09/16
05/17/2018 10:27

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:28
02 Std. Gas	0.080	10:28
03 Room Air	0.000	10:29
04 Std. Gas	0.081	10:29
05 Room Air	0.000	10:30
06 Std. Gas	0.081	10:30
07 Room Air	0.000	10:30

Lot No = 34917080A3
Cyl No = 5
Exp Date = 02/05/2020
County = 35 Oper No. = 666666

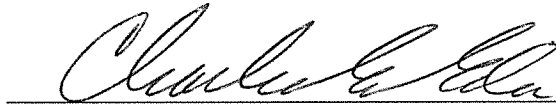

Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
0.080AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006690
Location = RUGB 8164.14.00 09/16
05/17/2018 10:49

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:50
02 Std. Gas	0.081	10:50
03 Room Air	0.000	10:51
04 Std. Gas	0.081	10:51
05 Room Air	0.000	10:52
06 Std. Gas	0.080	10:52
07 Room Air	0.000	10:52

Lot No = 34917080A3
Cyl No = 5
Exp Date = 02/05/2020
County = 35 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
0.080 AC