



# ETHANOL BREATH STANDARD CYLINDER REPORT

OFFICE OF ATTORNEY GENERAL  
CRIME LABORATORY DIVISION  
Toxicology Section/Breath Alcohol Program  
SFN 59282 (5-2011)

Chemical Test Operator Name (Print) <u>Pinski, Gregory E.</u>		
Location <u>Burke County SO.</u>	Intoxilyzer® Serial Number <u>80-006673</u>	
Gas Lot Number <u>13518080AL6</u>	Gas Cylinder Number <u>30</u>	Gas Expiration Date <u>08/05/2020</u>

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C).  
Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
  - A. 0.083 AC
  - B. 0.082 AC
  - C. 0.082 AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 5. Send the Following to the Crime Laboratory:
  - A. Ethanol Breath Standard Cylinder Report (SFN59282, Form 105-G).
  - B. ACA Test Record.

Chemical Test Operator Signature <u>[Signature]</u>	Date <u>12/16/2019</u>
Reviewed By (Crime Laboratory Use Only) <u>[Signature]</u>	Date <u>12.18.2019</u>


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006673  
Location = BURK      8164.14.00 09/16  
12/16/2019      14:19

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:20
02 Std. Gas	0.083	14:21
03 Room Air	0.000	14:21
04 Std. Gas	0.082	14:21
05 Room Air	0.000	14:22
06 Std. Gas	0.082	14:22
07 Room Air	0.000	14:23

Lot No = 13518080A6  
Cyl No = 30  
Exp Date = 08/05/2020  
County = 07      Oper No. = 132567

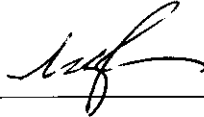

Operator Signature  
GREGORY PINSKI

Remarks:

Form 106-I8000