



ETHANOL GAS STANDARD CYLINDER REPORT
 NORTH DAKOTA OFFICE OF ATTORNEY GENERAL
 CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM
 SFN 59282 (08/2018)

Chemical Test Operator Name (Print) <i>Rustebakke, Dona</i>		
Location <i>VALL</i>	Intoxilyzer® Serial Number <i>80-006679</i>	
Gas Lot Number <i>02621080A1</i>	Gas Cylinder Number <i>36</i>	Gas Expiration Date <i>3/5/2023</i>

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C). *RECEIVED ACA TEST (DRY CAL CHECK)*
- CEE*
02-24-21 Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
 - A. *0.0&L* AC
 - B. *0.0&L* AC
 - C. *0.0&L* AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 5. Send the Following to the Crime Laboratory:
 - A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).
 - B. ACA Test Record.

Chemical Test Operator Signature <i>Dona Rustebakke</i>	Date <i>2/24/21</i>
Reviewed By (Crime Laboratory Use Only) <i>Charles E. Edr</i>	Date <i>02-24-2021</i>


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006679
Location = VALL 8164.14.00 09/16
02/24/2021 14:36

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:37
02 Std. Gas	0.081	14:37
03 Room Air	0.000	14:38
04 Std. Gas	0.081	14:38
05 Room Air	0.000	14:39
06 Std. Gas	0.081	14:39
07 Room Air	0.000	14:40

Lot No = 02621080A1
Cyl No = 36
Exp Date = 03/05/2023
County = 02 Oper No. = 100012

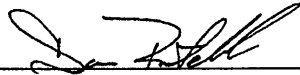

Operator Signature
DANA RUSTEBAKKE

Remarks:

Form 106-I8000

Replaced Gas Cylinder