



ETHANOL GAS STANDARD CYLINDER REPORT
 NORTH DAKOTA OFFICE OF ATTORNEY GENERAL
 CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM
 SFN 59282 (08/2018)

Chemical Test Operator Name (Print) Jeffrey Bauske		
Location Devils Lake		Intoxilyzer® Serial Number 80-006671
Gas Lot Number 02621080A1	Gas Cylinder Number 29	Gas Expiration Date 03/05/2023

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C).
 Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
 A. 0.081 AC
 B. 0.081 AC
 C. 0.082 AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 5. Send the Following to the Crime Laboratory:
 - A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).
 - B. ACA Test Record.

Chemical Test Operator Signature <i>Jeffrey Bauske</i>	Date 04/08/2021
Reviewed By (Crime Laboratory Use Only) <i>Charles E. Edr</i>	Date 5.7.21

SGT. BAUSKE NOTICED A DISCREPANCY WITH HIS ORIGINAL FORM 105-G SUBMISSION PERTAINING TO THE DATE. HE FORWARDED A CORRECTED FORM 105-G TO BE POSTED.

Charles E. Edr
5.7.21



ETHANOL GAS STANDARD CYLINDER REPORT
 NORTH DAKOTA OFFICE OF ATTORNEY GENERAL
 CRIME LABORATORY DIVISION-TOXICOLOGY SECTION/BREATH ALCOHOL PROGRAM
 SFN 59282 (08/2018)

Chemical Test Operator Name (Print) Jeffrey Bauske		
Location Devils Lake	Intoxilyzer® Serial Number 80-006671	
Gas Lot Number 02621080A1	Gas Cylinder Number 29	Gas Expiration Date 03/05/2023

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C).
Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
A. 0.581 AC
B. 0.581 AC
C. 0.582 AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 5. Send the Following to the Crime Laboratory:
 - A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).
 - B. ACA Test Record.

Chemical Test Operator Signature <i>Jeffrey Bauske</i>	Date 4/9/2021
Reviewed By (Crime Laboratory Use Only) <i>[Signature]</i>	Date 4.12.2021

FORM 105-G

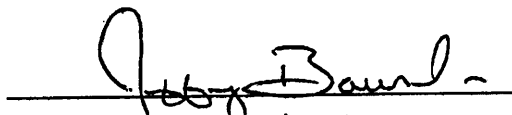
Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006671
Location = DEVL 8164.14.00 09/16
04/08/2021 18:09

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	18:10
02 Std. Gas	0.081	18:10
03 Room Air	0.000	18:11
04 Std. Gas	0.081	18:11
05 Room Air	0.000	18:12
06 Std. Gas	0.082	18:12
07 Room Air	0.000	18:12

Lot No = 02621080A1
Cyl No = 29
Exp Date = 03/05/2023
County = 36 Oper No. = 120501


Operator Signature
JEFFREY BAUSKE

Remarks: AQA TEST

Form 106-I8000