

INTOXILYZER® 8000 INSPECTION

Intoxilyzer® 8000 Serial Number: 80-006685 Location: TOXL

A. General Setup and Checks:

1. Diagnostic Tests Pass and Instrument in Ready Mode
2. Breath tube heated
3. Date, time and location code (Level 2,E). Re-set if necessary.
4. Print test (Level 1,P). Sign and attach test record.
5. Tank monitor (Level 3,D,G). Display and Regulator within 50 psi.
 - a. Display: 675 psi Regulator: 950 psi

Note: Regulator replaced then display read 958 and regulator read 950 psi.

B. Tests (Sign and attach test records):

- Configure simulator for the following tests (Level 1,S).
1. Wet Calibration Check
 - a. Low AC (Level 1,C): Known Value ≤ 0.03 AC: 0.020 AC
Sim. Ser #: MP6036 Lot #: 21340 Exp. Date: 8/24/23
 - b. High AC (Level 1,C): Known Value ≥ 0.25 AC: 0.300 AC
Sim. Ser #: MP3068 Lot #: 207201F Exp. Date: 1/18/24

Configure dry gas standard for the following tests (Level 1,S).

2. Interferent Check (Level 1,B): Known Value: 0.10 AC + 0.05% Acetone
Sim. Ser #: DR3847 Lot #: 1CS6 Exp. Date: NA
 Display reads "Interferent Detect"
3. RFI Check (CMS Mode) Display reads "RFI Detect"
4. Dry Calibration Check (Level 1,C): Known Value 0.080 AC
Gas Cyl Lot #: 34921080A1 Cyl #: 47 Exp. Date: 2/15/24

Test 1 <u>0.082</u>	Test 4 <u>0.081</u>	Test 7 <u>0.081</u>
Test 2 <u>0.080</u>	Test 5 <u>0.081</u>	Test 8 <u>0.081</u>
Test 3 <u>0.081</u>	Test 6 <u>0.082</u>	Test 9 <u>0.082</u>
Average <u>0.081</u>		

C. Remarks/Maintenance: Annual inspection -

Replaced breath tube o-ring and simulator return o-ring.
Replaced regulator after 90 difference between display and regulator was > 50 psi. After regulator was 90 replaced the display + regulator readings were less than 50 psi apart.

Instrument is acceptable to be used in the field.

Janelle Pertschiller
Breath Analyst Signature

6/22/22
Date

[Signature]
Reviewed by

062422
Date

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006685
Location = TOXL 8164.14.00 09/16
06/22/2022 15:51

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890- = |
ABCDEFGHIJKLMNopqrstuvwxyz!@#%\$^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890- = |
ABCDEFGHIJKLMNopqrstuvwxyz!@#%\$^&*()_+?

Current Instrument Setup


Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 26021080A1
Standard Cyl #? 32
Standard Expiration? 10/05/2023
Oper No? 132128

Flow Cal. Date: 09/06/2016
Slope 705
Intercept -687370

IR Calibration Date: 09/06/2016
 3um 9um

	3um	9um
0th Coef(*100):	-14155	-22915
1st Coef(*100):	271436	135296
2nd Coef(*100):	2016	1251
H2O adj(mg/l*10k):	372	505

***** Printer Test End *****

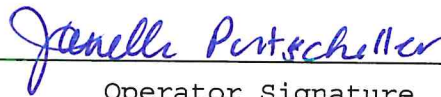

Operator Signature
JANELLE PORTSCHELLER

Remarks: *Print test*

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006685
Location = TOXL 8164.14.00 09/16
06/22/2022 17:10

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	17:11
02 Std. Sol.	0.019	17:13
03 Room Air	0.000	17:13
04 Std. Sol.	0.020	17:14
05 Room Air	0.000	17:14
06 Std. Sol.	0.019	17:15
07 Room Air	0.000	17:16


08 Sim Temp = 34.0°C

Simul Ser No = MP6036

Std Sol No = 21340

County = 08

Oper No. = 132128



Operator Signature
JANELLE PORTSCHELLER

Remarks: *Low AC Check 0.020*

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006685
Location = TOXL 8164.14.00 09/16
06/22/2022 17:19

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	17:20
02 Std. Sol.	0.299	17:21
03 Room Air	0.000	17:22
04 Std. Sol.	0.298	17:23
05 Room Air	0.000	17:23
06 Std. Sol.	0.298	17:24
07 Room Air	0.000	17:25


08 Sim Temp = 34.0°C

Simul Ser No = MP3068

Std Sol No = 202201F

County = 08

Oper No. = 132128

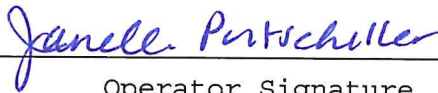

Operator Signature
JANELLE PORTSCHELLER

Remarks: *High AC check 0.300*


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006685
Location = TOXL 8164.14.00 09/16
06/22/2022 17:28

Test	AC	Time
01 Room Air	0.000	17:29
02 *Subject Test	INT*	17:30
03 Room Air	0.000	17:31

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = INTERFERENT CHK
Dr. Lic. = ND/DIS921456
Lot No = 34921080A1
Cyl No = 47
Expiration Date = 02/05/2024
County = 08 Oper No. = 132128

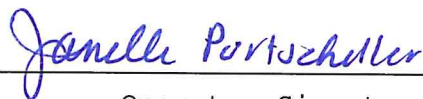

Operator Signature
JANELLE PORTSCHELLER

Remarks: *Interferent check*

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006685
Location = TOXL 8164.14.00 09/16
06/22/2022 17:31

Test	AC	Time
01 Diagnostic	OK	17:32
02 Room Air	0.000	17:32
03 *Subject Test	RFI*	17:32
04 Room Air	0.000	17:33

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = RFI CHK
Dr. Lic. = ND/DIS921456
Lot No = 34921080A1
Cyl No = 47
Expiration Date = 02/05/2024
County = 08 Oper No. = 132128

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature
JANELLE PORTSCHELLER

Remarks: *RFI check*


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006685
Location = TOXL 8164.14.00 09/16
06/22/2022 17:33

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	17:34
02 Std. Gas	0.082	17:34
03 Room Air	0.000	17:35
04 Std. Gas	0.080	17:35
05 Room Air	0.000	17:36
06 Std. Gas	0.081	17:36
07 Room Air	0.000	17:36

Lot No = 34921080A1
Cyl No = 47
Exp Date = 02/05/2024
County = 08 Oper No. = 132128

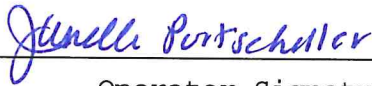

Operator Signature
JANELLE PORTSCHELLER

Remarks:  ^{JP}
Dry calibration check #1-3

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006685
Location = TOXL 8164.14.00 09/16
06/22/2022 17:44

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	17:44
02 Std. Gas	0.081	17:44
03 Room Air	0.000	17:45
04 Std. Gas	0.081	17:45
05 Room Air	0.000	17:46
06 Std. Gas	0.082	17:46
07 Room Air	0.000	17:47


Lot No = 34921080A1

Cyl No = 47

Exp Date = 02/05/2024

County = 08

Oper No. = 132128

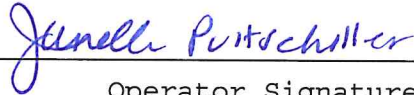

Operator Signature
JANELLE PORTSCHELLER

Remarks: *Dry Calibration check #4-6*


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006685
Location = TOXL 8164.14.00 09/16
06/22/2022 17:48

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	17:49
02 Std. Gas	0.081	17:49
03 Room Air	0.000	17:50
04 Std. Gas	0.081	17:50
05 Room Air	0.000	17:50
06 Std. Gas	0.082	17:51
07 Room Air	0.000	17:51

Lot No = 34921080A1
Cyl No = 47
Exp Date = 02/05/2024
County = 08 Oper No. = 132128

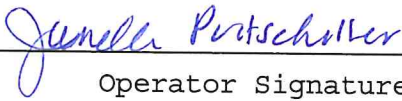

Operator Signature
JANELLE PORTSCHELLER

Remarks: *Dry Calibration Check #7-9*