

INTOXILYZER® 8000 INSPECTION

Intoxilyzer® 8000 Serial Number: 80-006666

Location: TOXL

A. General Setup and Checks:

1. ☒ Diagnostic Tests Pass and Instrument in Ready Mode
2. ☒ Breath tube heated
3. ☒ Date, time and location code (Level 2,E). Re-set if necessary.
4. ☒ Print test (Level 1,P). Sign and attach test record.
5. ☒ Tank monitor (Level 3,D,G). Display and Regulator within 50 psi.
 - a. Display: 944 psi Regulator: 950 psi

B. Tests (Sign and attach test records):

☒ Configure simulator for the following tests (Level 1,S).

1. Wet Calibration Check

- a. ☒ Low AC (Level 1,C): Known Value ≤ 0.03 AC: 0.020 AC
Sim. Ser #: MP3060 Lot #: 8820201A Exp. Date: 1/11/24
- b. ☒ High AC (Level 1,C): Known Value ≥ 0.25 AC: 0.300 AC
Sim. Ser #: MP3064 Lot #: 202012A Exp. Date: 12/8/22

☒ Configure dry gas standard for the following tests (Level 1,S).

2. ☒ Interferent Check (Level 1,B): Known Value: 0.10 AC + 0.05% Acetone
Sim. Ser #: DR5114 Lot #: 1CS6 Exp. Date: NA
☒ Display reads "Interferent Detect"
3. ☒ RFI Check (CMS Mode) ☒ Display reads "RFI Detect"
4. ☒ Dry Calibration Check (Level 1,C): Known Value 0.080 AC
Gas Cyl Lot #: 34421080A1 Cyl #: 49 Exp. Date: 2/5/24

Test 1 <u>0.080</u>	Test 4 <u>0.081</u>	Test 7 <u>0.081</u>
Test 2 <u>0.080</u>	Test 5 <u>0.081</u>	Test 8 <u>0.082</u>
Test 3 <u>0.081</u>	Test 6 <u>0.081</u>	Test 9 <u>0.082</u>
Average <u>0.081</u>		

C. Remarks/Maintenance: Annual inspection.

Replaced breath tube O-ring + simulator return O-ring.

☒ Instrument is acceptable to be used in the field.

Janelle Portscheller
Breath Analyst Signature

06/23/2022
Date

[Signature]
Reviewed by

062422
Date

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006666
Location = TOXL 8164.14.00 09/16
06/23/2022 13:53

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=|
ABCDEFGHIJKLMNOPQRSTUVWXYZ!@#%^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=|
ABCDEFGHIJKLMNOPQRSTUVWXYZ!@#%^&*()_+?

Current Instrument Setup

Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Wet
Standard Value? 0.080
Standard Solution #? 20330
Sim Serial #? MP5318
Oper No? 132128

Flow Cal. Date: 07/07/2021
Slope 636
Intercept -574580

IR Calibration Date: 07/07/2021
 3um 9um

0th Coef(*100):	-14158	-18690
1st Coef(*100):	262999	134281
2nd Coef(*100):	2469	1349
H2O adj(mg/l*10k):	331	452

***** Printer Test End *****



Operator Signature
JANELLE PORTSCHELLER

Remarks: *Print test*

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006666
Location = TOXL 8164.14.00 09/16
06/23/2022 13:58

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:59
02 Std. Sol.	0.019	14:00
03 Room Air	0.000	14:00
04 Std. Sol.	0.019	14:01
05 Room Air	0.000	14:02
06 Std. Sol.	0.019	14:02
07 Room Air	0.000	14:03


08 Sim Temp = 34.0°C

Simul Ser No = MP3060

Std Sol No = 202201A

County = 08

Oper No. = 132128

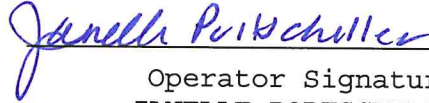

Operator Signature
JANELLE PORTSCHELLER

Remarks:

Low Ac Check 0.020

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006666
Location = TOXL 8164.14.00 09/16
06/23/2022 14:07

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:08
02 Std. Sol.	0.298	14:09
03 Room Air	0.000	14:09
04 Std. Sol.	0.300	14:10
05 Room Air	0.000	14:10
06 Std. Sol.	0.300	14:11
07 Room Air	0.000	14:12

08 Sim Temp = 34.0°C

Simul Ser No = MP3064

Std Sol No = 202012A

County = 08

Oper No. = 132128



Operator Signature
JANELLE PORTSCHELLER

Remarks: *High Hc check 0.300*


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006666
Location = TOXL 8164.14.00 09/16
06/23/2022 14:15

Test	AC	Time
01 Room Air	0.000	14:16
02 *Subject Test	INT*	14:17
03 Room Air	0.000	14:17

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = INTERFERENT CK
Dr. Lic. = ND/DIS821456
Lot No = 34921080A1
Cyl No = 49
Expiration Date = 02/05/2024
County = 08 Oper No. = 132128

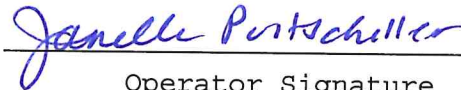

Operator Signature
JANELLE PORTSCHELLER

Remarks: *Interferent Check*

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006666
Location = TOXL 8164.14.00 09/16
06/23/2022 14:18

Test	AC	Time
01 Diagnostic	OK	14:19
02 Room Air	0.000	14:20
03 *Subject Test	RFI*	14:20
04 Room Air	0.000	14:20

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = RFI CHK
Dr. Lic. = ND/DIS821456
Lot No = 34921080A1
Cyl No = 49
Expiration Date = 02/05/2024
County = 08 Oper No. = 132128

I followed the Approved Method and the
instructions displayed by the Intoxilyzer
in conducting this test.

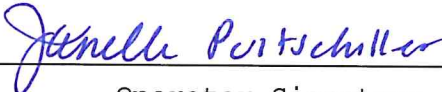

Operator Signature
JANELLE PORTSCHELLER

Remarks: *RFI Check*

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006666
Location = TOXL 8164.14.00 09/16
06/23/2022 14:21

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:21
02 Std. Gas	0.080	14:22
03 Room Air	0.000	14:22
04 Std. Gas	0.080	14:22
05 Room Air	0.000	14:23
06 Std. Gas	0.081	14:23
07 Room Air	0.000	14:24

Lot No = 34921080A1

Cyl No = 49

Exp Date = 02/05/2024

County = 08

Oper No. = 132128


Operator Signature
JANELLE PORTSCHELLER

Remarks: *Dry Calibration Check #1-3*

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006666
Location = TOXL 8164.14.00 09/16
06/23/2022 14:26

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:27
02 Std. Gas	0.081	14:27
03 Room Air	0.000	14:28
04 Std. Gas	0.081	14:28
05 Room Air	0.000	14:29
06 Std. Gas	0.081	14:29
07 Room Air	0.000	14:30


Lot No = 34921080A1

Cyl No = 49

Exp Date = 02/05/2024

County = 08

Oper No. = 132128

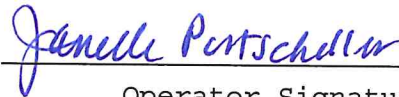

Operator Signature
JANELLE PORTSCHELLER

Remarks: *Dry Calibration Check #4-6*

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-006666
Location = TOXL 8164.14.00 09/16
06/23/2022 16:40

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	16:41
02 Std. Gas	0.081	16:41
03 Room Air	0.000	16:42
04 Std. Gas	0.082	16:42
05 Room Air	0.000	16:42
06 Std. Gas	0.082	16:43
07 Room Air	0.000	16:43


Lot No = 34921080A1

Cyl No = 49

Exp Date = 02/05/2024

County = 08

Oper No. = 132128

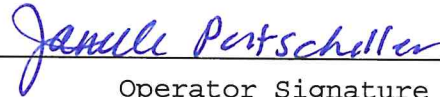

Operator Signature
JANELLE PORTSCHELLER

Remarks: *Dry Calibration Check # 7-9*

Form 106-I8000