



ETHANOL GAS STANDARD CYLINDER REPORT

OFFICE OF ATTORNEY GENERAL
CRIME LABORATORY DIVISION
Toxicology Section/Breath Alcohol Program
SFN 59282 (10-2013)

Chemical Test Operator Name (Print) <i>Conner J Monk</i>		
Location <i>NAPO</i>	Intoxilyzer® Serial Number <i>80-004947</i>	
Gas Lot Number <i>29614080A2</i>	Gas Cylinder Number <i>31</i>	Gas Expiration Date <i>12/5/2016</i>

Check When Done:

1. Scan/Enter Gas Cylinder Information (Level 1, Function S).

2. Perform an ACA Test (Level 1, Function C).

Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).

A. *0.019* AC

B. *0.079* AC

C. *0.080* AC

3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.

N/A 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.

5. Send the Following to the Crime Laboratory:

A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).

B. ACA Test Record.

Chemical Test Operator Signature <i>Conner Monk</i>	Date <i>2/11/2016</i>
Reviewed By (Crime Laboratory Use Only) <i>Deb Shanaver</i>	Date <i>11 Feb 16</i>


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004947
Location = NAPO 8164.13.00 06/09
02/11/2016 12:07

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:07
02 Std. Gas	0.079	12:08
03 Room Air	0.000	12:08
04 Std. Gas	0.079	12:08
05 Room Air	0.000	12:09
06 Std. Gas	0.080	12:09
07 Room Air	0.000	12:10

Lot No = 29614080A2
Cyl No = 31
Exp Date = 12/05/2016
County = 24 Oper No. = 131150

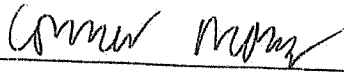

Operator Signature
CONNER MONK

Remarks: