

**INTOXILYZER 5000 G and 8000 INSPECTION**

Serial No.: 80-004943

Location: Tox LAB

Check When Done

- A. Is the warm-up time less than 20 minutes? ✓
- B. Is a three-pronged grounded outlet used? ✓
- C. Is the breath tube heated? ✓
- D. Is the diagnostic check complete? ✓
- E. Is the time, date, and year correct (re-set if necessary)? ✓
- F. Print test. (Attach test record.) ✓
- G. Low AC. Use <0.03 AC in ACA mode. (Attach test record.)  
Sim SN DR5190 Lot No. — Rep AC 0.019 ✓
- H. Linearity Test. Use ≥ 0.25 AC in ACA mode. (Attach test record.)  
Sim SN DR5113 Lot No. — Rep AC 0.285 ✓
- I. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode.  
(Attach test record.)  
Sim SN DR5135 Lot No. — Rep AC 0.100 + 0.05% Acetone ✓
- J. RFI check. Run CMS mode. Key radio on first room air.  
(Attach test record.) ✓
- K. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.  
(Attach test record.)

Lot No. <u>709352</u>	Cyl No. <u>34</u>	Exp. Date <u>23 Dec 2011</u>
Test 1 <u>0.081</u>	Test 4 <u>0.081</u>	Test 7 <u>0.081</u>
Test 2 <u>0.082</u>	Test 5 <u>0.081</u>	Test 8 <u>0.081</u>
Test 3 <u>0.081</u>	Test 6 <u>0.081</u>	Test 9 <u>0.080</u>
Average <u>0.081</u>		

L. Remarks/Maintenance record: Annual Inspection

<u>10 Mar 10</u> Date	<u>Mary Pearson</u> Field Inspector's Signature
<u>15 Mar '10</u> Date	<u>Deb Kashu</u> Reviewed by

80-004943  
23/31 PM

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004943  
Location = TOXL      8164.13.00 06/09  
03/10/2010      09:05

\*\*\*\*\* Printer Test \*\*\*\*\*

abcdefghijklmnopqrstuvwxyz1234567890-=\_|  
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#%&^&\*()\_+?

abcdefghijklmnopqrstuvwxyz1234567890-=\_|  
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#%&^&\*()\_+?

Current Instrument Setup

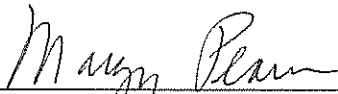
Data Entry Mode:      Enabled  
Start Test Sequence:      DABACABA  
Display Prelim Rslt?      Yes  
Display Third Digit?      Yes  
Inhib Printer(Y/N)?      No  
Display Volume?      No  
Disable On Memfull?      No  
# of Print Copies?      1  
Select Std (D/W/I)?      Dry  
Standard Value?      0.080  
Standard Lot #?      709352  
Standard Cyl #?      34  
Standard Expiration?      12/23/2011  
Oper No?      107501

Flow Cal. Date:      01/18/2010  
Slope      693  
Intercept      -665387

IR Calibration Date:      01/18/2010

	3um	9um
0th Coef(*100):	-20351	-24446
1st Coef(*100):	248542	131200
2nd Coef(*100):	2180	1190
H2O adj(mg/l*10k):	594	450

\*\*\*\*\* Printer Test End \*\*\*\*\*



Operator Signature  
MARGY PEARSON

Remarks:

*Print Test*

Form 106-18000

80-004943  
24/31 FKL

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004943  
Location = TOXL      8164.13.00 06/09  
03/10/2010      09:39

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:40
02 Std. Sol.	0.016	09:41
03 Room Air	0.000	09:41
04 Std. Sol.	0.017	09:42
05 Room Air	0.000	09:43
06 Std. Sol.	0.018	09:43
07 Room Air	0.000	09:44

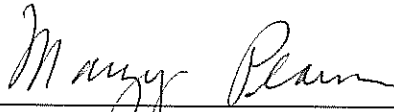
08 Sim Temp = 34.0°C

Simul Ser No = 000000

Std Sol No = 0000

County = 08

Oper No. = 107501



Operator Signature  
MARGY PEARSON

Remarks:

Low AC 0.019AC

Form 106-I8000

80-004943  
25/31 FM

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004943  
Location = TOXL      8164.13.00 06/09  
03/10/2010      09:25

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:26
02 Std. Sol.	0.288	09:27
03 Room Air	0.000	09:27
04 Std. Sol.	0.289	09:28
05 Room Air	0.000	09:28
06 Std. Sol.	0.290	09:29
07 Room Air	0.000	09:30

08 Sim Temp = 34.0°C

Simul Ser No = 000000

Std Sol No = 0000

County = 08

Oper No. = 107501



Operator Signature  
MARGY PEARSON

Remarks:

*Linearity*

Form 106-I8000

80-004943  
26/31/10


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004943  
Location = TOXL      8164.13.00 06/09  
03/10/2010      09:33

Test	AC	Time
01 Room Air	0.000	09:34
02 *Subject Test	INT*	09:35
03 Room Air	0.000	09:35

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female      Weight = 150  
Test = DUI      Cit = NA  
Dr. Lic. = ND/DIS821456  
Simul Ser No = 000000  
Std Sol No = 0000  
County = 08      Oper No. = 107501

  
\_\_\_\_\_  
Operator Signature  
MARGY PEARSON

Remarks:

*Inter - AIC + Acetone*

Form 106-18000

80-004943  
2/31/14

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004943  
Location = TOXL      8164.13.00 06/09  
03/10/2010      09:31

Test	AC	Time
01 Room Air	RFI*	09:31
02 Room Air	0.000	09:32

\*Invalid Test  
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female      Weight = 150  
Test = DUI      Cit = NA  
Dr. Lic. = ND/DIS821456  
Simul Ser No = 000000  
Std Sol No = 0000  
County = 08      Oper No. = 107501



Operator Signature  
MARGY PEARSON

Remarks:

RFI

Form 106-I8000

80-004943  
28/31 HZ

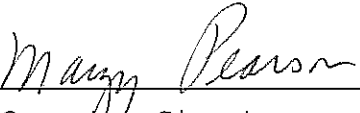
Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004943  
Location = TOXL      8164.13.00 06/09  
03/10/2010      10:32

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:32
02 Std. Gas	0.081	10:33
03 Room Air	0.000	10:33
04 Std. Gas	0.082	10:33
05 Room Air	0.000	10:34
06 Std. Gas	0.081	10:34
07 Room Air	0.000	10:35

Lot No = 709352  
Cyl No = 34  
Exp Date = 12/23/2011  
County = 08      Oper No. = 107501

  
\_\_\_\_\_  
Operator Signature  
MARGY PEARSON

Remarks:

*ACA Cal Check*

Form 106-I8000

80-004943  
39/31 KHL

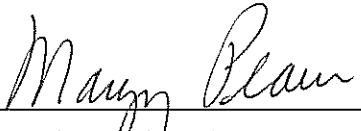
Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004943  
Location = TOXL      8164.13.00 06/09  
03/10/2010      10:41

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:42
02 Std. Gas	0.081	10:42
03 Room Air	0.000	10:42
04 Std. Gas	0.081	10:43
05 Room Air	0.000	10:43
06 Std. Gas	0.081	10:44
07 Room Air	0.000	10:44

Lot No = 709352  
Cyl No = 34  
Exp Date = 12/23/2011  
County = 08      Oper No. = 107501

  
\_\_\_\_\_  
Operator Signature  
MARGY PEARSON

Remarks:

*ACA Cal Check*

Form 106-I8000

*80-004943  
30/31/10*



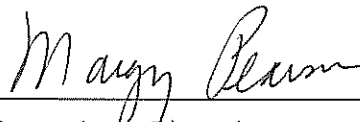
Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004943  
Location = TOXL      8164.13.00 06/09  
03/10/2010      10:45

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:46
02 Std. Gas	0.081	10:46
03 Room Air	0.000	10:47
04 Std. Gas	0.081	10:47
05 Room Air	0.000	10:48
06 Std. Gas	0.080	10:48
07 Room Air	0.000	10:48

Lot No = 709352  
Cyl No = 34  
Exp Date = 12/23/2011  
County = 08      Oper No. = 107501

  
\_\_\_\_\_  
Operator Signature  
MARGY PEARSON

Remarks:

*ACA Cal Check*

Form 106-I8000

80-004943  
3/31 JFA