



INTOXILYZER® 8000 INSTALLATION AND REPAIR CHECKOUT
 OFFICE OF ATTORNEY GENERAL
 CRIME LABORATORY DIVISION
 Toxicology Section/Breath Alcohol Program
 SFN 59281 (5-2011)

Serial Number 80-005939	Instrument Location Cass County Jail
Reason for Install/Repair <input checked="" type="checkbox"/> Install After Receiving From Crime Laboratory <input type="checkbox"/> Install After Location Change <input type="checkbox"/> Other (Specify) _____	

Check When Done:

- 1. Surge Protector Installed/Properly Grounded.
- 2. Telephone Line Connected to Intoxilyzer® 8000.
- 3. Breath Tube Heated.
- 4. Enter Preliminary Data (i.e. Date, Time, DST (Y), and Location; Level 2, Function E).
- 5. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 6. Run Tests:
 - A. Print Test (Level 1, Function P). ✓
 - B. ACA Test (Level 1, Function C). ✓
 - C. RFI Test (Use CMS Mode; Key Radio During First Room Air). ✓
- 7. Repair and/or Maintenance Performed (if Any): N/A
- 8. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 9. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 10. Send the Following to the Crime Laboratory: Completed Intoxilyzer®8000 Installation and Repair Checkout (SFN59281, Form 104-G), Print Test, ACA Test, and RFI Test.

Chad Thompson

Field Inspector Signature <i>[Signature]</i>	Date 3/13/14 ✓
Reviewed By (Crime Laboratory Use Only) <i>Deb Kashner</i>	Date <i>DK</i> 26 Mar 2014

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005939
Location = CASS 8164.13.00 06/09
03/13/2014 19:41

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	19:42
02 Std. Gas	0.081	19:42
03 Room Air	0.000	19:43
04 Std. Gas	0.081	19:43
05 Room Air	0.000	19:43
06 Std. Gas	0.081	19:44
07 Room Air	0.000	19:44

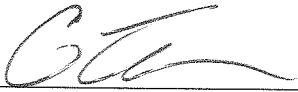
Lot No = 16512080A2

Cyl No = 66

Exp Date = 08/01/2014

County = 09

Oper No. = 050910

 3524

Operator Signature
CHAD THOMPSON

Remarks:

ACA Test

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005939
Location = CASS 8164.13.00 06/09
03/13/2014 19:56

Test	AC	Time
01 Diagnostic	OK	19:57
02 Room Air	RFI*	19:57
03 Room Air	0.000	19:58

*Invalid Test
Inhibited - RFI

Sub Name = RFI TEST, NA NA
Sub DOB = 01/01/1990
Sub Sex = Male Weight = 200
Test = DUI Cit = NA
Dr. Lic. = ND/NA
Lot No = 16512080A2
Cyl No = 66
Expiration Date = 08/01/2014
County = 09 Oper No. = 050910

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

 3524

Operator Signature
CHAD THOMPSON

Remarks:

RFI Test