



ETHANOL BREATH STANDARD CYLINDER REPORT

OFFICE OF ATTORNEY GENERAL

CRIME LABORATORY DIVISION

Toxicology Section/Breath Alcohol Program

SFN 59282 (5-2011)

Chemical Test Operator Name (Print) <i>Matthew W. Anderson</i>		
Location <i>Richland County Jail</i>	Intoxilyzer® Serial Number <i>80-05944</i>	
Gas Lot Number <i>16512080 A2</i>	Gas Cylinder Number <i>58</i>	Gas Expiration Date <i>08/01/2014</i>

Check When Done:

1. Scan/Enter Gas Cylinder Information (Level 1, Function S).

2. Perform an ACA Test (Level 1, Function C).

Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).

A. 0.082 AC

B. 0.081 AC

C. 0.082 AC

3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.

4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.

5. Send the Following to the Crime Laboratory:

A. Ethanol Breath Standard Cylinder Report (SFN59282, Form 105-G).

B. ACA Test Record.

Chemical Test Operator Signature <i>Matthew W. Anderson</i>	Date <i>5-26-2014</i>
Reviewed By (Crime Laboratory Use Only) <i>Deb Kashner</i>	Date <i>28 May 14</i>


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005944
Location = WAHP 8164.13.00 06/09
03/26/2014 21:09

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	21:10
02 Std. Gas	0.082	21:10
03 Room Air	0.000	21:11
04 Std. Gas	0.081	21:11
05 Room Air	0.000	21:12
06 Std. Gas	0.082	21:12
07 Room Air	0.000	21:12

Lot No = 16512080A2
Cyl No = 58
Exp Date = 08/01/2014
County = 39 Oper No. = 020401

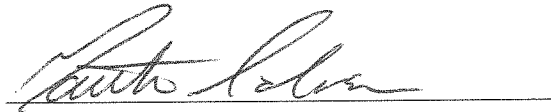

Operator Signature
MATTHEW ANDERSON

Remarks:

Ethanol Brock Standard cylinder report

Form 106-I8000