



ETHANOL BREATH STANDARD CYLINDER REPORT

OFFICE OF ATTORNEY GENERAL

CRIME LABORATORY DIVISION

Toxicology Section/Breath Alcohol Program

SFN 59282 (5-2011)

Chemical Test Operator Name (Print) <i>Johnny Z Lawson</i>		
Location <i>Bowman Co SO</i>	Intoxilyzer® Serial Number <i>80-004191</i>	
Gas Lot Number <i>05911080A1</i>	Gas Cylinder Number <i>4</i>	Gas Expiration Date <i>02-01-13</i>

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C).
Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
 - A. 0. 076 AC
 - B. 0. 077 AC
 - C. 0. 077 AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 5. Send the Following to the Crime Laboratory:
 - A. Ethanol Breath Standard Cylinder Report (SFN59282, Form 105-G).
 - B. ACA Test Record.

Chemical Test Operator Signature <i>Johnny Z Lawson</i>	Date <i>09-15-11</i>
Reviewed By (Crime Laboratory Use Only) <i>Deb Kashner</i>	Date <i>02 Apr 2012</i>

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = BOWM 8164.13.00 06/09
09/15/2011 13:14

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:14
02 Std. Gas	0.076	13:15
03 Room Air	0.000	13:15
04 Std. Gas	0.077	13:16
05 Room Air	0.000	13:16
06 Std. Gas	0.077	13:17
07 Room Air	0.000	13:17


Lot No = 05911080A1

Cyl No = 4

Exp Date = 02/01/2013

County = 06

Oper No. = 020814

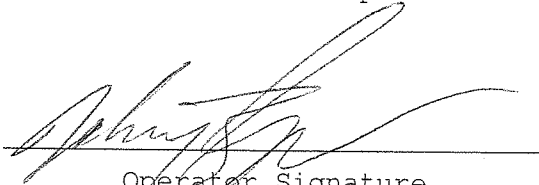

Operator Signature
JOHNNY Z. LAWSON

Remarks: