

INTOXILYZER 8000 INSPECTION

Serial No.: 80-004191

Location: TOXL

- | | Check When Done |
|--|-----------------|
| A. Is the warm-up time less than 20 minutes? | <u>✓</u> |
| B. Is a three-pronged grounded outlet used? | <u>✓</u> |
| C. Is the breath tube heated? | <u>✓</u> |
| D. Is the diagnostic check complete? | <u>✓</u> |
| E. Is the time, date, and year correct (re-set if necessary)? | <u>✓</u> |
| F. Print test. (Attach test record.) | <u>✓</u> |
| G. Low AC. Use <0.03 AC in ACA mode. (Attach test record.)
Sim SN <u>DR5144</u> Lot No. <u>—</u> Rep AC <u>0.021</u> | <u>✓</u> |
| H. Linearity Test. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN <u>DR3847</u> Lot No. <u>—</u> Rep AC <u>0.294</u> | <u>✓</u> |
| I. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode.
(Attach test record.)
Sim SN <u>D83452</u> Lot No. <u>—</u> Rep AC <u>0.10</u> | <u>✓</u> |
| J. RFI check. Run CMS mode. Key radio on first room air.
(Attach test record.) | <u>✓</u> |
| K. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) | |

Lot No. <u>15410080A1</u>	Cyl No. <u>12</u>	Exp. Date <u>7-1-12</u>
Test 1 <u>0.079</u>	Test 4 <u>0.080</u>	Test 7 <u>0.079</u>
Test 2 <u>0.079</u>	Test 5 <u>0.080</u>	Test 8 <u>0.080</u>
Test 3 <u>0.079</u>	Test 6 <u>0.080</u>	Test 9 <u>0.079</u>
Average <u>0.079</u>		

L. Tank pressure, Level 3 Func D Sub G P: 375 psi Reg 380 psi

M. Remarks/Maintenance record: Annual Inspection

4-5-11
Date

[Signature]
Field Inspector's Signature

05 Apr 2011
Date

Deb Kashner
Reviewed by

80-004191
P. 77 of 85
FM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = TOXL 8164.13.00 06/09
04/05/2011 11:25

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

Current Instrument Setup


Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? ~~05911080A1~~¹² 15410080A1
Standard Cyl #? ~~17~~¹²
Standard Expiration? ~~02/01/2013~~¹² 7-1-2012
Oper No? 555555

Flow Cal. Date: 01/28/2010
Slope 738
Intercept -777574

IR Calibration Date: 01/28/2010

	3um	9um
0th Coef(*100):	-23128	-33609
1st Coef(*100):	272485	131724
2nd Coef(*100):	1490	1056
H2O adj(mg/l*10k):	432	614

***** Printer Test End *****

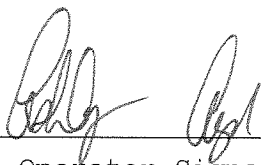

Operator Signature
ASHLEY AYD

Remarks: *Print test*

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = TOXL 8164.13.00 06/09
04/05/2011 11:26

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:27
02 Std. Sol.	0.018	11:27
03 Room Air	0.000	11:28
04 Std. Sol.	0.018	11:29
05 Room Air	0.000	11:29
06 Std. Sol.	0.019	11:30
07 Room Air	0.000	11:31


08 Sim Temp = 34.0°C

Simul Ser No = 0

Std Sol No = 0

County = 08

Oper No. = 555555

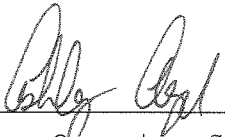

Operator Signature
ASHLEY AYD

Remarks: Low AC 0.021 AC

Form 106-I8000

80-004191
P. 79 of 85
JAN

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = TOXL 8164.13.00 06/09
04/05/2011 11:32

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:32
02 Std. Sol.	0.294	11:33
03 Room Air	0.000	11:34
04 Std. Sol.	0.293	11:35
05 Room Air	0.000	11:35
06 Std. Sol.	0.295	11:36
07 Room Air	0.000	11:37


08 Sim Temp = 34.0°C

Simul Ser No = 0

Std Sol No = 0

County = 08

Oper No. = 555555

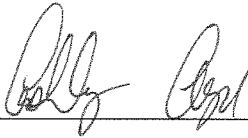

Operator Signature
ASHLEY AYD

Remarks: High AC ~~0.0~~⁰⁰ 0.294 AC

Form 106-I8000

80-004191
P. 80 of 85
AM


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = TOXL 8164.13.00 06/09
04/05/2011 11:42

Test	AC	Time
01 Room Air	0.000	11:43
02 *Subject Test	INT*	11:43
03 Room Air	0.000	11:44

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Lot No = 15410080A1
Cyl No = 12
Expiration Date = 07/01/2012
County = 08 Oper No. = 555555

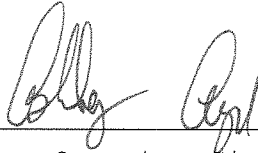

Operator Signature
ASHLEY AYD

Remarks: *Interferent check 0.10 AC + Acetone*

Form 106-I8000

80-004191
P. 81-485
JAM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = TOXL 8164.13.00 06/09
04/05/2011 11:45

Test	AC	Time
01 Diagnostic	OK	11:46
02 Room Air	RFI*	11:46
03 Room Air	0.000	11:47

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Lot No = 15410080A1
Cyl No = 12
Expiration Date = 07/01/2012
County = 08 Oper No. = 555555

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature
ASHLEY AYD

Remarks: RFI

Form 106-I8000

80-004191
P. 82 of 85
KTM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = TOXL 8164.13.00 06/09
04/05/2011 11:48

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:48
02 Std. Gas	0.079	11:49
03 Room Air	0.000	11:49
04 Std. Gas	0.079	11:50
05 Room Air	0.000	11:50
06 Std. Gas	0.079	11:51
07 Room Air	0.000	11:51

Lot No = 15410080A1
Cyl No = 12
Exp Date = 07/01/2012
County = 08 Oper No. = 555555



Operator Signature
ASHLEY AYD

Remarks: Calibration checks 0.080 AC

Form 106-I8000

80-004191
P.83485
FM


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = TOXL 8164.13.00 06/09
04/05/2011 11:52

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:53
02 Std. Gas	0.080	11:53
03 Room Air	0.000	11:53
04 Std. Gas	0.080	11:54
05 Room Air	0.000	11:54
06 Std. Gas	0.080	11:55
07 Room Air	0.000	11:55

Lot No = 15410080A1
Cyl No = 12
Exp Date = 07/01/2012
County = 08 Oper No. = 555555

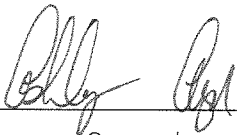

Operator Signature
ASHLEY AYD

Remarks: Calibration check 0.080 AC

Form 106-I8000

80-004191
P. 84 of 85
KM


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004191
Location = TOXL 8164.13.00 06/09
04/05/2011 11:56

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:57
02 Std. Gas	0.079	11:57
03 Room Air	0.000	11:58
04 Std. Gas	0.080	11:58
05 Room Air	0.000	11:58
06 Std. Gas	0.079	11:59
07 Room Air	0.000	11:59

Lot No = 15410080A1
Cyl No = 12
Exp Date = 07/01/2012
County = 08 Oper No. = 555555

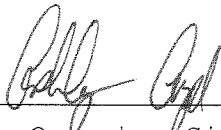

Operator Signature
ASHLEY AYD

Remarks: Calibration check 0.080 AC

Form 106-I8000

80-004191
P. 85 of 85
KTM