

**INTOXILYZER 8000 INSPECTION**

Serial No.: 80-004184

Location: TOXL

Check When Done

- A. Is the warm-up time less than 20 minutes? ✓
- B. Is a three-pronged grounded outlet used? ✓
- C. Is the breath tube heated? ✓
- D. Is the diagnostic check complete? ✓
- E. Is the time, date, and year correct (re-set if necessary)? ✓
- F. Print test. (Attach test record.) ✓
- G. Low AC. Use <0.03 AC in ACA mode. (Attach test record.)  
Sim SN DR5144 Lot No. — Rep AC 0.021 ✓
- H. Linearity Test. Use ≥ 0.25 AC in ACA mode. (Attach test record.)  
Sim SN DR3847 Lot No. — Rep AC 0.294 ✓
- I. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode.  
(Attach test record.)  
Sim SN DR3452 Lot No. — Rep AC 0.10 AC + Acetone ✓
- J. RFI check. Run CMS mode. Key radio on first room air.  
(Attach test record.) ✓
- K. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.  
(Attach test record.)

Lot No. <u>15410080A1</u>	Cyl No. <u>12</u>	Exp. Date <u>7-1-12</u>
Test 1 <u>0.081</u>	Test 4 <u>0.081</u>	Test 7 <u>0.081</u>
Test 2 <u>0.081</u>	Test 5 <u>0.081</u>	Test 8 <u>0.080</u>
Test 3 <u>0.080</u>	Test 6 <u>0.080</u>	Test 9 <u>0.080</u>
Average <u>0.081</u>		

L. Tank pressure, Level 3 Func D Sub G P: 202 psi Reg 220 psi

M. Remarks/Maintenance record: Annual Inspection

4-7-11  
Date

[Signature]  
Field Inspector's Signature

07 Apr 2011  
Date

Deb Kashner  
Reviewed by

80-004184  
P. 109 of 117  
FM

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-004184  
Location = TOXL                    8164.13.00 06/09  
04/07/2011                            11:12

\*\*\*\*\* Printer Test \*\*\*\*\*

abcdefghijklmnopqrstuvwxyz1234567890-=\_|  
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&\*()\_+?

abcdefghijklmnopqrstuvwxyz1234567890-=\_|  
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&\*()\_+?

Current Instrument Setup


Data Entry Mode:            Enabled  
Start Test Sequence:        DABACABA  
Display Prelim Rslt?        Yes  
Display Third Digit?        Yes  
Inhib Printer(Y/N)?        No  
Display Volume?            No  
Disable On Memfull?        Yes  
# of Print Copies?         1  
Select Std (D/W/I)?        Dry  
Standard Value?            0.080  
Standard Lot #?            15410080A1  
Standard Cyl #?            12  
Standard Expiration?       07/01/2012  
Oper No?                    555555

Flow Cal. Date:            07/30/2010  
Slope                        804  
Intercept                   -925008

IR Calibration Date:        07/30/2010  
                                    3um                            9um

-----  
0th Coef(\*100):            -19612                        -18975  
1st Coef(\*100):            270159                        137378  
2nd Coef(\*100):            1174                            963  
H2O adj(mg/l\*10k):        427                            372

\*\*\*\*\* Printer Test End \*\*\*\*\*

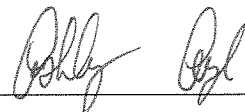

Operator Signature  
ASHLEY AYD

Remarks: *Print test*

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004184  
Location = TOXL      8164.13.00 06/09  
04/07/2011      11:16

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:16
02 Std. Sol.	0.021	11:17
03 Room Air	0.000	11:18
04 Std. Sol.	0.020	11:19
05 Room Air	0.000	11:19
06 Std. Sol.	0.021	11:20
07 Room Air	0.000	11:20


08 Sim Temp = 34.0°C

Simul Ser No = 0

Std Sol No = 0

County = 08

Oper No. = 555555

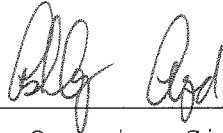

Operator Signature  
ASHLEY AYD

Remarks:

Low AC 0.021 AC

Form 106-I8000

80-004184  
p. 111 of 117  
YFM

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004184  
Location = TOXL      8164.13.00 06/09  
04/07/2011      11:21

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:22
02 Std. Sol.	0.297	11:22
03 Room Air	0.000	11:23
04 Std. Sol.	0.299	11:24
05 Room Air	0.000	11:24
06 Std. Sol.	0.299	11:25
07 Room Air	0.000	11:26


08 Sim Temp = 34.0°C

Simul Ser No = 0

Std Sol No = 0

County = 08

Oper No. = 555555

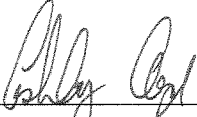  
\_\_\_\_\_  
Operator Signature  
ASHLEY AYD

Remarks: *Linearity test 0.294 AC*

Form 106-I8000

*80-004184  
P. 112 of 117  
AM*


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-004184  
Location = TOXL                8164.13.00 06/09  
04/07/2011                        11:29

Test	AC	Time
01 Room Air	0.000	11:30
02 *Subject Test	INT*	11:30
03 Room Air	0.000	11:31

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = NA  
Dr. Lic. = ND/DIS821456  
Lot No = 15410080A1  
Cyl No = 12  
Expiration Date = 07/01/2012  
County = 08                          Oper No. = 555555

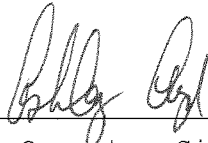

Operator Signature  
ASHLEY AYD

Remarks: *Interferent 0.10 AC+Acetone*

Form 106-I8000

80-004184  
P.113d/117  
ym

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-004184  
Location = TOXL            8164.13.00 06/09  
04/07/2011            11:31

Test	AC	Time
01 Diagnostic	OK	11:32
02 Room Air	RFI*	11:32
03 Room Air	0.000	11:33

\*Invalid Test  
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female            Weight = 150  
Test = DUI            Cit = NA  
Dr. Lic. = ND/DIS821456  
Lot No = 15410080A1  
Cyl No = 12  
Expiration Date = 07/01/2012  
County = 08            Oper No. = 555555

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

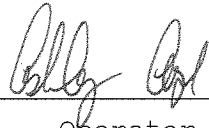

Operator Signature  
ASHLEY AYD

Remarks: RFI

Form 106-I8000

80-009184  
P.114 of 117  
VAM


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004184  
Location = TOXL      8164.13.00 06/09  
04/07/2011      11:33

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:34
02 Std. Gas	0.081	11:34
03 Room Air	0.000	11:35
04 Std. Gas	0.081	11:35
05 Room Air	0.000	11:36
06 Std. Gas	0.080	11:36
07 Room Air	0.000	11:37

Lot No = 15410080A1  
Cyl No = 12  
Exp Date = 07/01/2012  
County = 08      Oper No. = 555555

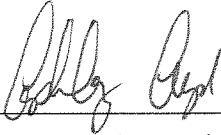  
\_\_\_\_\_  
Operator Signature  
ASHLEY AYD

Remarks: Calibration check 0.080 AC

Form 106-I8000

80-004184  
P. 1154/117  
AM


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004184  
Location = TOXL      8164.13.00 06/09  
04/07/2011      11:37

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:38
02 Std. Gas	0.081	11:38
03 Room Air	0.000	11:39
04 Std. Gas	0.081	11:39
05 Room Air	0.000	11:40
06 Std. Gas	0.080	11:40
07 Room Air	0.000	11:40

Lot No = 15410080A1  
Cyl No = 12  
Exp Date = 07/01/2012  
County = 08      Oper No. = 555555

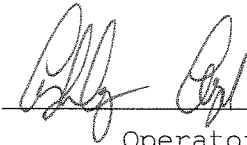  
\_\_\_\_\_  
Operator Signature  
ASHLEY AYD

Remarks: Calibration check 0.080 AC

Form 106-I8000

80-004184  
P. 116 of 117  
EM




Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004184  
Location = TOXL      8164.13.00 06/09  
04/07/2011      12:06

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:06
02 Std. Gas	0.081	12:07
03 Room Air	0.000	12:07
04 Std. Gas	0.080	12:08
05 Room Air	0.000	12:08
06 Std. Gas	0.080	12:09
07 Room Air	0.000	12:09

Lot No = 15410080A1  
Cyl No = 12  
Exp Date = 07/01/2012  
County = 08      Oper No. = 555555

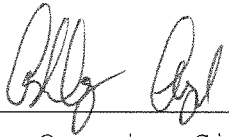  
\_\_\_\_\_  
Operator Signature  
ASHLEY AYD

Remarks: Calibration check 0.080 AC

Form 106-I8000

80-004184  
p. 117 of 117  
KLN