



INTOXILYZER® 8000 INSTALLATION AND REPAIR CHECKOUT
 OFFICE OF ATTORNEY GENERAL
 CRIME LABORATORY DIVISION
 Toxicology Section/Breath Alcohol Program
 SFN 59281 (1-2010)

Serial Number: <u>80-004198</u>	Instrument Location: <u>Stutsman Co.</u>
Reason for Install/Repair:	
<input checked="" type="checkbox"/> Install After Lab Repair/Inspection <input type="checkbox"/> Other (Specify) _____	

Check When Done:

- A. Surge Protector Installed/Properly Grounded.
- B. Telephone Line Connected to Intoxilyzer® 8000.
- C. Breath Tube Heated.
- D. Review/Enter Preliminary Data Entry (Level 2, Function E).
- E. Review/Enter Gas Setup (Level 1, Function S).
- F. Print Test (Level 1, Function P). (Attach Test Record.)
- G. ACA Test (Level 1, Function C). (Attach Test Record.)
- H. RFI Test (Use CMS Mode-First Room Air). (Attach Test Record.)
- I. Repair and/or Maintenance Performed: _____

[Signature]
 Field Inspector Signature

Deb Kashner
 Reviewed By

04-08-11
 Date

19 Apr 11
 Date

80-004198
 P. 133 of 137
 VM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004198
Location = ~~TOX~~ JAME 8164.13.00 06/09
04/08/2011 19 April DK 18:11

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

Current Instrument Setup

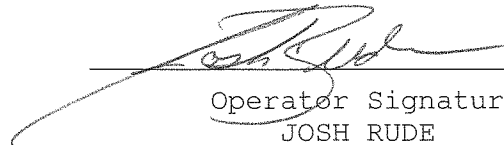
Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 05911080A1
Standard Cyl #? 21
Standard Expiration? 02/01/2013
Oper No? 100414

Flow Cal. Date: 02/02/2011
Slope 658
Intercept -563002

IR Calibration Date: 02/03/2011

	3um	9um
0th Coef(*100):	-24997	-30437
1st Coef(*100):	271010	135136
2nd Coef(*100):	2084	890
H2O adj(mg/l*10k):	431	443

***** Printer Test End *****


Operator Signature
JOSH RUDE

Remarks:

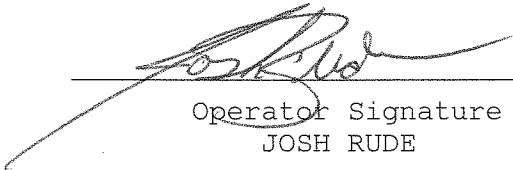
Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004198
Location = ~~TOXIL JAME~~ 8164.13.00 06/09
04/08/2011 19 April DK 18:33

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	18:34
02 Std. Gas	0.077	18:34
03 Room Air	0.000	18:35
04 Std. Gas	0.076	18:35
05 Room Air	0.000	18:35
06 Std. Gas	0.076	18:36
07 Room Air	0.000	18:36

Lot No = 05911080A1
Cyl No = 21
Exp Date = 02/01/2013
County = 47 Oper No. = 100414


Operator Signature
JOSH RUDE

Remarks:

Form 106-I8000

80-004198
P.135 of 137
JRM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004198
Location = ~~FOXL~~ JAME 8164.13.00 06/09
04/08/2011 19 Apr 11 18:38
 DK

DRY CAL CHECK

Test	AC	Time
01 Room Air	RFI*	18:38
02 Room Air	0.000	18:39

*Invalid Test
 Inhibited - RFI

Lot No = 05911080A1
Cyl No = 21
Exp Date = 02/01/2013
County = 47 Oper No. = 100414


Operator Signature
JOSH RUDE

Remarks:

Form 106-I8000

80-004198
P. 136 of 137
FR

19 April 2011

I called Josh Rude informing him that he forgot to change the Location Code when he installed Intoxilyzer 80-004198. I told him that I can change it on the paperwork he sent in, but it needs to be changed *in* the 8000. He told me that he will go to the jail today and enter the correct Location Code.

Deb Kasher

80-004198
P. 137 of 137
FTM