

INTOXILYZER 8000 INSPECTION

Serial No.: 80-004198 Location: TOXL

- | | Check When Done |
|--|-----------------|
| A. Is the warm-up time less than 20 minutes? | <u>✓</u> |
| B. Is a three-pronged grounded outlet used? | <u>✓</u> |
| C. Is the breath tube heated? | <u>✓</u> |
| D. Is the diagnostic check complete? | <u>✓</u> |
| E. Is the time, date, and year correct (re-set if necessary)? | <u>✓</u> |
| F. Print test. (Attach test record.) | <u>✓</u> |
| G. Low AC. Use <0.03 AC in ACA mode. (Attach test record.)
Sim SN <u>DR5132</u> Lot No. <u>10290</u> Rep AC <u>0.020</u> | <u>✓</u> |
| H. Linearity Test. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN <u>DR3803</u> Lot No. <u>11110</u> Rep AC <u>0.300</u> | <u>✓</u> |
| I. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode.
(Attach test record.)
Sim SN <u>DR3846</u> Lot No. <u> </u> Rep AC <u>0.100</u> | <u>✓</u> |
| J. RFI check. Run CMS mode. Key radio on first room air.
(Attach test record.) | <u>✓</u> |
| K. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) | |

Lot No. <u>05911080A1</u>	Cyl No. <u>14</u>	Exp. Date <u>02/01/2013</u>
Test 1 <u>0.080</u>	Test 4 <u>0.079</u>	Test 7 <u>0.080</u>
Test 2 <u>0.080</u>	Test 5 <u>0.079</u>	Test 8 <u>0.080</u>
Test 3 <u>0.080</u>	Test 6 <u>0.080</u>	Test 9 <u>0.080</u>
Average <u>0.080</u>		

L. Tank pressure, Level 3 Func D Sub G P: 617 psi Reg 600 psi

M. Remarks/Maintenance record: Annual inspection

05 Apr 12
Date

09 APR 12
Date

Deb Kashur
Field Inspector's Signature

Charles E. [Signature]
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004198
Location = TOXL 8164.13.00 06/09
04/05/2012 14:50

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890- = |
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&^*()_+?

abcdefghijklmnopqrstuvwxyz1234567890- = |
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&^*()_+?

Current Instrument Setup

Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 05911080A1
Standard Cyl #? 14
Standard Expiration? 02/01/2013
Oper No? 777777

Flow Cal. Date: 01/09/2012
Slope 637
Intercept -467208

IR Calibration Date: 01/09/2012


	3um	9um
0th Coef(*100):	-31777	-32014
1st Coef(*100):	271021	132277
2nd Coef(*100):	1454	964
H2O adj(mg/l*10k):	558	576

***** Printer Test End *****



Operator Signature
DEB KASHUR

Remarks:

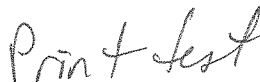

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004198
Location = TOXL 8164.13.00 06/09
04/05/2012 14:51

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:51
02 Std. Sol.	0.019	14:52
03 Room Air	0.000	14:53
04 Std. Sol.	0.019	14:53
05 Room Air	0.000	14:54
06 Std. Sol.	0.019	14:55
07 Room Air	0.000	14:55

08 Sim Temp = 34.0°C

Simul Ser No = 0

Std Sol No = 0

County = 08

Oper No. = 777777



Operator Signature
DEB KASHUR

Remarks:

Low AC 0.020 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004198
Location = TOXL 8164.13.00 06/09
04/05/2012 15:00

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:01
02 Std. Sol.	0.298	15:02
03 Room Air	0.000	15:02
04 Std. Sol.	0.296	15:03
05 Room Air	0.000	15:03
06 Std. Sol.	0.296	15:04
07 Room Air	0.000	15:05

08 Sim Temp = 34.0°C

Simul Ser No = 0

Std Sol No = 0

County = 08

Oper No. = 77777

Deb Kashur

Operator Signature
DEB KASHUR

Remarks:

Linearity test 0.300AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004198
Location = TOXL 8164.13.00 06/09
04/05/2012 16:30

Test	AC	Time
01 Room Air	0.000	16:30
02 *Subject Test	INT*	16:31
03 Room Air	0.000	16:31

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS921456
Lot No = 05911080A1
Cyl No = 14
Expiration Date = 02/01/2013
County = 08 Oper No. = 777777

Deb Kasher

Operator Signature
DEB KASHUR

Remarks:

Interferent check 0.10 AC + acetone

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004198
Location = TOXL 8164.13.00 06/09
04/05/2012 16:32

Test	AC	Time
01 Diagnostic	OK	16:33
02 Room Air	RFI*	16:33
03 Room Air	0.000	16:33

*Invalid Test
 Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS921456
Lot No = 05911080A1
Cyl No = 14
Expiration Date = 02/01/2013
County = 08 Oper No. = 777777

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature
DEB KASHUR

Remarks:

RFI check

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004198
Location = TOXL 8164.13.00 06/09
04/04/2012 15:54

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:55
02 Std. Gas	0.080	15:55
03 Room Air	0.000	15:56
04 Std. Gas	0.080	15:56
05 Room Air	0.000	15:57
06 Std. Gas	0.080	15:57
07 Room Air	0.000	15:57

Lot No = 05911080A1
Cyl No = 14
Exp Date = 02/01/2013
County = 08 Oper No. = 777777

Deb Kasher

Operator Signature
DEB KASHUR

Remarks:

Calibration check 0.080 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004198
Location = TOXL 8164.13.00 06/09
04/04/2012 16:00

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	16:00
02 Std. Gas	0.079	16:01
03 Room Air	0.000	16:01
04 Std. Gas	0.079	16:01
05 Room Air	0.000	16:02
06 Std. Gas	0.080	16:02
07 Room Air	0.000	16:03

Lot No = 05911080A1

Cyl No = 14

Exp Date = 02/01/2013

County = 08

Oper No. = 777777

Deb Kashur

Operator Signature
DEB KASHUR

Remarks:

Calibration check 0.080 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004198
Location = TOXL 8164.13.00 06/09
04/04/2012 16:05

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	16:06
02 Std. Gas	0.080	16:06
03 Room Air	0.000	16:07
04 Std. Gas	0.080	16:07
05 Room Air	0.000	16:07
06 Std. Gas	0.080	16:08
07 Room Air	0.000	16:08

Lot No = 05911080A1
Cyl No = 14
Exp Date = 02/01/2013
County = 08 Oper No. = 777777

Deb Kashur
Operator Signature
DEB KASHUR

Remarks:
Calibration check 0.080 AC