

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-006499

Location: TOXL

Check When Done

- A. Is the breath tube heated? ✓
- B. Is the diagnostic check complete? ✓
- C. Is the time, date, and year correct (re-set if necessary)? ✓
- D. Print test. (Attach test record.) ✓
- E. Low AC. Use  $\leq 0.03$  AC in ACA mode. (Attach test record.)  
Sim SN DR5114 Lot No. 201505B Rep AC 0.030 ✓
- F. High AC. Use  $\geq 0.25$  AC in ACA mode. (Attach test record.)  
Sim SN MP3068 Lot No. 14102 Rep AC 0.400 ✓
- G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.  
(Attach test record.)  
Sim SN DR5135 Lot No. ICS4 Solution 0.10AC + 0.05% ACETONE ✓
- H. RFI check. Run CMS mode. Key radio during CMS test.  
(Attach test record.) ✓
- I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.  
(Attach test record.) 0.080 AC

Lot No. 13314080A3 Cyl No. 028 Exp. Date 6/5/2016

Test 1 0.082 Test 4 0.081 Test 7 0.081

Test 2 0.082 Test 5 0.082 Test 8 0.082

Test 3 0.080 Test 6 0.082 Test 9 0.082

Average 0.082

J. Tank pressure, Level 3 Func D Sub G P: 644 psi Reg 625 psi

K. Remarks/Maintenance record: ANNUAL INSPECTION

4.12.16  
Date

Charles E. Edr  
Field Inspector's Signature

12 Apr 16  
Date

Deb Shanaver  
Reviewed by

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-006499  
Location = TOXL                 8164.13.00 06/09  
04/12/2016                        09:30

\*\*\*\*\* Printer Test \*\*\*\*\*

abcdefghijklmnopqrstuvwxyz1234567890- = |  
ABCDEFGHIJKLMNopqrstuvwxyz!@#\$\$%^&\*()\_+?

abcdefghijklmnopqrstuvwxyz1234567890- = |  
ABCDEFGHIJKLMNopqrstuvwxyz!@#\$\$%^&\*()\_+?

Current Instrument Setup

Data Entry Mode:                Enabled  
Start Test Sequence:            DABACABA  
Display Prelim Rslt?            Yes  
Display Third Digit?            Yes  
Inhib Printer(Y/N)?            No  
Display Volume?                No  
Disable On Memfull?            Yes  
# of Print Copies?              1  
Select Std (D/W/I)?            Dry  
Standard Value?                0.080  
Standard Lot #?                 13314080A3  
Standard Cyl #?                28  
Standard Expiration?            06/05/2016  
Oper No?                        666666

Flow Cal. Date:                 08/18/2015  
Slope                            678  
Intercept                        -685002

IR Calibration Date:            08/18/2015  
                                  3um                        9um

-----  
0th Coef(\*100):                -22123                    -16014  
1st Coef(\*100):                264933                    138505  
2nd Coef(\*100):                3281                      1621  
H2O adj(mg/l\*10k):            704                        456

\*\*\*\*\* Printer Test End \*\*\*\*\*



Operator Signature  
CHARLES EDER

Remarks:

PRINT TEST

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006499  
Location = TOXL      8164.13.00 06/09  
04/12/2016      09:35

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:36
02 Std. Sol.	0.028	09:36
03 Room Air	0.000	09:37
04 Std. Sol.	0.029	09:38
05 Room Air	0.000	09:38
06 Std. Sol.	0.028	09:39
07 Room Air	0.000	09:39

08 Sim Temp = 34.0°C

Simul Ser No = DR5114  
Std Sol No = 201505B  
County = 08      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: *Low AC*

~~0.050~~ *0.030 AC*  
*CEE*

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006499  
Location = TOXL      8164.13.00 06/09  
04/12/2016      09:41

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:42
02 Std. Sol.	0.399	09:43
03 Room Air	0.000	09:43
04 Std. Sol.	0.401	09:44
05 Room Air	0.000	09:45
06 Std. Sol.	0.400	09:45
07 Room Air	0.000	09:46

08 Sim Temp = 34.0°C

Simul Ser No = MP3068  
Std Sol No = 14102  
County = 08      Oper No. = 666666

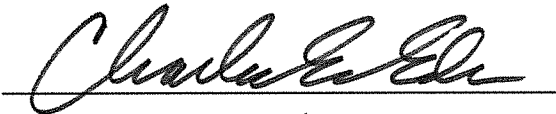

Operator Signature  
CHARLES EDER

Remarks:

HIGH AC  
0.400AC

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006499  
Location = TOXL      8164.13.00 06/09  
04/12/2016      09:48

Test	AC	Time
01 Room Air	0.000	09:49
02 *Subject Test	INT*	09:49
03 Room Air	0.000	09:50

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female      Weight = 150  
Test = DUI      Cit = INTERFERENT  
Dr. Lic. = ND/DIS821456  
Lot No = 13314080A3  
Cyl No = 28  
Expiration Date = 06/05/2016  
County = 08      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks:

*INTERFERENT CHECK*  
*0.10 AC + 0.05% ACETONE*

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006499  
Location = TOXL      8164.13.00 06/09  
04/12/2016      09:50

Test	AC	Time
01 Diagnostic	OK	09:51
02 Room Air	0.000	09:51
03 *Subject Test	RFI*	09:52
04 Room Air	0.000	09:52

\*Invalid Test  
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female      Weight = 150  
Test = DUI      Cit = RFI CHECK  
Dr. Lic. = ND/DIS821456  
Lot No = 13314080A3  
Cyl No = 28  
Expiration Date = 06/05/2016  
County = 08      Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature  
CHARLES EDER

Remarks:

*RFI CHECK*

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-006499  
Location = TOXL                    8164.13.00 06/09  
04/12/2016                            09:54

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:54
02 Std. Gas	0.082	09:55
03 Room Air	0.000	09:55
04 Std. Gas	0.082	09:56
05 Room Air	0.000	09:56
06 Std. Gas	0.080	09:56
07 Room Air	0.000	09:57

Lot No = 13314080A3  
Cyl No = 28  
Exp Date = 06/05/2016  
County = 08                            Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: *CALIBRATION CHECK*  
*0.080 AC*


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006499  
Location = TOXL      8164.13.00 06/09  
04/12/2016      09:57

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:58
02 Std. Gas	0.081	09:58
03 Room Air	0.000	09:59
04 Std. Gas	0.082	09:59
05 Room Air	0.000	10:00
06 Std. Gas	0.082	10:00
07 Room Air	0.000	10:01

Lot No = 13314080A3  
Cyl No = 28  
Exp Date = 06/05/2016  
County = 08      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
0.080 AC




Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-006499  
Location = TOXL      8164.13.00 06/09  
04/12/2016      10:01

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:01
02 Std. Gas	0.081	10:02
03 Room Air	0.000	10:02
04 Std. Gas	0.082	10:03
05 Room Air	0.000	10:03
06 Std. Gas	0.082	10:03
07 Room Air	0.000	10:04

Lot No = 13314080A3  
Cyl No = 28  
Exp Date = 06/05/2016  
County = 08      Oper No. = 666666

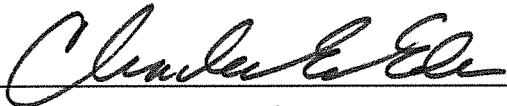  
Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
0.080 AC