

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-004197

Location: TOXL

- Check When Done
- A. Is the breath tube heated? ✓
- B. Is the diagnostic check complete? ✓
- C. Is the time, date, and year correct (re-set if necessary)? ✓
- D. Print test. (Attach test record.) ✓
- E. Low AC. Use <0.03 AC in ACA mode. (Attach test record.)
Sim SN DR3803 Lot No. 85 Rep AC 0.020
13070 ✓
- F. High AC. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN DR3374 Lot No. 13080 Rep AC 0.300 ✓
- G. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode.
(Attach test record.)
Sim SN DR3375 Lot No. FCS 2 Rep AC 0.10 AC + 0.05% acetone ✓
- H. RFI check. Run CMS mode. Key radio on first room air.
(Attach test record.) ✓
- I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.)

Lot No. 1651208042 Cyl No. 22 Exp. Date 08/01/2014

Test 1 <u>0.080</u>	Test 4 <u>0.080</u>	Test 7 <u>0.082</u>
Test 2 <u>0.081</u>	Test 5 <u>0.081</u>	Test 8 <u>0.081</u>
Test 3 <u>0.080</u>	Test 6 <u>0.081</u>	Test 9 <u>0.080</u>

Average 0.081

J. Tank pressure, Level 3 Func D Sub G P: 411 psi Reg 400 psi

K. Remarks/Maintenance record: Annual inspection, replaced simulator
return o-ring.

11 Apr 2014
Date

14 APR 14 ✓
Date

Deb Kashner
Field Inspector's Signature

Charles E. Ed...
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004197
Location = TOXL 8164.13.00 06/09
04/11/2014 12:46

***** Printer Test *****

abcdefghijklmnopqrstuvwxy1234567890-=_|
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&^*()_+?

abcdefghijklmnopqrstuvwxy1234567890-=_|
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&^*()_+?

Current Instrument Setup


Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 16512080A2
Standard Cyl #? 22
Standard Expiration? 08/01/2014
Oper No? 777777

Flow Cal. Date: 06/04/2009
Slope 697
Intercept -597078

IR Calibration Date: 06/03/2011

	3um	9um
0th Coef(*100):	2548	-942
1st Coef(*100):	242333	119472
2nd Coef(*100):	7141	2296
H2O adj(mg/l*10k):	665	778

***** Printer Test End *****

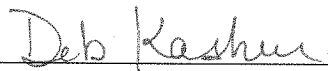

Operator Signature
DEB KASHUR

Remarks:



Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004197
Location = TOXL 8164.13.00 06/09
04/11/2014 12:57

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:57
02 Std. Sol.	0.023	12:58
03 Room Air	0.000	12:58
04 Std. Sol.	0.022	12:59
05 Room Air	0.000	13:00
06 Std. Sol.	0.023	13:00
07 Room Air	0.000	13:01

08 Sim Temp = 34.0°C

Simul Ser No = DR3803

Std Sol No = 13070

County = 08

Oper No. = 777777

Deb Kashur

Operator Signature
DEB KASHUR

Remarks:

Low AC 0.020AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004197
Location = TOXL 8164.13.00 06/09
04/11/2014 12:51

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:51
02 Std. Sol.	0.305	12:52
03 Room Air	0.000	12:53
04 Std. Sol.	0.303	12:53
05 Room Air	0.000	12:54
06 Std. Sol.	0.302	12:55
07 Room Air	0.000	12:55

08 Sim Temp = 34.0°C

Simul Ser No = DR3374

Std Sol No = 13080

County = 08

Oper No. = 777777

Deb Kasher

Operator Signature
DEB KASHUR

Remarks:

High AC 0.300AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004197
Location = TOXL 8164.13.00 06/09
04/11/2014 13:36

Test	AC	Time
01 Room Air	0.000	13:37
02 *Subject Test	INT*	13:39
03 Room Air	0.000	13:39

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS921456
Lot No = 16512080A2
Cyl No = 22
Expiration Date = 08/01/2014
County = 08 Oper No. = 777777

Deb Kashur

Operator Signature
DEB KASHUR

Remarks:

DK
Interferent Interferent check
0.10AC + 0.05% acetone

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004197
Location = TOXL 8164.13.00 06/09
04/11/2014 13:40

Test	AC	Time
01 Diagnostic	OK	13:41
02 Room Air	RFI*	13:41
03 Room Air	0.000	13:41

*Invalid Test
 Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS921456
Lot No = 16512080A2
Cyl No = 22
Expiration Date = 08/01/2014
County = 08 Oper No. = 777777

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

Deb Kasher

Operator Signature
DEB KASHUR

Remarks:

RFI test

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004197
Location = TOXL 8164.13.00 06/09
04/11/2014 10:21

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	10:22
02 Std. Gas	0.080	10:22
03 Room Air	0.000	10:23
04 Std. Gas	0.081	10:23
05 Room Air	0.000	10:23
06 Std. Gas	0.080	10:24
07 Room Air	0.000	10:24

Lot No = 16512080A2

Cyl No = 22

Exp Date = 08/01/2014

County = 08

Oper No. = 777777



Operator Signature
DEB KASHUR

Remarks:

Calibration check 0.080AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004197
Location = TOXL 8164.13.00 06/09
04/11/2014 11:28

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:29
02 Std. Gas	0.080	11:29
03 Room Air	0.000	11:30
04 Std. Gas	0.081	11:30
05 Room Air	0.000	11:30
06 Std. Gas	0.081	11:31
07 Room Air	0.000	11:31


Lot No = 16512080A2

Cyl No = 22

Exp Date = 08/01/2014

County = 08

Oper No. = 777777

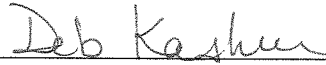

Operator Signature
DEB KASHUR

Remarks:

Calibration check 0.080 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004197
Location = TOXL 8164.13.00 06/09
04/11/2014 12:46

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:46
02 Std. Gas	0.082	12:47
03 Room Air	0.000	12:47
04 Std. Gas	0.081	12:48
05 Room Air	0.000	12:48
06 Std. Gas	0.080	12:49
07 Room Air	0.000	12:49

Lot No = 16512080A2
Cyl No = 22
Exp Date = 08/01/2014
County = 08 Oper No. = 777777

Deb Kasher
Operator Signature
DEB KASHUR

Remarks:

Calibration check 0.080 Ac

Form 106-I8000