



**INTOXILYZER® 8000 INSTALLATION AND REPAIR CHECKOUT**  
 OFFICE OF ATTORNEY GENERAL  
 CRIME LABORATORY DIVISION  
 Toxicology Section/Breath Alcohol Program  
 SFN 59281 (03-2014)

Serial Number <b>80-004956</b>	Instrument Location <b>WASH (McLEAN County SD)</b>
Reason for Install/Repair <input type="checkbox"/> Install After Receiving From Crime Laboratory <input type="checkbox"/> Install After Location Change <input checked="" type="checkbox"/> Other (Specify) <b>INSTALL DURING ANNUAL INSPECTION</b>	

Check When Done:

- 1. Surge Protector Installed/Properly Grounded.
- 2. Telephone Line Connected to Intoxilyzer® 8000.
- 3. Breath Tube Heated.
- 4. Enter Preliminary Data (i.e. Date, Time, DST (Y), and Location; Level 2, Function E).
- 5. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 6. Run Tests:
  - A. Print Test (Level 1, Function P). ✓
  - B. ACA Test (Level 1, Function C). ✓
  - C. RFI Test (CMS Mode or Level 1, Function B or C; Key Radio During Test). ✓
- 7. Repair and/or Maintenance Performed (if Any): **NA - INSTALL DURING ANNUAL INSPECTION**
- 8. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 9. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 10. Send the Following to the Crime Laboratory: Completed Intoxilyzer®8000 Installation and Repair Checkout (SFN59281, Form 104-G), Print Test, ACA Test, and RFI Test.

Field Inspector Signature <i>Charles E. Edler</i>	Date <b>4.14.15 ✓</b>
Reviewed-By (Crime Laboratory Use Only) <i>Deb Kashner</i>	Date <b>17 Apr 15</b>

This record on file at the Office of Attorney General, Crime Laboratory Division, in the County of Burleigh, North Dakota, is certified to be a true and correct copy of the documents received.

Certified By: *Charles E. Edler* Certified Date: **30 Nov. 2016**




Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004956  
Location = WASH      8164.13.00 06/09  
04/14/2015      14:11

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:11
02 Std. Gas	0.080	14:12
03 Room Air	0.000	14:12
04 Std. Gas	0.079	14:12
05 Room Air	0.000	14:13
06 Std. Gas	0.080	14:13
07 Room Air	0.000	14:14

Lot No = 13314080A3  
Cyl No = 30  
Exp Date = 06/05/2016  
County = 28      Oper No. = 666666

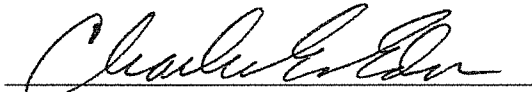

Operator Signature  
CHARLES EDER

Remarks:      *ACA TEST*  
                    *0.080 AC*

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004956  
Location = WASH      8164.13.00 06/09  
04/14/2015      14:14

Test	AC	Time
01 Diagnostic	OK	14:15
02 Room Air	0.000	14:15
03 *Subject Test	RFI*	14:16
04 Room Air	0.000	14:16

\*Invalid Test  
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female      Weight = 150  
Test = DUI      Cit = RFI CHECK  
Dr. Lic. = ND/DIS821456  
Lot No = 13314080A3  
Cyl No = 30  
Expiration Date = 06/05/2016  
County = 28      Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature  
CHARLES EDER

Remarks:

RFI CHECK