

INTOXILYZER 8000 INSPECTION

Serial No.: 80-004937

Location: TOX LAB

- | | Check When Done |
|---|-----------------|
| A. Is the warm-up time less than 20 minutes? | <u>✓</u> |
| B. Is a three-pronged grounded outlet used? | <u>✓</u> |
| C. Is the breath tube heated? | <u>✓</u> |
| D. Is the diagnostic check complete? | <u>✓</u> |
| E. Is the time, date, and year correct (re-set if necessary)? | <u>✓</u> |
| F. Print test. (Attach test record.) | <u>✓</u> |
| G. Low AC. Use <0.03 AC in ACA mode. (Attach test record.)
Sim SN <u>DR5140</u> Lot No. <u>—</u> Rep AC <u>0.019 AC</u> | <u>✓</u> |
| H. Linearity Test. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN <u>DR5113</u> Lot No. <u>—</u> Rep AC <u>0.285 AC</u> | <u>✓</u> |
| I. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode.
(Attach test record.)
Sim SN <u>DR5135</u> Lot No. <u>—</u> Rep AC <u>0.10 AC + 0.05% Acetone</u> | <u>✓</u> |
| J. RFI check. Run CMS mode. Key radio on first room air.
(Attach test record.) | <u>✓</u> |
| K. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) | |

Lot No. <u>709352</u>	Cyl No. <u>18</u>	Exp. Date <u>12/23/2011</u>
Test 1 <u>0.078</u>	Test 4 <u>0.079</u>	Test 7 <u>0.079</u>
Test 2 <u>0.079</u>	Test 5 <u>0.079</u>	Test 8 <u>0.078</u>
Test 3 <u>0.079</u>	Test 6 <u>0.078</u>	Test 9 <u>0.079</u>
Average <u>0.079</u>		

L. Tank pressure, Level 3 Func D Sub G P: 303 psi Reg 300 psi

M. Remarks/Maintenance record: Annual Inspection

16 April 2010
Date

Mary Pearson
Field Inspector's Signature

19 Apr 2010
Date

Deb Kashner
Reviewed by

80-004937
3/4/2 km

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004937
Location = TOXL 8164.13.00 06/09
04/16/2010 11:40

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&^*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&^*()_+?

Current Instrument Setup


Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Wet
Standard Value? 0.000
Standard Solution #? 000000
Sim Serial #? 00000
Oper No? 107501

Flow Cal. Date: 01/18/2010
Slope 730
Intercept -779991

IR Calibration Date: 02/10/2010
 3um 9um

0th Coef(*100): -23590 -32037
1st Coef(*100): 262965 133872
2nd Coef(*100): -566 311
H2O adj(mg/l*10k): 569 537

***** Printer Test End *****

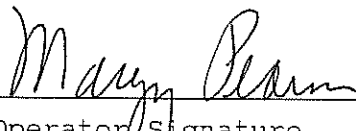

Operator Signature
MARGY PEARSON

Remarks:

Print

Form 106-I8000

80-004937
35/12 KM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004937
Location = TOXL 8164.13.00 06/09
04/16/2010 11:25

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:26
02 Std. Sol.	0.016	11:27
03 Room Air	0.000	11:28
04 Std. Sol.	0.016	11:29
05 Room Air	0.000	11:29
06 Std. Sol.	0.017	11:30
07 Room Air	0.000	11:30


08 Sim Temp = 34.0°C

Simul Ser No = 00000

Std Sol No = 000000

County = 08

Oper No. = 107501

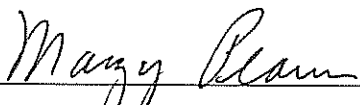

Operator Signature
MARGY PEARSON

Remarks:

LOW AC - 0.019 AC

Form 106-I8000

80-004937
3/6/12 PR

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004937
Location = TOXL 8164.13.00 06/09
04/16/2010 11:32

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:33
02 Std. Sol.	0.277	11:35
03 Room Air	0.000	11:35
04 Std. Sol.	0.278	11:36
05 Room Air	0.000	11:37
06 Std. Sol.	0.278	11:37
07 Room Air	0.000	11:38


08 Sim Temp = 34.0°C

Simul Ser No = 00000

Std Sol No = 000000

County = 08

Oper No. = 107501

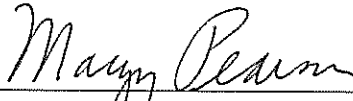

Operator Signature
MARGY PEARSON

Remarks:

Linearity 0.285 Ac

Form 106-I8000

80-004937
37/42 KM


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004937
Location = TOXL 8164.13.00 06/09
04/16/2010 11:41

Test	AC	Time
01 Room Air	0.000	11:41
02 *Subject Test	INT*	11:42
03 Room Air	0.000	11:42

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = 00000
Std Sol No = 000000
County = 08 Oper No. = 107501

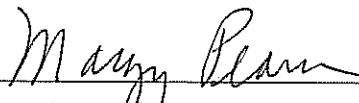

Operator Signature
MARGY PEARSON

Remarks:

Aic + Acetone

Form 106-I8000

80-004937
38/42 *MP*

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004937
Location = TOXL 8164.13.00 06/09
04/16/2010 11:44

Test	AC	Time
01 Room Air	RFI*	11:45
02 Room Air	0.000	11:45

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = 00000
Std Sol No = 000000
County = 08 Oper No. = 107501



Operator Signature
MARGY PEARSON

Remarks:

RFI

Form 106-I8000

80-004937
3/4/2 KH


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004937
Location = TOXL 8164.13.00 06/09
04/16/2010 11:09

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:10
02 Std. Gas	0.078	11:10
03 Room Air	0.000	11:11
04 Std. Gas	0.079	11:11
05 Room Air	0.000	11:12
06 Std. Gas	0.079	11:12
07 Room Air	0.000	11:12

Lot No = 709352
Cyl No = 18
Exp Date = 12/23/2011
County = 08 Oper No. = 107501

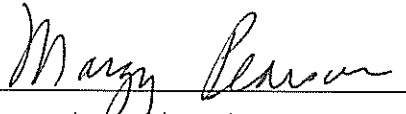

Operator Signature
MARGY PEARSON

Remarks:

Cal Check

Form 106-I8000

80-004937
4/16/10


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004937
Location = TOXL 8164.13.00 06/09
04/16/2010 11:13

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:14
02 Std. Gas	0.079	11:14
03 Room Air	0.000	11:14
04 Std. Gas	0.079	11:15
05 Room Air	0.000	11:15
06 Std. Gas	0.078	11:16
07 Room Air	0.000	11:16

Lot No = 709352
Cyl No = 18
Exp Date = 12/23/2011
County = 08 Oper No. = 107501

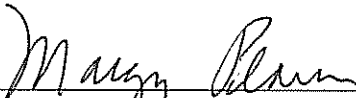

Operator Signature
MARGY PEARSON

Remarks:

Cal Check

Form 106-I8000

80-004937
4/16/12 FM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004937
Location = TOXL 8164.13.00 06/09
04/16/2010 11:17

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:18
02 Std. Gas	0.079	11:18
03 Room Air	0.000	11:19
04 Std. Gas	0.078	11:19
05 Room Air	0.000	11:20
06 Std. Gas	0.079	11:20
07 Room Air	0.000	11:21

Lot No = 709352
Cyl No = 18
Exp Date = 12/23/2011
County = 08 Oper No. = 107501



Operator Signature
MARGY PEARSON

Remarks:

Cal Check

Form 106-I8000

80-004937
4/2/12 SM