

INTOXILYZER 8000 INSPECTION

Serial No.: 80-002669

Location: TOXL

- | | Check When Done |
|--|-----------------|
| A. Is the warm-up time less than 20 minutes? | <u>✓</u> |
| B. Is a three-pronged grounded outlet used? | <u>✓</u> |
| C. Is the breath tube heated? | <u>✓</u> |
| D. Is the diagnostic check complete? | <u>✓</u> |
| E. Is the time, date, and year correct (re-set if necessary)? | <u>✓</u> |
| F. Print test. (Attach test record.) | <u>✓</u> |
| G. Low AC. Use <0.03 AC in ACA mode. (Attach test record.)
Sim SN <u>DR3455</u> Lot No. <u>10290</u> Rep AC <u>0.020</u> | <u>✓</u> |
| H. Linearity Test. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN <u>DR3847</u> Lot No. <u>11110</u> Rep AC <u>0.300</u> | <u>✓</u> |
| I. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode.
(Attach test record.)
Sim SN <u>DR3377</u> Lot No. <u> </u> Rep AC <u>0.100</u> | <u>✓</u> |
| J. RFI check. Run CMS mode. Key radio on first room air.
(Attach test record.) | <u>✓</u> |
| K. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) | |

Lot No. <u>23411080A1</u>	Cyl No. <u>55</u>	Exp. Date <u>10/01/2013</u>
Test 1 <u>0.082</u>	Test 4 <u>0.082</u>	Test 7 <u>0.082</u>
Test 2 <u>0.081</u>	Test 5 <u>0.082</u>	Test 8 <u>0.082</u>
Test 3 <u>0.081</u>	Test 6 <u>0.082</u>	Test 9 <u>0.082</u>
Average <u>0.082</u>		

L. Tank pressure, Level 3 Func D Sub G P: 958 psi Reg 950 psi

M. Remarks/Maintenance record: ANNUAL INSPECTION
13 APR 12 Date Charles E. [Signature] Field Inspector's Signature
16 Apr 12 Date Deb Kashner Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.13.00 06/09
04/13/2012 16:18

***** Printer Test *****

abcdefghijklmnopqrstuvwxy1234567890-==|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxy1234567890-==|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

Current Instrument Setup

Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Wet
Standard Value? 0.000
Standard Solution #? 0000
Sim Serial #? 000000
Oper No? 666666

Flow Cal. Date: 09/09/2011
Slope 662
Intercept -510765

IR Calibration Date: 09/09/2011
 3um 9um

0th Coef(*100): -15623 -9015
1st Coef(*100): 263092 133867
2nd Coef(*100): 3130 1285
H2O adj(mg/l*10k): 642 392

***** Printer Test End *****



Operator Signature
CHARLES EDER

Remarks:

PRINT TEST

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.13.00 06/09
04/13/2012 16:18

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	16:19
02 Std. Sol.	0.020	16:20
03 Room Air	0.000	16:20
04 Std. Sol.	0.020	16:21
05 Room Air	0.000	16:21
06 Std. Sol.	0.020	16:22
07 Room Air	0.000	16:23

08 Sim Temp = 34.0°C

Simul Ser No = 000000

Std Sol No = 0000

County = 08

Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

Low AC CHECK
0.020 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.13.00 06/09
04/13/2012 16:31

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	16:32
02 Std. Sol.	0.298	16:33
03 Room Air	0.000	16:33
04 Std. Sol.	0.300	16:34
05 Room Air	0.000	16:35
06 Std. Sol.	0.300	16:35
07 Room Air	0.000	16:36

08 Sim Temp = 34.0°C

Simul Ser No = 000000

Std Sol No = 0000

County = 08

Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

LINEARITY TEST
0.300 AC

Form 106-I8000


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.13.00 06/09
04/13/2012 16:38

Test	AC	Time
01 Room Air	0.000	16:39
02 *Subject Test	INT*	16:40
03 Room Air	0.000	16:40

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = 000000
Std Sol No = 0000
County = 08 Oper No. = 666666

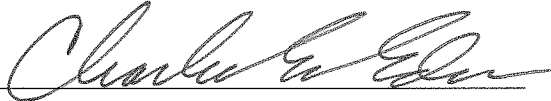

Operator Signature
CHARLES EDER

Remarks:

*INTERFERENT CHECK
0.100 AC + ACETONE*

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.13.00 06/09
04/13/2012 16:42

Test	AC	Time
01 Room Air	RFI*	16:43
02 Room Air	0.000	16:43

*Invalid Test
 Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Lot No = 23411080A1
Cyl No = 55
Expiration Date = 10/01/2013
County = 08 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

RFI CHECK

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.13.00 06/09
04/13/2012 16:44

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	16:45
02 Std. Gas	0.082	16:45
03 Room Air	0.000	16:46
04 Std. Gas	0.081	16:46
05 Room Air	0.000	16:47
06 Std. Gas	0.081	16:47
07 Room Air	0.000	16:48


Lot No = 23411080A1

Cyl No = 55

Exp Date = 10/01/2013

County = 08

Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK

0.080 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.13.00 06/09
04/13/2012 16:48

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	16:49
02 Std. Gas	0.082	16:49
03 Room Air	0.000	16:50
04 Std. Gas	0.082	16:50
05 Room Air	0.000	16:51
06 Std. Gas	0.082	16:51
07 Room Air	0.000	16:52

Lot No = 23411080A1
Cyl No = 55
Exp Date = 10/01/2013
County = 08 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: *CALIBRATION CHECK*
 0.080 AC


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-002669
Location = TOXL 8164.13.00 06/09
04/13/2012 16:52

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	16:52
02 Std. Gas	0.082	16:53
03 Room Air	0.000	16:53
04 Std. Gas	0.082	16:54
05 Room Air	0.000	16:54
06 Std. Gas	0.082	16:54
07 Room Air	0.000	16:55

Lot No = 23411080A1
Cyl No = 55
Exp Date = 10/01/2013
County = 08 Oper No. = 666666

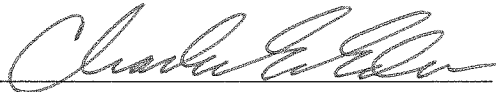

Operator Signature
CHARLES EDER

Remarks:

CALIBRATION CHECK
0.080A-C

Form 106-I8000