

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-004952

Location: HETT (ADAMS Co. S.O.)

- Check When Done
- A. Is the breath tube heated? ✓
 - B. Is the diagnostic check complete? ✓
 - C. Is the time, date, and year correct (re-set if necessary)? ✓
 - D. Print test. (Attach test record.) ✓
 - E. Low AC. Use ≤ 0.03 AC in ACA mode. (Attach test record.)
Sim SN DR5143 Lot No. 201501B Rep AC 0.020 ✓
 - F. High AC. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN DR3455 Lot No. 14270 Rep AC 0.300 ✓
 - G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.
(Attach test record.)
Sim SN DR5132 Lot No. ICS # 2 Solution 0.10AC + 0.05% ACE ✓
 - H. RFI check. Run CMS mode. Key radio during CMS test.
(Attach test record.) ✓
 - I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) 0.080 AC

Lot No. 33913080A4 Cyl No. 13 Exp. Date 1/1/2016
Test 1 0.079 Test 4 0.080 Test 7 0.079
Test 2 0.079 Test 5 0.080 Test 8 0.079
Test 3 0.080 Test 6 0.080 Test 9 0.080
Average 0.080

J. Tank pressure, Level 3 Func D Sub G P: 1006 psi Reg 1000+ psi

K. Remarks/Maintenance record: ANNUAL INSPECTION

4.22.15
Date

Charles E. Edwards
Field Inspector's Signature

28 Apr 2015 ✓
Date

Deb Kashner
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004952
Location = HETT 8164.13.00 06/09
04/22/2015 11:20

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&^*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNopqrstuvwxyz!@#%&^*()_+?

Current Instrument Setup

Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 33913080A4
Standard Cyl #? 13
Standard Expiration? 01/01/2016
Oper No? 666666

Flow Cal. Date: 01/19/2010
Slope 694
Intercept -638616

IR Calibration Date: 01/19/2010
 3um 9um

0th Coef(*100): -32832 -34603
1st Coef(*100): 262566 131011
2nd Coef(*100): 1418 975
H2O adj(mg/l*10k): 686 550

***** Printer Test End *****



Operator Signature
CHARLES EDER

Remarks:

PRINT TEST

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004952
Location = HETT 8164.13.00 06/09
04/22/2015 11:41

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:42
02 Std. Sol.	0.017	11:43
03 Room Air	0.000	11:43
04 Std. Sol.	0.018	11:44
05 Room Air	0.000	11:44
06 Std. Sol.	0.017	11:45
07 Room Air	0.000	11:45


08 Sim Temp = 34.0°C

Simul Ser No = DR5143

Std Sol No = 201501B

County = 01

Oper No. = 666666

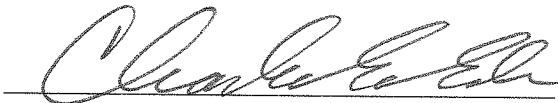

Operator Signature
CHARLES EDER

Remarks:

Low AC

0.020 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004952
Location = HETT 8164.13.00 06/09
04/22/2015 11:51

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:52
02 Std. Sol.	0.293	11:53
03 Room Air	0.000	11:53
04 Std. Sol.	0.294	11:54
05 Room Air	0.000	11:55
06 Std. Sol.	0.294	11:55
07 Room Air	0.000	11:56

08 Sim Temp = 34.0°C

Simul Ser No = DR3455

Std Sol No = 14270

County = 01

Oper No. = 666666



Operator Signature

CHARLES EDER

Remarks:

HIGH AC

0.300 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004952
Location = HETT 8164.13.00 06/09
04/22/2015 11:33

Test	AC	Time
01 Room Air	0.000	11:34
02 *Subject Test	INT*	11:35
03 Room Air	0.000	11:36

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = INTERFERENT
Dr. Lic. = ND/DIS821456
Lot No = 33913080A4
Cyl No = 13
Expiration Date = 01/01/2016
County = 01 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

Interferent Check
PER ~~0.0~~ 0.10AC + 0.05% ACETONE

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004952
Location = HETT 8164.13.00 06/09
04/22/2015 11:22

Test	AC	Time
01 Diagnostic	OK	11:24
02 Room Air	RFI*	11:24
03 Room Air	0.000	11:24

*Invalid Test
 Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = RFI CHECK
Dr. Lic. = ND/DIS821456
Lot No = 33913080A4
Cyl No = 13
Expiration Date = 01/01/2016
County = 01 Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

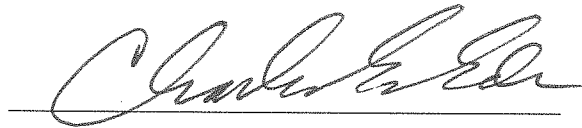

Operator Signature
CHARLES EDER

Remarks: *RFI CHECK*


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004952
Location = HETT 8164.13.00 06/09
04/22/2015 11:26

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:27
02 Std. Gas	0.079	11:27
03 Room Air	0.000	11:28
04 Std. Gas	0.079	11:28
05 Room Air	0.000	11:29
06 Std. Gas	0.080	11:29
07 Room Air	0.000	11:29

Lot No = 33913080A4
Cyl No = 13
Exp Date = 01/01/2016
County = 01 Oper No. = 666666

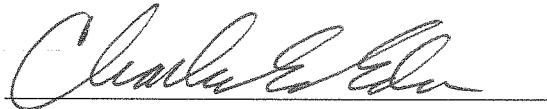

Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
0.080 AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004952
Location = HETT 8164.13.00 06/09
04/22/2015 11:30

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:31
02 Std. Gas	0.080	11:31
03 Room Air	0.000	11:31
04 Std. Gas	0.080	11:32
05 Room Air	0.000	11:32
06 Std. Gas	0.080	11:33
07 Room Air	0.000	11:33

Lot No = 33913080A4
Cyl No = 13
Exp Date = 01/01/2016
County = 01 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
0.080 AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004952
Location = HETT 8164.13.00 06/09
04/22/2015 11:57

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:57
02 Std. Gas	0.079	11:58
03 Room Air	0.000	11:58
04 Std. Gas	0.079	11:59
05 Room Air	0.000	11:59
06 Std. Gas	0.080	11:59
07 Room Air	0.000	12:00

Lot No = 33913080A4

Cyl No = 13

Exp Date = 01/01/2016

County = 01

Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

CALIBRATION CHECK
0.080 AC

Form 106-I8000