

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-005360

Location: ASHL (McINTOSH Co. S.O.)

- Check When Done
- A. Is the breath tube heated? ✓
  - B. Is the diagnostic check complete? ✓
  - C. Is the time, date, and year correct (re-set if necessary)? ✓
  - D. Print test. (Attach test record.) ✓
  - E. Low AC. Use  $\leq 0.03$  AC in ACA mode. (Attach test record.)  
Sim SN DR5143 Lot No. 201501 B Rep AC 0.020 ✓
  - F. High AC. Use  $\geq 0.25$  AC in ACA mode. (Attach test record.)  
Sim SN DR3455 Lot No. 14270 Rep AC 0.300 ✓
  - G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.  
(Attach test record.)  
Sim SN DR5132 Lot No. Ics#2 Solution 0.10AC+0.05%ACE ✓
  - H. RFI check. Run CMS mode. Key radio during CMS test.  
(Attach test record.) ✓
  - I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.  
(Attach test record.) 0.080 AC

Lot No. 33913080A4 Cyl No. 41 Exp. Date 1/1/16  
Test 1 0.083 Test 4 0.082 Test 7 0.084  
Test 2 0.082 Test 5 0.082 Test 8 0.083  
Test 3 0.081 Test 6 0.083 Test 9 0.083  
Average 0.083

J. Tank pressure, Level 3 Func D Sub G P: 1015 psi Reg 1000+ psi

K. Remarks/Maintenance record: ANNUAL INSPECTION

4.21.15

Date

28 Apr 2015 ✓

Date

Charles E. Edr

Field Inspector's Signature

Deb Kashner

Reviewed by

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-005360  
Location = ASHL                8164.13.00 06/09  
04/21/2015                      12:18

\*\*\*\*\* Printer Test \*\*\*\*\*

abcdefghijklmnopqrstuvwxyz1234567890-==|  
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#%&^\*()\_+?

abcdefghijklmnopqrstuvwxyz1234567890-==|  
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#%&^\*()\_+?

Current Instrument Setup

Data Entry Mode:            Enabled  
Start Test Sequence:        DABACABA  
Display Prelim Rslt?        Yes  
Display Third Digit?        Yes  
Inhib Printer(Y/N)?        No  
Display Volume?            No  
Disable On Memfull?        Yes  
# of Print Copies?         1  
Select Std (D/W/I)?        Dry  
Standard Value?            0.080  
Standard Lot #?            33913080A4  
Standard Cyl #?            41  
Standard Expiration?       01/01/2016  
Oper No?                    666666

Flow Cal. Date:            08/17/2011  
Slope                      639  
Intercept                  -343114

IR Calibration Date:       08/17/2011  
                              3um                    9um

-----  
0th Coef(\*100):        -17266                -25747  
1st Coef(\*100):        263773                142804  
2nd Coef(\*100):        2606                  1406  
H2O adj(mg/l\*10k):    499                    577

\*\*\*\*\* Printer Test End \*\*\*\*\*



Operator Signature  
CHARLES EDER

Remarks:

PRINT TEST

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-005360  
Location = ASHL      8164.13.00 06/09  
04/21/2015      12:32

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:32
02 Std. Sol.	0.017	12:33
03 Room Air	0.000	12:34
04 Std. Sol.	0.017	12:34
05 Room Air	0.000	12:35
06 Std. Sol.	0.018	12:36
07 Room Air	0.000	12:36

08 Sim Temp = 34.0°C

Simul Ser No = DR5143  
Std Sol No = ~~14270~~ *CEE 201501B*  
County = 26      Oper No. = 666666

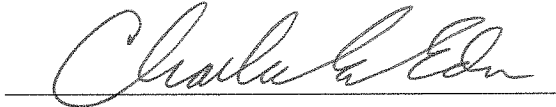

Operator Signature  
CHARLES EDER

Remarks:

*Low AC*  
*0.020 AC*

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-005360  
Location = ASHL                    8164.13.00 06/09  
04/21/2015                            12:42

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:43
02 Std. Sol.	0.301	12:44
03 Room Air	0.000	12:44
04 Std. Sol.	0.301	12:45
05 Room Air	0.000	12:45
06 Std. Sol.	0.302	12:46
07 Room Air	0.000	12:47


08 Sim Temp = 34.0°C

Simul Ser No = DR3455

Std Sol No = 14270

County = 26

Oper No. = 666666

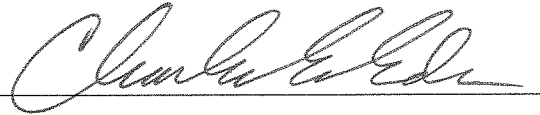

Operator Signature  
CHARLES EDER

Remarks:

HIGH AC  
0.300 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-005360  
Location = ASHL                 8164.13.00 06/09  
04/21/2015                         12:24

Test	AC	Time
01 Room Air	0.000	12:25
02 *Subject Test	INT*	12:26
03 Room Air	0.000	12:26

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                     Weight = 150  
Test = DUI                             Cit = INTERFERENT  
Dr. Lic. = ND/DIS821456  
Lot No = 33913080A4  
Cyl No = 41  
Expiration Date = 01/01/2016  
County = 26                             Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks:

*Interferent Check*  
*0.10 AC + 0.05% ACETONE*

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-005360  
Location = ASHL                 8164.13.00 06/09  
04/21/2015                        12:10

Test	AC	Time
01 Diagnostic	OK	12:11
02 Room Air	RFI*	12:12
03 Room Air	0.000	12:12

\*Invalid Test  
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                 Weight = 150  
Test = DUI                         Cit = RFI CHECK  
Dr. Lic. = ND/DIS821456  
Lot No = 33913080A4  
Cyl No = 41  
Expiration Date = 01/01/2016  
County = 26                        Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

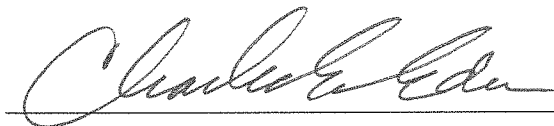

Operator Signature  
CHARLES EDER

Remarks:

RFI CHECK


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-005360  
Location = ASHL      8164.13.00 06/09  
04/21/2015      12:20

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:21
02 Std. Gas	0.083	12:21
03 Room Air	0.000	12:22
04 Std. Gas	0.082	12:22
05 Room Air	0.000	12:22
06 Std. Gas	0.081	12:23
07 Room Air	0.000	12:23

Lot No = 33913080A4  
Cyl No = 41  
Exp Date = 01/01/2016  
County = 26      Oper No. = 666666

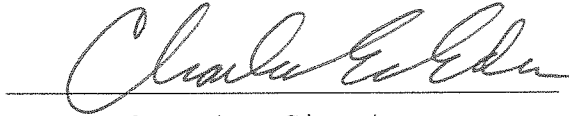

Operator Signature  
CHARLES EDER

Remarks:

*Calibration Check*  
*0.080 AC*

Form 106-I8000


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-005360  
Location = ASHL                8164.13.00 06/09  
04/21/2015                        12:47

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:48
02 Std. Gas	0.082	12:48
03 Room Air	0.000	12:49
04 Std. Gas	0.082	12:49
05 Room Air	0.000	12:50
06 Std. Gas	0.083	12:50
07 Room Air	0.000	12:50

Lot No = 33913080A4  
Cyl No = 41  
Exp Date = 01/01/2016  
County = 26                            Oper No. = 666666

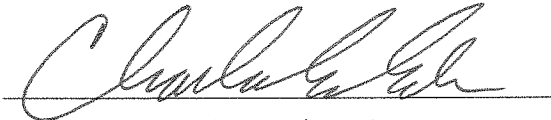

Operator Signature  
CHARLES EDER

Remarks:

*Calibration Check*  
*0.080 AC*

Form 106-I8000




Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-005360  
Location = ASHL                    8164.13.00 06/09  
04/21/2015                            12:51

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:51
02 Std. Gas	0.084	12:52
03 Room Air	0.000	12:52
04 Std. Gas	0.083	12:52
05 Room Air	0.000	12:53
06 Std. Gas	0.083	12:53
07 Room Air	0.000	12:54

Lot No = 33913080A4  
Cyl No = 41  
Exp Date = 01/01/2016  
County = 26                            Oper No. = 666666

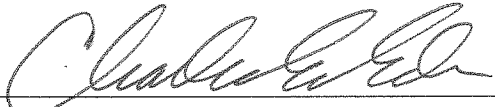  
Operator Signature  
CHARLES EDER

Remarks: *Geto CEE CALIBRATION Check*  
*0.080AC*