

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-005952

Location: BCSO

- Check When Done
- A. Is the breath tube heated? ✓
 - B. Is the diagnostic check complete? ✓
 - C. Is the time, date, and year correct (re-set if necessary)? ✓
 - D. Print test. (Attach test record.) ✓
 - E. Low AC. Use ≤ 0.03 AC in ACA mode. (Attach test record.)
Sim SN DR3803 Lot No. 15160 Rep AC 0.020 ✓
 - F. High AC. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN DR3847 Lot No. 16060 Rep AC 0.300 ✓
 - G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.
(Attach test record.)
Sim SN DR5190 Lot No. ICS4 Solution 0.10 AC + 0.05% ACETONE ✓
 - H. RFI check. Run CMS mode. Key radio during CMS test.
(Attach test record.) ✓
 - I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) 0.080 AC

Lot No. 13314080A3 Cyl No. 025 Exp. Date 6/5/2016

Test 1 0.080 Test 4 0.080 Test 7 0.081

Test 2 0.080 Test 5 0.081 Test 8 0.081

Test 3 0.081 Test 6 0.081 Test 9 0.081

Average 0.081

J. Tank pressure, Level 3 Func D Sub G P: 600 psi Reg 600 psi

K. Remarks/Maintenance record: ANNUAL INSPECTION

4-26-2016
Date

28 Apr 16
Date

Charles E. Edwards
Field Inspector's Signature

Deb Shanauer
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005952
Location = BCSO 8164.13.00 06/09
04/26/2016 15:40

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNopqrstuvwxyz!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNopqrstuvwxyz!@#\$\$%^&*()_+?

Current Instrument Setup

Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 13314080A3
Standard Cyl #? 25
Standard Expiration? 06/05/2016
Oper No? 666666

Flow Cal. Date: 06/21/2013
Slope 734
Intercept -865305

IR Calibration Date: 06/21/2013
 3um 9um

	3um	9um
0th Coef(*100):	-21281	-23057
1st Coef(*100):	292132	137092
2nd Coef(*100):	2334	1231
H2O adj(mg/l*10k):	497	482

***** Printer Test End *****



Operator Signature
CHARLES EDER

Remarks:

PRINT TEST

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005952
Location = BCSO 8164.13.00 06/09
04/26/2016 15:46

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:47
02 Std. Sol.	0.019	15:48
03 Room Air	0.000	15:48
04 Std. Sol.	0.018	15:49
05 Room Air	0.000	15:49
06 Std. Sol.	0.018	15:50
07 Room Air	0.000	15:51

08 Sim Temp = 34.0°C

Simul Ser No = DR3803
Std Sol No = 15160
County = 08 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: *Low AC*
0.020 AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005952
Location = BCSO 8164.13.00 06/09
04/26/2016 15:52

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:52
02 Std. Sol.	0.305	15:53
03 Room Air	0.000	15:54
04 Std. Sol.	0.306	15:55
05 Room Air	0.000	15:55
06 Std. Sol.	0.306	15:56
07 Room Air	0.000	15:56

08 Sim Temp = 34.0°C

Simul Ser No = DR3847

Std Sol No = 16060

County = 08

Oper No. = 666666



Operator Signature

CHARLES EDER

Remarks:

HIGH AC
0.300 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005952
Location = BCSO 8164.13.00 06/09
04/26/2016 15:41

Test	AC	Time
01 Room Air	0.000	15:42
02 *Subject Test	INT*	15:43
03 Room Air	0.000	15:43

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = INTERFERENT
Dr. Lic. = ND/DIS821456
Lot No = 13314080A3
Cyl No = 25
Expiration Date = 06/05/2016
County = 08 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: *INTERFERENT CHECK*
0.10AC + 0.05% ACETONE

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005952
Location = BCSO 8164.13.00 06/09
04/26/2016 15:38

Test	AC	Time
01 Diagnostic	OK	15:39
02 Room Air	0.000	15:39
03 *Subject Test	RFI*	15:39
04 Room Air	0.000	15:40

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = RFI CHECK
Dr. Lic. = ND/DIS821456
Lot No = 13314080A3
Cyl No = 25
Expiration Date = 06/05/2016
County = 08 Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature
CHARLES EDER

Remarks:

RFI CHECK

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005952
Location = BCSO 8164.13.00 06/09
04/26/2016 15:59

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:59
02 Std. Gas	0.080	16:00
03 Room Air	0.000	16:00
04 Std. Gas	0.080	16:00
05 Room Air	0.000	16:01
06 Std. Gas	0.081	16:01
07 Room Air	0.000	16:02

Lot No = 13314080A3
Cyl No = 25
Exp Date = 06/05/2016
County = 08 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
0.080 AC

Form 106-I8000


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005952
Location = BCSO 8164.13.00 06/09
04/26/2016 16:02

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	16:03
02 Std. Gas	0.080	16:03
03 Room Air	0.000	16:03
04 Std. Gas	0.081	16:04
05 Room Air	0.000	16:04
06 Std. Gas	0.081	16:05
07 Room Air	0.000	16:05

Lot No = 13314080A3
Cyl No = 25
Exp Date = 06/05/2016
County = 08 Oper No. = 666666

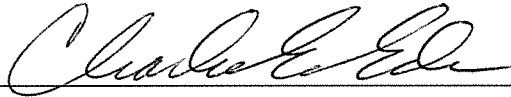

Operator Signature
CHARLES EDER

Remarks: *CALIBRATION CHECK*
0.080 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005952
Location = BCSO 8164.13.00 06/09
04/26/2016 16:05

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	16:06
02 Std. Gas	0.081	16:06
03 Room Air	0.000	16:07
04 Std. Gas	0.081	16:07
05 Room Air	0.000	16:08
06 Std. Gas	0.081	16:08
07 Room Air	0.000	16:09

Lot No = 13314080A3
Cyl No = 25
Exp Date = 06/05/2016
County = 08 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
0.080 AC

Form 106-I8000