

INTOXILYZER 8000 INSPECTION

Serial No.: 80-004189

Location: Williston

- |  | Check When Done |
|--|-----------------|
| A. Is the warm-up time less than 20 minutes?   | <u>✓</u>        |
| B. Is a three-pronged grounded outlet used?  | <u>✓</u>        |
| C. Is the breath tube heated?  | <u>✓</u>        |
| D. Is the diagnostic check complete?   | <u>✓</u>        |
| E. Is the time, date, and year correct (re-set if necessary)?  | <u>✓</u>        |
| F. Print test. (Attach test record.)   | <u>✓</u>        |
| G. Low AC. Use <0.03 AC in ACA mode. (Attach test record.)<br>Sim SN <u>DR5130</u> Lot No. <u>    </u> Rep AC <u>0.021</u>   | <u>✓</u>        |
| H. Linearity Test. Use ≥ 0.25 AC in ACA mode. (Attach test record.)<br>Sim SN <u>DR5115</u> Lot No. <u>    </u> Rep AC <u>0.298</u>                                | <u>✓</u>        |
| I. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode.<br>(Attach test record.)<br>Sim SN <u>DR5186</u> Lot No. <u>    </u> Rep AC <u>0.100</u> | <u>✓</u>        |
| J. RFI check. Run CMS mode. Key radio on first room air.<br>(Attach test record.)  | <u>✓</u>        |
| K. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.<br>(Attach test record.)  |                 |

|                       |                     |                             |
|-----------------------|---------------------|-----------------------------|
| Lot No. <u>659357</u> | Cyl No. <u>45</u>   | Exp. Date <u>09/03/2011</u> |
| Test 1 <u>0.078</u>   | Test 4 <u>0.078</u> | Test 7 <u>0.078</u>         |
| Test 2 <u>0.078</u>   | Test 5 <u>0.078</u> | Test 8 <u>0.078</u>         |
| Test 3 <u>0.078</u>   | Test 6 <u>0.077</u> | Test 9 <u>0.079</u>         |
| Average <u>0.078</u>  |                     |                             |

L. Tank pressure, Level 3 Func D Sub G P: 527 psi Reg 520 psi

M. Remarks/Maintenance record: Annual inspection

|                            |  |
|----------------------------|--|
| <u>05 May 2010</u><br>Date | <u>Deb Kashur</u><br>Field Inspector's Signature |
| <u>15 May 2010</u><br>Date | <u>Margy Pearson</u><br>Reviewed by              |

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-004189  
Location = WILL                8164.13.00 06/09  
05/05/2010                      11:08

\*\*\*\*\* Printer Test \*\*\*\*\*

abcdefghijklmnopqrstuvwxyz1234567890-=  
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#%^&\*()\_+?

abcdefghijklmnopqrstuvwxyz1234567890-=  
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#%^&\*()\_+?

Current Instrument Setup


Data Entry Mode:            Enabled  
Start Test Sequence:        DABACABA  
Display Prelim Rslt?        Yes  
Display Third Digit?        Yes  
Inhib Printer(Y/N)?        No  
Display Volume?            No  
Disable On Memfull?        Yes  
# of Print Copies?          1  
Select Std (D/W/I)?        Dry  
Standard Value?            0.080  
Standard Lot #?            659357  
Standard Cyl #?            45  
Standard Expiration?        09/03/2011  
Oper No?                    777777

Flow Cal. Date:            06/03/2009  
Slope                        666  
Intercept                  -550210

IR Calibration Date:        01/05/2010  
                                    3um                      9um


-----  
0th Coef(\*100):            -13528                      -21726  
1st Coef(\*100):            252059                      130963  
2nd Coef(\*100):            270                         387  
H2O adj(mg/l\*10k):        558                         521

\*\*\*\*\* Printer Test End \*\*\*\*\*

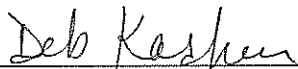

Operator Signature  
DEB KASHUR

Remarks:

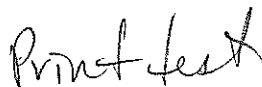

Form 106-18000

80-004189  
P. 75/82 194

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004189  
Location = WILL      8164.13.00 06/09  
05/05/2010      11:13

WET CAL CHECK

| Test         | AC    | Time  |
|--------------|-------|-------|
| 01 Room Air  | 0.000 | 11:13 |
| 02 Std. Sol. | 0.020 | 11:14 |
| 03 Room Air  | 0.000 | 11:15 |
| 04 Std. Sol. | 0.020 | 11:15 |
| 05 Room Air  | 0.000 | 11:16 |
| 06 Std. Sol. | 0.020 | 11:17 |
| 07 Room Air  | 0.000 | 11:17 |

08 Sim Temp = 34.0°C

Simul Ser No = 00000

Std Sol No = 000

County = 53

Oper No. = 777777

Deb Kasher

Operator Signature  
DEB KASHUR

Remarks:

Low AC 0.021 AC

Form 106-I8000

80-004189  
P. 76/82 YK

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004189  
Location = WILL      8164.13.00 06/09  
05/05/2010      11:18

WET CAL CHECK

| Test         | AC    | Time  |
|--------------|-------|-------|
| 01 Room Air  | 0.000 | 11:18 |
| 02 Std. Sol. | 0.288 | 11:19 |
| 03 Room Air  | 0.000 | 11:20 |
| 04 Std. Sol. | 0.285 | 11:20 |
| 05 Room Air  | 0.000 | 11:21 |
| 06 Std. Sol. | 0.283 | 11:22 |
| 07 Room Air  | 0.000 | 11:22 |

08 Sim Temp = 34.0°C

Simul Ser No = 00000  
Std Sol No = 000  
County = 53

Oper No. = 777777

Deb Kashur

Operator Signature  
DEB KASHUR

Remarks:

*Linearity test 0.298 AC*

Form 106-18000

80-004189  
P. 77/82 KA

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-004189  
Location = WILL                8164.13.00 06/09  
05/05/2010                        11:22

| Test             | AC    | Time  |
|------------------|-------|-------|
| 01 Room Air      | 0.000 | 11:23 |
| 02 *Subject Test | INT*  | 11:23 |
| 03 Room Air      | 0.000 | 11:24 |

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 02/01/1992  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = NA  
Dr. Lic. = ND/DIS921456  
Simul Ser No = 00000  
Std Sol No = 000  
County = 53                            Oper No. = 777777

Deb Kasher

Operator Signature  
DEB KASHUR

Remarks:

*Interferent check*

Form 106-I8000

*0.100 AC  
+ acetone*

*80-004189  
P. 7/8/02 1A*

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-004189  
Location = WILL                8164.13.00 06/09  
05/05/2010                        11:29

| Test          | AC    | Time  |
|---------------|-------|-------|
| 01 Diagnostic | OK    | 11:30 |
| 02 Room Air   | RFI*  | 11:30 |
| 03 Room Air   | 0.000 | 11:31 |

\*Invalid Test  
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 02/01/1992  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = NA  
Dr. Lic. = ND/DIS921456  
Lot No = 659357  
Cyl No = 45  
Expiration Date = 09/03/2011  
County = 53                            Oper No. = 777777

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

Deb Kashur

Operator Signature  
DEB KASHUR

Remarks:

RFI check

Form 106-18000

80-004189  
P. 79/82 KM

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004189  
Location = WILL      8164.13.00 06/09  
05/05/2010      11:09

DRY CAL CHECK

| Test        | AC    | Time  |
|-------------|-------|-------|
| 01 Room Air | 0.000 | 11:10 |
| 02 Std. Gas | 0.078 | 11:10 |
| 03 Room Air | 0.000 | 11:10 |
| 04 Std. Gas | 0.078 | 11:11 |
| 05 Room Air | 0.000 | 11:11 |
| 06 Std. Gas | 0.078 | 11:12 |
| 07 Room Air | 0.000 | 11:12 |

Lot No = 659357  
Cyl No = 45  
Exp Date = 09/03/2011  
County = 53      Oper No. = 777777

*Deb Kasher*

Operator Signature  
DEB KASHUR

Remarks:

*Calibration check 0.080 AC*

Form 106-I8000

80-004189  
P. 86/82 KRL

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004189  
Location = WILL      8164.13.00 06/09  
05/05/2010      11:25

DRY CAL CHECK

| Test        | AC    | Time  |
|-------------|-------|-------|
| 01 Room Air | 0.000 | 11:25 |
| 02 Std. Gas | 0.078 | 11:26 |
| 03 Room Air | 0.000 | 11:26 |
| 04 Std. Gas | 0.078 | 11:27 |
| 05 Room Air | 0.000 | 11:27 |
| 06 Std. Gas | 0.077 | 11:27 |
| 07 Room Air | 0.000 | 11:28 |

Lot No = 659357  
Cyl No = 45  
Exp Date = 09/03/2011  
County = 53      Oper No. = 777777

Deb Kashur  
Operator Signature  
DEB KASHUR

Remarks:

*calibration check 0.080Ac*

Form 106-I8000

80-004189  
P 8/82 KH



Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004189  
Location = WILL      8164.13.00 06/09  
05/05/2010      11:31

DRY CAL CHECK

| Test        | AC    | Time  |
|-------------|-------|-------|
| 01 Room Air | 0.000 | 11:32 |
| 02 Std. Gas | 0.078 | 11:32 |
| 03 Room Air | 0.000 | 11:32 |
| 04 Std. Gas | 0.078 | 11:33 |
| 05 Room Air | 0.000 | 11:33 |
| 06 Std. Gas | 0.079 | 11:34 |
| 07 Room Air | 0.000 | 11:34 |

Lot No = 659357  
Cyl No = 45  
Exp Date = 09/03/2011  
County = 53      Oper No. = 777777

*Deb Kashur*

Operator Signature  
DEB KASHUR

Remarks:

*Calibration check 0.080 AC*

Form 106-I8000

80-004189  
P. 8/8/12 *JK*