

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-005942

Location: PLPD (POWERS LAKE P.D.)

- Check When Done
- A. Is the breath tube heated? ✓
 - B. Is the diagnostic check complete? ✓
 - C. Is the time, date, and year correct (re-set if necessary)? ✓
 - D. Print test. (Attach test record.) ✓
 - E. Low AC. Use ≤ 0.03 AC in ACA mode. (Attach test record.)
Sim SN DR3803 Lot No. 15160 Rep AC 0.020 ✓
 - F. High AC. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN DR3847 Lot No. 16060 Rep AC 0.300 ✓
 - G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.
(Attach test record.)
Sim SN DR5190 Lot No. ICS4 Solution 0.10AC + 0.05% ACETONE ✓
 - H. RFI check. Run CMS mode. Key radio during CMS test.
(Attach test record.) ✓
 - I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) 0.080 AC

Lot No. 1421408046 Cyl No. 057 Exp. Date 07/05/2016

Test 1 <u>0.078</u>	Test 4 <u>0.079</u>	Test 7 <u>0.078</u>
Test 2 <u>0.078</u>	Test 5 <u>0.078</u>	Test 8 <u>0.078</u>
Test 3 <u>0.079</u>	Test 6 <u>0.078</u>	Test 9 <u>0.078</u>

Average 0.078

J. Tank pressure, Level 3 Func D Sub G P: 825 psi Reg ~825 psi

K. Remarks/Maintenance record: ANNUAL INSPECTION

5/3/2016
Date

Charles E. [Signature]
Field Inspector's Signature

05 May 16
Date

Deb Shanaver
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005942
Location = PLPD 8164.13.00 06/09
05/03/2016 10:53

***** Printer Test *****

abcdefghijklmnopqrstuvwxy1234567890-=|
ABCDEFGHIJKLMNQPQRSTUVWXYZ!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxy1234567890-=|
ABCDEFGHIJKLMNQPQRSTUVWXYZ!@#\$\$%^&*()_+?

Current Instrument Setup

Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 14214080A6
Standard Cyl #? 57
Standard Expiration? 07/05/2016
Oper No? 666666

Flow Cal. Date: 06/18/2013
Slope 715
Intercept -740740

IR Calibration Date: 05/20/2015
3um 9um

	3um	9um
0th Coef(*100):	-32049	-32893
1st Coef(*100):	272377	139703
2nd Coef(*100):	2023	1232
H2O adj(mg/l*10k):	468	428

***** Printer Test End *****



Operator Signature
CHARLES EDER

Remarks:

PRINT TEST

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005942
Location = PLPD 8164.13.00 06/09
05/03/2016 11:03

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:03
02 Std. Sol.	0.018	11:04
03 Room Air	0.000	11:05
04 Std. Sol.	0.018	11:05
05 Room Air	0.000	11:06
06 Std. Sol.	0.018	11:07
07 Room Air	0.000	11:07


08 Sim Temp = 34.0°C

Simul Ser No = DR3803

Std Sol No = 15160

County = 07

Oper No. = 666666

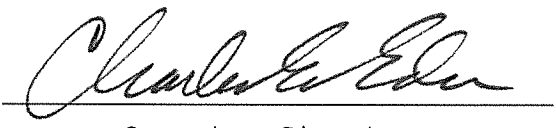

Operator Signature
CHARLES EDER

Remarks:

Low AC
0.020 AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005942
Location = PLPD 8164.13.00 06/09
05/03/2016 11:15

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:15
02 Std. Sol.	0.311	11:16
03 Room Air	0.000	11:17
04 Std. Sol.	0.310	11:17
05 Room Air	0.000	11:18
06 Std. Sol.	0.310	11:19
07 Room Air	0.000	11:19

08 Sim Temp = 34.0°C

Simul Ser No = DR3847

Std Sol No = 16060

County = 07

Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

HIGH AC
0.300 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005942
Location = PLPD 8164.13.00 06/09
05/03/2016 11:24

Test	AC	Time
01 Room Air	0.000	11:25
02 *Subject Test	INT*	11:26
03 Room Air	0.000	11:26

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = INTERFERENT
Dr. Lic. = ND/DIS821456
Lot No = 14214080A6
Cyl No = 57
Expiration Date = 07/05/2016
County = 07 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

Interferent Check
0.10AC + 0.05% ACETONE

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005942
Location = PLPD 8164.13.00 06/09
05/03/2016 11:28

Test	AC	Time
01 Diagnostic	OK	11:29
02 Room Air	0.000	11:29
03 *Subject Test	RFI*	11:29
04 Room Air	0.000	11:30

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = RFI CHECK
Dr. Lic. = ND/DIS821456
Lot No = 14214080A6
Cyl No = 57
Expiration Date = 07/05/2016
County = 07 Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature
CHARLES EDER

Remarks:

RFI CHECK

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005942
Location = PLPD 8164.13.00 06/09
05/03/2016 11:31

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:32
02 Std. Gas	0.078	11:32
03 Room Air	0.000	11:32
04 Std. Gas	0.078	11:33
05 Room Air	0.000	11:33
06 Std. Gas	0.079	11:34
07 Room Air	0.000	11:34


Lot No = 14214080A6

Cyl No = 57

Exp Date = 07/05/2016

County = 07

Oper No. = 666666

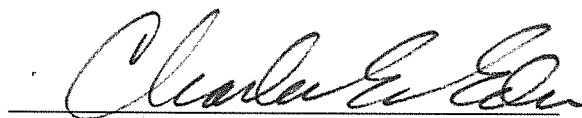

Operator Signature
CHARLES EDER

Remarks:

CALIBRATION CHECK
0.080 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005942
Location = PLPD 8164.13.00 06/09
05/03/2016 11:34

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:35
02 Std. Gas	0.079	11:35
03 Room Air	0.000	11:36
04 Std. Gas	0.078	11:36
05 Room Air	0.000	11:37
06 Std. Gas	0.078	11:37
07 Room Air	0.000	11:38


Lot No = 14214080A6

Cyl No = 57

Exp Date = 07/05/2016

County = 07

Oper No. = 666666

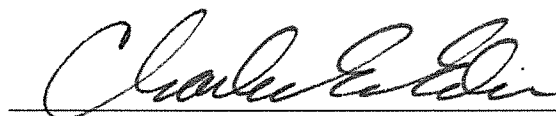

Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK

0.080AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005942
Location = PLPD 8164.13.00 06/09
05/03/2016 11:38

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:39
02 Std. Gas	0.078	11:39
03 Room Air	0.000	11:39
04 Std. Gas	0.078	11:40
05 Room Air	0.000	11:40
06 Std. Gas	0.078	11:41
07 Room Air	0.000	11:41

Lot No = 14214080A6

Cyl No = 57

Exp Date = 07/05/2016

County = 07

Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

CALIBRATION CHECK

0.080AC

Form 106-I8000