

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-003072

Location: GARR

- Check When Done
- A. Is the breath tube heated? ✓
  - B. Is the diagnostic check complete? ✓
  - C. Is the time, date, and year correct (re-set if necessary)? ✓
  - D. Print test. (Attach test record.) ✓
  - E. Low AC. Use  $\leq 0.03$  AC in ACA mode. (Attach test record.)  
Sim SN DR5143 Lot No. 201501B Rep AC 0.020 ✓
  - F. High AC. Use  $\geq 0.25$  AC in ACA mode. (Attach test record.)  
Sim SN DR3455 Lot No. 14270 Rep AC 0.300 ✓
  - G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.  
(Attach test record.)  
Sim SN DR5132 Lot No. ICS#2 Solution 0.10AC + 0.05% ACE ✓
  - H. RFI check. Run CMS mode. Key radio during CMS test.  
(Attach test record.) ✓
  - I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.  
(Attach test record.) 0.080 AC

Lot No. 33913080A4 Cyl No. 24 Exp. Date 1/1/2016

Test 1 <u>0.079</u>	Test 4 <u>0.079</u>	Test 7 <u>0.079</u>
Test 2 <u>0.079</u>	Test 5 <u>0.079</u>	Test 8 <u>0.079</u>
Test 3 <u>0.080</u>	Test 6 <u>0.080</u>	Test 9 <u>0.080</u>

Average 0.079

J. Tank pressure, Level 3 Func D Sub G P: 783 psi Reg 800 psi

K. Remarks/Maintenance record: ANNUAL INSPECTION

5.6.15  
Date

06 May 2015  
Date

Charles E. Eda  
Field Inspector's Signature

Deb Kashner  
Reviewed by

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-003072  
Location = GARR                    8164.13.00 06/09  
05/06/2015                            11:27

\*\*\*\*\* Printer Test \*\*\*\*\*

abcdefghijklmnopqrstuvwxyz1234567890-=\_|  
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&\*()\_+?

abcdefghijklmnopqrstuvwxyz1234567890-=\_|  
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&\*()\_+?

Current Instrument Setup

Data Entry Mode:            Enabled  
Start Test Sequence:        DABACABA  
Display Prelim Rslt?        Yes  
Display Third Digit?        Yes  
Inhib Printer(Y/N)?        No  
Display Volume?            No  
Disable On Memfull?        Yes  
# of Print Copies?          1  
Select Std (D/W/I)?        Dry  
Standard Value?            0.080  
Standard Lot #?            33913080A4  
Standard Cyl #?            24  
Standard Expiration?        01/01/2016  
Oper No?                    666666

Flow Cal. Date:            03/09/2011  
Slope                        682  
Intercept                   -620512

IR Calibration Date:        03/09/2011  
                                    3um                            9um

-----  
0th Coef(\*100):            -24766                    -17499  
1st Coef(\*100):            264559                    138630  
2nd Coef(\*100):            1876                      1114  
H2O adj(mg/l\*10k):        497                        331

\*\*\*\*\* Printer Test End \*\*\*\*\*



Operator Signature  
CHARLES EDER

Remarks:

PRINT TEST

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-003072  
Location = GARR      8164.13.00 06/09  
05/06/2015      11:52

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:53
02 Std. Sol.	0.017	11:53
03 Room Air	0.000	11:54
04 Std. Sol.	0.018	11:55
05 Room Air	0.000	11:55
06 Std. Sol.	0.018	11:56
07 Room Air	0.000	11:56

08 Sim Temp = 34.0°C

Simul Ser No = DR5143  
Std Sol No = 201501B  
County = 28      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: *Low AC*  
*0.020 AC*

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-003072  
Location = GARR      8164.13.00 06/09  
05/06/2015      11:59

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:00
02 Std. Sol.	0.294	12:01
03 Room Air	0.000	12:01
04 Std. Sol.	0.300	12:02
05 Room Air	0.000	12:02
06 Std. Sol.	0.300	12:03
07 Room Air	0.000	12:03

08 Sim Temp = 34.0°C

Simul Ser No = DR3455  
Std Sol No = 14270  
County = 28      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks:      HIGH AC  
                    0.300 AC


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-003072  
Location = GARR                8164.13.00 06/09  
05/06/2015                        11:43

Test	AC	Time
01 Room Air	0.000	11:44
02 *Subject Test	INT*	11:45
03 Room Air	0.000	11:45

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = INTERFERENT  
Dr. Lic. = ND/DIS821456  
Lot No = 33913080A4  
Cyl No = 24  
Expiration Date = 01/01/2016  
County = 28                            Oper No. = 666666

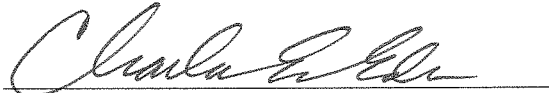

Operator Signature  
CHARLES EDER

Remarks: *Interferent Check*  
*0.10 AC + 0.05% ACETONE*

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-003072  
Location = GARR                8164.13.00 06/09  
05/06/2015                      11:30

Test	AC	Time
01 Diagnostic	OK	11:32
02 Room Air	0.000	11:32
03 *Subject Test	RFI*	11:33
04 Room Air	0.000	11:33

\*Invalid Test  
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female                      Weight = 150  
Test = DUI                              Cit = RFI CHECK  
Dr. Lic. = ND/DIS821456  
Lot No = 33913080A4  
Cyl No = 24  
Expiration Date = 01/01/2016  
County = 28                              Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature  
CHARLES EDER

Remarks:

RFI CHECK

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-003072  
Location = GARR      8164.13.00 06/09  
05/06/2015      11:34

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:35
02 Std. Gas	0.079	11:35
03 Room Air	0.000	11:36
04 Std. Gas	0.079	11:36
05 Room Air	0.000	11:37
06 Std. Gas	0.080	11:37
07 Room Air	0.000	11:38

Lot No = 33913080A4  
Cyl No = 24  
Exp Date = 01/01/2016  
County = 28      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
0.080 AC

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000        SN 80-003072  
Location = GARR                8164.13.00 06/09  
05/06/2015                        11:39

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	11:39
02 Std. Gas	0.079	11:40
03 Room Air	0.000	11:40
04 Std. Gas	0.079	11:41
05 Room Air	0.000	11:41
06 Std. Gas	0.080	11:41
07 Room Air	0.000	11:42

Lot No = 33913080A4  
Cyl No = 24  
Exp Date = 01/01/2016  
County = 28                      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks:

CALIBRATION CHECK  
0.080 AC

Form 106-I8000



Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-003072  
Location = GARR      8164.13.00 06/09  
05/06/2015      12:07

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:07
02 Std. Gas	0.079	12:08
03 Room Air	0.000	12:08
04 Std. Gas	0.079	12:08
05 Room Air	0.000	12:09
06 Std. Gas	0.080	12:09
07 Room Air	0.000	12:10

Lot No = 33913080A4  
Cyl No = 24  
Exp Date = 01/01/2016  
County = 28      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
0.080 AC