



INTOXILYZER® 8000 INSTALLATION AND REPAIR CHECKOUT
 OFFICE OF ATTORNEY GENERAL
 CRIME LABORATORY DIVISION
 Toxicology Section/Breath Alcohol Program
 SFN 59281 (1-2010)

| | |
|---|------------------------------------|
| Serial Number: <u>80-004934</u> | Instrument Location: <u>UND PD</u> |
| Reason for Install/Repair: | |
| <input checked="" type="checkbox"/> Install After Lab Repair/Inspection | |
| <input type="checkbox"/> Other (Specify) _____ | |

Check When Done:

- A. Surge Protector Installed/Properly Grounded.
- B. Telephone Line Connected to Intoxilyzer® 8000.
- C. Breath Tube Heated.
- D. Review/Enter Preliminary Data Entry (Level 2, Function E).
- E. Review/Enter Gas Setup (Level 1, Function S).
- F. Print Test (Level 1, Function P). (Attach Test Record.)
- G. ACA Test (Level 1, Function C). (Attach Test Record.)
- H. RFI Test (Use CMS Mode-First Room Air). (Attach Test Record.)
- I. Repair and/or Maintenance Performed: Install

Mary Pearson
 Field Inspector Signature

10 May 2010
 Date

Deb Kashner
 Reviewed By

15
DK 17 May 2010
 Date

80-004934
 P. 44/47 KSK

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004934
Location = UNDP 8164.13.00 06/09
05/10/2010 18:53

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-==|
ABCDEFGHIJKLMNopqrstuvwxyz!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-==|
ABCDEFGHIJKLMNopqrstuvwxyz!@#\$\$%^&*()_+?

Current Instrument Setup


Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 659358
Standard Cyl #? 31
Standard Expiration? 09/03/2011
Oper No? 107501

Flow Cal. Date: 01/20/2010
Slope 698
Intercept -694859

IR Calibration Date: 03/12/2010

| | 3um | 9um |
|--------------------|--------|--------|
| 0th Coef(*100): | -18922 | -21325 |
| 1st Coef(*100): | 252748 | 128359 |
| 2nd Coef(*100): | 509 | 713 |
| H2O adj(mg/l*10k): | 732 | 642 |

***** Printer Test End *****

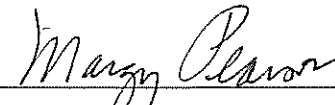

Operator Signature
MARGY PEARSON

Remarks:

Print Test

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004934
Location = UNDP 8164.13.00 06/09
05/10/2010 19:12

DRY CAL CHECK

| Test | AC | Time |
|-------------|-------|-------|
| 01 Room Air | 0.000 | 19:12 |
| 02 Std. Gas | 0.080 | 19:13 |
| 03 Room Air | 0.000 | 19:13 |
| 04 Std. Gas | 0.079 | 19:14 |
| 05 Room Air | 0.000 | 19:14 |
| 06 Std. Gas | 0.080 | 19:14 |
| 07 Room Air | 0.000 | 19:15 |

Lot No = 659358
Cyl No = 31
Exp Date = 09/03/2011
County = 18 Oper No. = 107501



Operator Signature
MARGY PEARSON

Remarks:

Install ACA

Form 106-I8000

80-004934
P. 46/47 1912

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004934
Location = UNDP 8164.13.00 06/09
05/10/2010 19:05

| Test | AC | Time |
|---------------|-------|-------|
| 01 Diagnostic | OK | 19:06 |
| 02 Room Air | RFI* | 19:06 |
| 03 Room Air | 0.000 | 19:07 |

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Lot No = 659358
Cyl No = 31
Expiration Date = 09/03/2011
County = 08 Oper No. = 107501

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

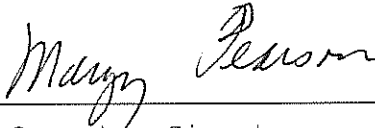

Operator Signature
MARGY PEARSON

Remarks:

RFI

Form 106-I8000

80-004934
P. 47/47 VAM