

INTOXILYZER 8000 INSPECTION

Serial No.: 80-004206

Location: UND PD

- | | Check When Done |
|---|-----------------|
| A. Is the warm-up time less than 20 minutes? | <u>✓</u> |
| B. Is a three-pronged grounded outlet used? | <u>✓</u> |
| C. Is the breath tube heated? | <u>✓</u> |
| D. Is the diagnostic check complete? | <u>✓</u> |
| E. Is the time, date, and year correct (re-set if necessary)? | <u>✓</u> |
| F. Print test. (Attach test record.) | <u>✓</u> |
| G. Low AC. Use <0.03 AC in ACA mode. (Attach test record.)
Sim SN <u>DR5132</u> Lot No. <u>—</u> Rep AC <u>0.021 AC</u> | <u>✓</u> |
| H. Linearity Test. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN <u>DR 3375</u> Lot No. <u>—</u> Rep AC <u>0.292 AC</u> | <u>✓</u> |
| I. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode.
(Attach test record.)
Sim SN <u>DR5186</u> Lot No. <u>—</u> Rep AC <u>0.10 AC + 0.05% Acetone</u> | <u>✓</u> |
| J. RFI check. Run CMS mode. Key radio on first room air.
(Attach test record.) | <u>✓</u> |
| K. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) | |

Lot No. <u>659358</u>	Cyl No. <u>31</u>	Exp. Date <u>09/03/2011</u>
Test 1 <u>0.083</u>	Test 4 <u>0.083</u>	Test 7 <u>0.082</u>
Test 2 <u>0.082</u>	Test 5 <u>0.082</u>	Test 8 <u>0.082</u>
Test 3 <u>0.083</u>	Test 6 <u>0.081</u>	Test 9 <u>0.082</u>
Average <u>0.082 AC</u>		

L. Tank pressure, Level 3 Func D Sub G P: 750 psi Reg 750 psi

M. Remarks/Maintenance record: Annual Inspection

<u>10 May 2010</u> Date	<u>Mary Pearson</u> Field Inspector's Signature
<u>15 May 2010</u> Date	<u>Deb Kashner</u> Reviewed by

80-004206
P. 75/83 1/2

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004206
Location = UNDP 8164.13.00 06/09
05/10/2010 18:17

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	18:18
02 Std. Sol.	0.020	18:19
03 Room Air	0.000	18:19
04 Std. Sol.	0.020	18:20
05 Room Air	0.000	18:21
06 Std. Sol.	0.020	18:21
07 Room Air	0.000	18:22

08 Sim Temp = 34.0°C

Simul Ser No = 0000000

Std Sol No = 0000000

County = 18

Oper No. = 107501



Operator Signature
MARGY PEARSON

Remarks:

Low AC 0.021 AC

Form 106-I8000

80-004206
P. 77/83 JEM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004206
Location = UNDP 8164.13.00 06/09
05/10/2010 18:11

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	18:12
02 Std. Sol.	0.285	18:13
03 Room Air	0.000	18:13
04 Std. Sol.	0.285	18:14
05 Room Air	0.000	18:14
06 Std. Sol.	0.286	18:15
07 Room Air	0.000	18:16


08 Sim Temp = 34.0°C

Simul Ser No = 0000000

Std Sol No = 0000000

County = 18

Oper No. = 107501

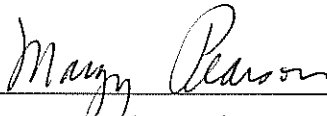

Operator Signature
MARGY PEARSON

Remarks:

Linearity 0.292 AC

Form 106-18000

80-004206
P. 7/8/83 JKH


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004206
Location = UNDP 8164.13.00 06/09
05/10/2010 17:55

Test	AC	Time
01 Room Air	0.000	17:56
02 *Subject Test	INT*	17:56
03 Room Air	0.000	17:57

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Lot No = 659358
Cyl No = 31
Expiration Date = 09/03/2011
County = 18 Oper No. = 107501

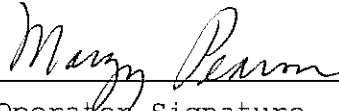

Operator Signature
MARGY PEARSON

Remarks:

Interferent - 0.10 AC + 0.05% acetone

Form 106-I8000

80-004206
P. 79/83 FK

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004206
Location = UNDP 8164.13.00 06/09
05/10/2010 17:58

Test	AC	Time
01 Diagnostic	OK	17:58
02 Room Air	RFI*	17:59
03 Room Air	0.000	17:59

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Lot No = 659358
Cyl No = 31
Expiration Date = 09/03/2011
County = 18 Oper No. = 107501

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

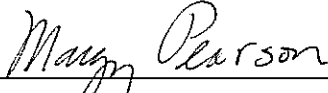

Operator Signature
MARGY PEARSON

Remarks:

RFI

Form 106-I8000

*80-004206
P. 80/83*

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004206
Location = UNDP 8164.13.00 06/09
05/10/2010 18:23

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	18:24
02 Std. Gas	0.083	18:24
03 Room Air	0.000	18:24
04 Std. Gas	0.082	18:25
05 Room Air	0.000	18:25
06 Std. Gas	0.083	18:26
07 Room Air	0.000	18:26

Lot No = 659358
Cyl No = 31
Exp Date = 09/03/2011
County = 18 Oper No. = 107501

Margy Pearson

Operator Signature
MARGY PEARSON

Remarks:

Cal Check OK ACA

Form 106-18000

80-004206
P. 8/1/83 KZ

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004206
Location = UNDP 8164.13.00 06/09
05/10/2010 18:27

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	18:27
02 Std. Gas	0.083	18:28
03 Room Air	0.000	18:28
04 Std. Gas	0.082	18:29
05 Room Air	0.000	18:29
06 Std. Gas	0.081	18:30
07 Room Air	0.000	18:30

Lot No = 659358

Cyl No = 31

Exp Date = 09/03/2011

County = 18

Oper No. = 107501

Margy Pearson

Operator Signature
MARGY PEARSON

Remarks:

Cal Check ACA

Form 106-I8000

80-004206
P. 8/2/83 Y/L


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004206
Location = UNDP 8164.13.00 06/09
05/10/2010 18:31

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	18:32
02 Std. Gas	0.082	18:32
03 Room Air	0.000	18:32
04 Std. Gas	0.082	18:33
05 Room Air	0.000	18:33
06 Std. Gas	0.082	18:34
07 Room Air	0.000	18:34

Lot No = 659358
Cyl No = 31
Exp Date = 09/03/2011
County = 18 Oper No. = 107501

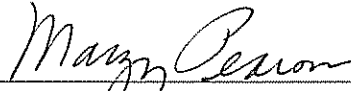

Operator Signature
MARGY PEARSON

Remarks:

Cal Check ACA

Form 106-I8000

80-004206
P. 83/83 1/2