

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-004188

Location: TOXL

- Check When Done
- A. Is the breath tube heated? ✓
 - B. Is the diagnostic check complete? ✓
 - C. Is the time, date, and year correct (re-set if necessary)? ✓
 - D. Print test. (Attach test record.) ✓
 - E. Low AC. Use ≤ 0.03 AC in ACA mode. (Attach test record.)
Sim SN MP3066 Lot No. 15160 Rep AC 0.020 ✓
 - F. High AC. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN MP3069 Lot No. 201603C Rep AC 0.400 ✓
 - G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.
(Attach test record.)
Sim SN DR5135 Lot No. ICS 4 Solution 0.10 AC + 0.05% ACETONE ✓
 - H. RFI check. Run CMS mode. Key radio during CMS test.
(Attach test record.) ✓
 - I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) 0.080 AC

Lot No. 19615080A5 Cyl No. 003 Exp. Date 09/09/2017
CEE

Test 1 <u>0.080</u>	Test 4 <u>0.081</u>	Test 7 <u>0.082</u>
Test 2 <u>0.081</u>	Test 5 <u>0.081</u>	Test 8 <u>0.081</u>
Test 3 <u>0.081</u>	Test 6 <u>0.082</u>	Test 9 <u>0.081</u>

Average 0.081

J. Tank pressure, Level 3 Func D Sub G P: 330 psi Reg ~325 psi

K. Remarks/Maintenance record: ANNUAL INSPECTION, REPLACED EXTERNAL

BATTERY
5/11/2016
Date
12 May 16 ✓
Date

Charles E. Edr
Field Inspector's Signature
Deb Shanaver
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = TOXL 8164.13.00 06/09
05/11/2016 13:58

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=-|
ABCDEFGHIJKLMNQPQRSTUVWXYZ!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=-|
ABCDEFGHIJKLMNQPQRSTUVWXYZ!@#\$\$%^&*()_+?

Current Instrument Setup

Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 19615080A5
Standard Cyl #? 3
Standard Expiration? 09/05/2017
Oper No? 666666

Flow Cal. Date: 07/05/2011
Slope 668
Intercept -597204

IR Calibration Date: 04/05/2012
 3um 9um

0th Coef(*100): -15090 -16406
1st Coef(*100): 263383 126966
2nd Coef(*100): 6533 2371
H2O adj(mg/l*10k): 801 809

***** Printer Test End *****



Operator Signature
CHARLES EDER

Remarks:

PRINT TEST

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = TOXL 8164.13.00 06/09
05/11/2016 14:07

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:07
02 Std. Sol.	0.019	14:08
03 Room Air	0.000	14:09
04 Std. Sol.	0.019	14:09
05 Room Air	0.000	14:10
06 Std. Sol.	0.019	14:11
07 Room Air	0.000	14:11

08 Sim Temp = 34.0°C

Simul Ser No = MP3066
Std Sol No = 15160
County = 08 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

Low AC
0.019 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = TOXL 8164.13.00 06/09
05/11/2016 14:30

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:31
02 Std. Sol.	0.391	14:32
03 Room Air	0.000	14:32
04 Std. Sol.	0.392	14:33
05 Room Air	0.000	14:33
06 Std. Sol.	0.391	14:34
07 Room Air	0.000	14:35

08 Sim Temp = 34.0°C

Simul Ser No = MP3069
Std Sol No = 201603C
County = 08 Oper No. = 666666

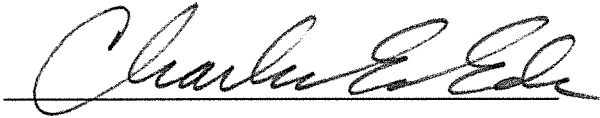

Operator Signature
CHARLES EDER

Remarks:

HIGH AC
0.400 AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = TOXL 8164.13.00 06/09
05/11/2016 14:41

Test	AC	Time
01 Room Air	0.000	14:42
02 *Subject Test	INT*	14:42
03 Room Air	0.000	14:43

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = INTERFERENT
Dr. Lic. = ND/DIS821456
Lot No = 19615080A5
Cyl No = 3
Expiration Date = 09/05/2017
County = 08 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

Interferent Check
0.10AC + 0.05% ACETONE

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = TOXL 8164.13.00 06/09
05/11/2016 14:44

Test	AC	Time
01 Diagnostic	OK	14:45
02 Room Air	0.000	14:46
03 *Subject Test	RFI*	14:46
04 Room Air	0.000	14:46

*Invalid Test
 Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = RFI CHECK
Dr. Lic. = ND/DIS821456
Lot No = 19615080A5
Cyl No = 3
Expiration Date = 09/05/2017
County = 08 Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature
CHARLES EDER

Remarks: *RFI CHECK*

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = TOXL 8164.13.00 06/09
05/11/2016 14:48

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:48
02 Std. Gas	0.080	14:49
03 Room Air	0.000	14:49
04 Std. Gas	0.081	14:50
05 Room Air	0.000	14:50
06 Std. Gas	0.081	14:50
07 Room Air	0.000	14:51

Lot No = 19615080A5
Cyl No = 3
Exp Date = 09/05/2017
County = 08 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
 0.080 AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = TOXL 8164.13.00 06/09
05/11/2016 14:51

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:52
02 Std. Gas	0.081	14:52
03 Room Air	0.000	14:53
04 Std. Gas	0.081	14:53
05 Room Air	0.000	14:54
06 Std. Gas	0.082	14:54
07 Room Air	0.000	14:55

Lot No = 19615080A5
Cyl No = 3
Exp Date = 09/05/2017
County = 08 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK

0.080 AC

Form 106-I8000


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004188
Location = TOXL 8164.13.00 06/09
05/11/2016 14:56

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:57
02 Std. Gas	0.082	14:57
03 Room Air	0.000	14:58
04 Std. Gas	0.081	14:58
05 Room Air	0.000	14:58
06 Std. Gas	0.081	14:59
07 Room Air	0.000	14:59

Lot No = 19615080A5
Cyl No = 3
Exp Date = 09/05/2017
County = 08 Oper No. = 666666

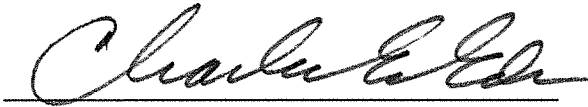

Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
0.080 AC