

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-004193

Location: Ward Co jail

- Check When Done
- A. Is the breath tube heated? ✓
 - B. Is the diagnostic check complete? ✓
 - C. Is the time, date, and year correct (re-set if necessary)? ✓
 - D. Print test. (Attach test record.) ✓
 - E. Low AC. Use ≤ 0.03 AC in ACA mode. (Attach test record.)
Sim SN DR5191 Lot No. 15160 Rep AC 0.020 ✓
 - F. High AC. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN DR 3847 Lot No. 16060 Rep AC 0.300 ✓
 - G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.
(Attach test record.)
Sim SN DR5190 Lot No. FCS4 Solution 0.10AC + 0.05% acetone ✓
 - H. RFI check. Run CMS mode. Key radio during CMS test.
(Attach test record.) ✓
 - I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) 0.080 AC

Lot No. <u>14214080A6</u>	Cyl No. <u>70</u>	Exp. Date <u>07/05/2016</u>
Test 1 <u>0.079</u>	Test 4 <u>0.079</u>	Test 7 <u>0.079</u>
Test 2 <u>0.079</u>	Test 5 <u>0.079</u>	Test 8 <u>0.078</u>
Test 3 <u>0.079</u>	Test 6 <u>0.079</u>	Test 9 <u>0.080</u>
Average <u>0.079</u>		

J. Tank pressure, Level 3 Func D Sub G P: 973 psi Reg 980 psi

K. Remarks/Maintenance record: Annual inspection, replaced o-rings on breath hose & simulator return

11 May 16
Date

5/12/2016 ✓
Date

Db. Shearer
Field Inspector's Signature

Charles E. Ed
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004193
Location = WARD 8164.13.00 06/09
05/11/2016 15:10

***** Printer Test *****

abcdefghijklmnopqrstuvwxy1234567890-=_|
ABCDEFGHIJKLMN0PQRSTUVWXYZ!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxy1234567890-=_|
ABCDEFGHIJKLMN0PQRSTUVWXYZ!@#\$\$%^&*()_+?

Current Instrument Setup

Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 14214080A6
Standard Cyl #? 70
Standard Expiration? 07/05/2016
Oper No? 777777

Flow Cal. Date: 06/04/2009
Slope 692
Intercept -611138

IR Calibration Date: 04/23/2015
 3um 9um


0th Coef(*100): -21817 -20283
1st Coef(*100): 254491 124899
2nd Coef(*100): 4701 1877
H2O adj(mg/l*10k): 502 469

***** Printer Test End *****



Operator Signature
DEB SHANAVER

Remarks:

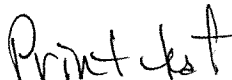

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004193
Location = WARD 8164.13.00 06/09
05/11/2016 15:23

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:23
02 Std. Sol.	0.016	15:24
03 Room Air	0.000	15:25
04 Std. Sol.	0.016	15:25
05 Room Air	0.000	15:26
06 Std. Sol.	0.016	15:27
07 Room Air	0.000	15:27

08 Sim Temp = 34.0°C

Simul Ser No = DR5191
Std Sol No = 15160
County = 51 Oper No. = 777777



Operator Signature
DEB SHANAVER

Remarks:

low alcohol 0.020 AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004193
Location = WARD 8164.13.00 06/09
05/11/2016 15:30

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:30
02 Std. Sol.	0.293	15:31
03 Room Air	0.000	15:32
04 Std. Sol.	0.293	15:32
05 Room Air	0.000	15:33
06 Std. Sol.	0.293	15:34
07 Room Air	0.000	15:34

08 Sim Temp = 34.0°C

Simul Ser No = DR3847

Std Sol No = 16060

County = 51

Oper No. = 777777

Deb Shanaver

Operator Signature
DEB SHANAVER

Remarks:

High alcohol 0.300 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004193
Location = WARD 8164.13.00 06/09
05/11/2016 15:36

Test	AC	Time
01 Room Air	0.000	15:36
02 *Subject Test	INT*	15:37
03 Room Air	0.000	15:37

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS921456
Lot No = 14214080A6
Cyl No = 70
Expiration Date = 07/05/2016
County = 51 Oper No. = 777777

Deb Shanaver

Operator Signature
DEB SHANAVER

Remarks:

Interferent detected 0.10AC + 0.05% acetone

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004193
Location = WARD 8164.13.00 06/09
05/11/2016 15:38

Test	AC	Time
01 Diagnostic	OK	15:38
02 Room Air	RFI*	15:38
03 Room Air	0.000	15:39

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS921456
Lot No = 14214080A6
Cyl No = 70
Expiration Date = 07/05/2016
County = 51 Oper No. = 777777

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature
DEB SHANAVER

Remarks:

RFI test


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004193
Location = WARD 8164.13.00 06/09
05/11/2016 15:10

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:11
02 Std. Gas	0.079	15:11
03 Room Air	0.000	15:12
04 Std. Gas	0.079	15:12
05 Room Air	0.000	15:13
06 Std. Gas	0.079	15:13
07 Room Air	0.000	15:14

Lot No = 14214080A6
Cyl No = 70
Exp Date = 07/05/2016
County = 51 Oper No. = 777777



Operator Signature
DEB SHANAVER

Remarks:

Calibration check 0.080AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004193
Location = WARD 8164.13.00 06/09
05/11/2016 15:16

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:16
02 Std. Gas	0.079	15:17
03 Room Air	0.000	15:17
04 Std. Gas	0.079	15:17
05 Room Air	0.000	15:18
06 Std. Gas	0.079	15:18
07 Room Air	0.000	15:19

Lot No = 14214080A6
Cyl No = 70
Exp Date = 07/05/2016
County = 51 Oper No. = 777777

Deb Shanaver

Operator Signature
DEB SHANAVER

Remarks:

Calibration check 0.080 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004193
Location = WARD 8164.13.00 06/09
05/11/2016 15:19

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:20
02 Std. Gas	0.079	15:20
03 Room Air	0.000	15:20
04 Std. Gas	0.078	15:21
05 Room Air	0.000	15:21
06 Std. Gas	0.080	15:22
07 Room Air	0.000	15:22

Lot No = 14214080A6
Cyl No = 70
Exp Date = 07/05/2016
County = 51 Oper No. = 777777

Deb Shanaver

Operator Signature
DEB SHANAVER

Remarks:

Calibration check 0.080 AC

Form 106-I8000