

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-004954

Location: Rolette Co SO

- Check When Done
- A. Is the breath tube heated? ✓
- B. Is the diagnostic check complete? ✓
- C. Is the time, date, and year correct (re-set if necessary)? ✓
- D. Print test. (Attach test record.) ✓
- E. Low AC. Use ≤ 0.03 AC in ACA mode. (Attach test record.)
Sim SN DR5191 Lot No. 15160 Rep AC 0.020 ✓
- F. High AC. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN DR 3847 Lot No. 16060 Rep AC 0.300 ✓
- G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.
(Attach test record.)
Sim SN DR5190 Lot No. ICS4 Solution 0.10ACT + 0.05% acetone ✓
- H. RFI check. Run CMS mode. Key radio during CMS test.
(Attach test record.) ✓
- I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) 0.080 AC

Lot No. <u>19615080AS</u>	Cyl No. <u>31</u>	Exp. Date <u>09/05/2017</u>
Test 1 <u>0.082</u>	Test 4 <u>0.083</u>	Test 7 <u>0.082</u>
Test 2 <u>0.082</u>	Test 5 <u>0.082</u>	Test 8 <u>0.082</u>
Test 3 <u>0.082</u>	Test 6 <u>0.082</u>	Test 9 <u>0.083</u>
Average <u>0.082</u>		

J. Tank pressure, Level 3 Func D Sub G P: 985 psi Reg ~1000 psi

K. Remarks/Maintenance record: DS Annual inspection, replaced internal printer paper

11 May 16
Date

12 MAY 2016 ✓
Date

Deb Shanaver
Field Inspector's Signature

Charles E. ...
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004954
Location = ROLL 8164.13.00 06/09
05/11/2016 08:54

***** Printer Test *****

abcdefghijklmnopqrstu vwxyz1234567890-=_|
ABCDEFGHIJKLMNopQRSTUvwxyz!@#\$\$%^&*()_+?

abcdefghijklmnopqrstu vwxyz1234567890-=_|
ABCDEFGHIJKLMNopQRSTUvwxyz!@#\$\$%^&*()_+?

Current Instrument Setup
Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? No
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 19615080A5
Standard Cyl #? 31
Standard Expiration? 09/05/2017
Oper No? 777777

Flow Cal. Date: 08/27/2012
Slope 670
Intercept -570610

IR Calibration Date: 08/27/2012		
	3um 9um	

0th Coef(*100):	-22689 -21316	
1st Coef(*100):	260406 132331	
2nd Coef(*100):	1084 919	
H2O adj(mg/l*10k):	570 470	

***** Printer Test End *****

Deb Shanaver

Operator Signature
DEB SHANAVER

Remarks:

Print Test

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004954
Location = ROLL 8164.13.00 06/09
05/11/2016 09:11

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:11
02 Std. Sol.	0.018	09:12
03 Room Air	0.000	09:13
04 Std. Sol.	0.017	09:13
05 Room Air	0.000	09:14
06 Std. Sol.	0.018	09:15
07 Room Air	0.000	09:15

08 Sim Temp = 34.0°C

Simul Ser No = DR5191

Std Sol No = 15160

County = 40

Oper No. = 777777

Deb Shanaver

Operator Signature
DEB SHANAVER

Remarks:

low alcohol 0.020 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004954
Location = ROLL 8164.13.00 06/09
05/11/2016 09:16

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:16
02 Std. Sol.	0.297	09:17
03 Room Air	0.000	09:18
04 Std. Sol.	0.298	09:18
05 Room Air	0.000	09:19
06 Std. Sol.	0.297	09:20
07 Room Air	0.000	09:20

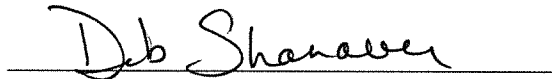
08 Sim Temp = 34.0°C

Simul Ser No = DR3847

Std Sol No = 16060

County = 40

Oper No. = 777777


Operator Signature
DEB SHANAVER

Remarks:

High alcohol 0.300 AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004954
Location = ROLL 8164.13.00 06/09
05/11/2016 09:07

Test	AC	Time
01 Room Air	0.000	09:08
02 *Subject Test	INT*	09:08
03 Room Air	0.000	09:09

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS921456
Lot No = 19615080A5
Cyl No = 31
Expiration Date = 09/05/2017
County = 40 Oper No. = 777777

Deb Shanaver

Operator Signature
DEB SHANAVER

Remarks:

Interferent detected 0.16AC + 0.05% acetone

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004954
Location = ROLL 8164.13.00 06/09
05/11/2016 09:09

Test	AC	Time
01 Diagnostic	OK	09:10
02 Room Air	RFI*	09:10
03 Room Air	0.000	09:10

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS921456
Lot No = 19615080A5
Cyl No = 31
Expiration Date = 09/05/2017
County = 40 Oper No. = 777777

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

Deb Shanaver

Operator Signature
DEB SHANAVER

Remarks:

RFI test

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004954
Location = ROLL 8164.13.00 06/09
05/11/2016 08:55

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	08:55
02 Std. Gas	0.082	08:56
03 Room Air	0.000	08:56
04 Std. Gas	0.082	08:57
05 Room Air	0.000	08:57
06 Std. Gas	0.082	08:58
07 Room Air	0.000	08:58

Lot No = 19615080A5
Cyl No = 31
Exp Date = 09/05/2017
County = 40 Oper No. = 777777

Deb Shanaver

Operator Signature
DEB SHANAVER

Remarks:

Calibration check 0.080 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004954
Location = ROLL 8164.13.00 06/09
05/11/2016 08:58

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	08:59
02 Std. Gas	0.083	08:59
03 Room Air	0.000	09:00
04 Std. Gas	0.082	09:00
05 Room Air	0.000	09:01
06 Std. Gas	0.082	09:01
07 Room Air	0.000	09:01

Lot No = 19615080A5
Cyl No = 31
Exp Date = 09/05/2017
County = 40 Oper No. = 777777

Deb Shanaver

Operator Signature
DEB SHANAVER

Remarks:

Calibration check 0.080 AC

Form 106-I8000


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004954
Location = ROLL 8164.13.00 06/09
05/11/2016 09:03

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:04
02 Std. Gas	0.082	09:04
03 Room Air	0.000	09:05
04 Std. Gas	0.082	09:05
05 Room Air	0.000	09:06
06 Std. Gas	0.083	09:06
07 Room Air	0.000	09:07

Lot No = 19615080A5
Cyl No = 31
Exp Date = 09/05/2017
County = 40 Oper No. = 777777


Operator Signature
DEB SHANAVER

Remarks:
Calibration check 0.080AC