

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-005940

Location: Griggs Co SO

- Check When Done
- A. Is the breath tube heated? ✓
- B. Is the diagnostic check complete? ✓
- C. Is the time, date, and year correct (re-set if necessary)? ✓
- D. Print test. (Attach test record.) ✓
- E. Low AC. Use ≤ 0.03 AC in ACA mode. (Attach test record.)
Sim SN DR591 Lot No. 15160 Rep AC 0.020 ✓
- F. High AC. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN DR3847 Lot No. 16060 Rep AC 0.300 ✓
- G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.
(Attach test record.)
Sim SN DR590 Lot No. ICS4 Solution 0.10 AC + 0.05% acetone ✓
- H. RFI check. Run CMS mode. Key radio during CMS test.
(Attach test record.) ✓
- I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) 0.080 AC

Lot No. 13314080 A2 Cyl No. 60 Exp. Date 06/05/16

Test 1 <u>0.079</u>	Test 4 <u>0.079</u>	Test 7 <u>0.079</u>
Test 2 <u>0.079</u>	Test 5 <u>0.079</u>	Test 8 <u>0.079</u>
Test 3 <u>0.079</u>	Test 6 <u>0.079</u>	Test 9 <u>0.078</u>

Average 0.079

J. Tank pressure, Level 3 Func D Sub G P: 371 psi Reg ~380 psi

K. Remarks/Maintenance record: Annual ^{DS} inspection

09 May 16
Date

12 MAY 2016 ✓
Date

Deb Shanaver
Field Inspector's Signature

Charles E. Edin
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005940
Location = COOP 8164.13.00 06/09
05/09/2016 14:29

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=_|
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&*()_+?

Current Instrument Setup


Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 13314080A2
Standard Cyl #? 60
Standard Expiration? 06/05/2016
Oper No? 777777

Flow Cal. Date: 06/21/2013
Slope 688
Intercept -659053

IR Calibration Date: 04/02/2015
 3um 9um


0th Coef(*100): -18917 -25088
1st Coef(*100): 267616 136158
2nd Coef(*100): 3374 1493
H2O adj(mg/l*10k): 427 417

***** Printer Test End *****



Operator Signature
DEB SHANAVER

Remarks:

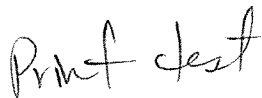

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005940
Location = COOP 8164.13.00 06/09
05/09/2016 14:31

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:31
02 Std. Sol.	0.019	14:32
03 Room Air	0.000	14:33
04 Std. Sol.	0.019	14:34
05 Room Air	0.000	14:34
06 Std. Sol.	0.019	14:35
07 Room Air	0.000	14:35

08 Sim Temp = 34.0°C

Simul Ser No = DR5191
Std Sol No = 15160
County = 20 Oper No. = 777777



Operator Signature
DEB SHANAVER

Remarks:

low AC 0.020AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005940
Location = COOP 8164.13.00 06/09
05/09/2016 14:37

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:37
02 Std. Sol.	0.314	14:38
03 Room Air	0.000	14:39
04 Std. Sol.	0.313	14:40
05 Room Air	0.000	14:40
06 Std. Sol.	0.313	14:41
07 Room Air	0.000	14:41

08 Sim Temp = 34.0°C

Simul Ser No = DR3847

Std Sol No = 16060

County = 20

Oper No. = 777777

Deb Shanaver

Operator Signature
DEB SHANAVER

Remarks:

High AC
Calibration check 0.350 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005940
Location = COOP 8164.13.00 06/09
05/09/2016 14:44

Test	AC	Time
01 Room Air	0.000	14:44
02 *Subject Test	INT*	14:45
03 Room Air	0.000	14:46

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS921456
Lot No = 13314080A2
Cyl No = 60
Expiration Date = 06/05/2016
County = 20 Oper No. = 777777

Deb Shanaver

Operator Signature
DEB SHANAVER

Remarks:

*Interferent detected 0.10 AC+
0.05% acetone*

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005940
Location = COOP 8164.13.00 06/09
05/09/2016 14:47

Test	AC	Time
01 Diagnostic	OK	14:48
02 Room Air	RFI*	14:48
03 Room Air	0.000	14:48

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS921456
Lot No = 13314080A2
Cyl No = 60
Expiration Date = 06/05/2016
County = 20 Oper No. = 777777

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

Deb Shanaver

Operator Signature
DEB SHANAVER

Remarks:

RFI test

Form 106-I8000


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005940
Location = COOP 8164.13.00 06/09
05/09/2016 14:11

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:12
02 Std. Gas	0.079	14:12
03 Room Air	0.000	14:12
04 Std. Gas	0.079	14:13
05 Room Air	0.000	14:13
06 Std. Gas	0.079	14:14
07 Room Air	0.000	14:14

Lot No = 13314080A2
Cyl No = 60
Exp Date = 06/05/2016
County = 20 Oper No. = 777777



Operator Signature
DEB SHANAVER

Remarks:
Calibration check 0.080 AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005940
Location = COOP 8164.13.00 06/09
05/09/2016 14:15

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:15
02 Std. Gas	0.079	14:16
03 Room Air	0.000	14:16
04 Std. Gas	0.079	14:16
05 Room Air	0.000	14:17
06 Std. Gas	0.079	14:17
07 Room Air	0.000	14:18

Lot No = 13314080A2
Cyl No = 60
Exp Date = 06/05/2016
County = 20 Oper No. = 777777

Deb Shanaver
Operator Signature
DEB SHANAVER

Remarks:

Calibration check 0.080 AC


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005940
Location = COOP 8164.13.00 06/09
05/09/2016 14:19

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:19
02 Std. Gas	0.079	14:20
03 Room Air	0.000	14:20
04 Std. Gas	0.079	14:20
05 Room Air	0.000	14:21
06 Std. Gas	0.078	14:21
07 Room Air	0.000	14:22

Lot No = 13314080A2
Cyl No = 60
Exp Date = 06/05/2016
County = 20 Oper No. = 777777



Operator Signature
DEB SHANAVER

Remarks:
Calibration check 0.080 AC