

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-005948

Location: Walsh Co So

- Check When Done
- A. Is the breath tube heated? ✓
 - B. Is the diagnostic check complete? ✓
 - C. Is the time, date, and year correct (re-set if necessary)? ✓
 - D. Print test. (Attach test record.) ✓
 - E. Low AC. Use ≤ 0.03 AC in ACA mode. (Attach test record.)
Sim SN DR5191 Lot No. 15160 Rep AC 0.020 ✓
 - F. High AC. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN DR3847 Lot No. 16060 Rep AC 0.300 ✓
 - G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.
(Attach test record.)
Sim SN DR5190 Lot No. IC54 Solution 0.10AC + 0.05% acetone ✓
 - H. RFI check. Run CMS mode. Key radio during CMS test.
(Attach test record.) ✓
 - I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) 0.080 AC

Lot No. 13314080A2 Cyl No. 55 Exp. Date 06/05/2016

Test 1 <u>0.079</u>	Test 4 <u>0.080</u>	Test 7 <u>0.080</u>
Test 2 <u>0.079</u>	Test 5 <u>0.079</u>	Test 8 <u>0.080</u>
Test 3 <u>0.079</u>	Test 6 <u>0.080</u>	Test 9 <u>0.079</u>

Average 0.079

J. Tank pressure, Level 3 Func D Sub G P: 101 psi Reg 100 psi

K. Remarks/Maintenance record: Annual inspection

10 May 16
Date

12 MAY 2016 ✓
Date

Deb Shanever
Field Inspector's Signature
Charles E. Ed
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005948
Location = GRAF 8164.13.00 06/09
05/10/2016 12:24

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-=
ABCDEFGHIJKLMNPNOPQRSTUVWXYZ!@#%&^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-=
ABCDEFGHIJKLMNPNOPQRSTUVWXYZ!@#%&^&*()_+?

Current Instrument Setup

Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 13314080A2
Standard Cyl #? 55
Standard Expiration? 06/05/2016
Oper No? 777777

Flow Cal. Date: 06/20/2013
Slope 701
Intercept -716431

IR Calibration Date: 04/02/2015
 3um 9um


0th Coef(*100): -25939 -20035
1st Coef(*100): 272687 134661
2nd Coef(*100): 3113 1576
H2O adj(mg/l*10k): 531 397

***** Printer Test End *****



Operator Signature
DEB SHANAVER

Remarks:

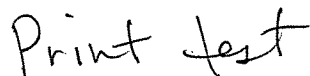

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005948
Location = GRAF 8164.13.00 06/09
05/10/2016 12:27

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:28
02 Std. Sol.	0.018	12:29
03 Room Air	0.000	12:29
04 Std. Sol.	0.018	12:30
05 Room Air	0.000	12:30
06 Std. Sol.	0.018	12:31
07 Room Air	0.000	12:32

08 Sim Temp = 34.0°C

Simul Ser No = DR5191
Std Sol No = 15160
County = 50 Oper No. = 777777

Deb Shanaver
Operator Signature
DEB SHANAVER

Remarks:

low alcohol 0.020 AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005948
Location = GRAF 8164.13.00 06/09
05/10/2016 12:35

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:36
02 Std. Sol.	0.314	12:36
03 Room Air	0.000	12:37
04 Std. Sol.	0.313	12:38
05 Room Air	0.000	12:38
06 Std. Sol.	0.312	12:39
07 Room Air	0.000	12:40

08 Sim Temp = 34.0°C

Simul Ser No = DR3847

Std Sol No = 16060

County = 50

Oper No. = 777777

Deb Shanaver

Operator Signature
DEB SHANAVER

Remarks:

^{DS}
~~Ju~~ High alcohol 0.300AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005948
Location = GRAF 8164.13.00 06/09
05/10/2016 12:41

Test	AC	Time
01 Room Air	0.000	12:41
02 *Subject Test	INT*	12:42
03 Room Air	RFI**	12:42

*Invalid Test
Interferent Detected

**Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS921456
Lot No = 13314080A2
Cyl No = 55
Expiration Date = 06/05/2016
County = 50 Oper No. = 777777

Deb Shanaver
Operator Signature
DEB SHANAVER

Remarks:

Interferent detected 0.10AC + 0.0590acetone
(plus RFI)

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005948
Location = GRAF 8164.13.00 06/09
05/10/2016 12:43

Test	AC	Time
01 Diagnostic	OK	12:44
02 Room Air	RFI*	12:44
03 Room Air	0.000	12:45

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS921456
Lot No = 13314080A2
Cyl No = 55
Expiration Date = 06/05/2016
County = 50 Oper No. = 777777

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

Deb Shanaver

Operator Signature
DEB SHANAVER

Remarks:

RFI test

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005948
Location = GRAF 8164.13.00 06/09
05/10/2016 12:11

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:11
02 Std. Gas	0.079	12:12
03 Room Air	0.000	12:12
04 Std. Gas	0.079	12:13
05 Room Air	0.000	12:13
06 Std. Gas	0.079	12:13
07 Room Air	0.000	12:14

Lot No = 13314080A2
Cyl No = 55
Exp Date = 06/05/2016
County = 50 Oper No. = 777777

Deb Shanaver
Operator Signature
DEB SHANAVER

Remarks:

Calibration check 0.080AC

Form 106-I8000


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005948
Location = GRAF 8164.13.00 06/09
05/10/2016 12:14

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:15
02 Std. Gas	0.080	12:15
03 Room Air	0.000	12:16
04 Std. Gas	0.079	12:16
05 Room Air	0.000	12:17
06 Std. Gas	0.080	12:17
07 Room Air	0.000	12:18

Lot No = 13314080A2
Cyl No = 55
Exp Date = 06/05/2016
County = 50 Oper No. = 777777



Operator Signature
DEB SHANAVER

Remarks:
Calibration check 0.080AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-005948
Location = GRAF 8164.13.00 06/09
05/10/2016 12:18

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:19
02 Std. Gas	0.080	12:19
03 Room Air	0.000	12:20
04 Std. Gas	0.080	12:20
05 Room Air	0.000	12:21
06 Std. Gas	0.079	12:21
07 Room Air	0.000	12:22

Lot No = 13314080A2
Cyl No = 55
Exp Date = 06/05/2016
County = 50 Oper No. = 777777

Deb Shanaver
Operator Signature
DEB SHANAVER

Remarks:
Calibration check 0.080 AC