

INTOXILYZER 8000 INSPECTION

Serial No.: 80-004935

Location: Harvey PD

- | | Check When Done |
|---|-----------------|
| A. Is the warm-up time less than 20 minutes? | <u>✓</u> |
| B. Is a three-pronged grounded outlet used? | <u>✓</u> |
| C. Is the breath tube heated? | <u>✓</u> |
| D. Is the diagnostic check complete? | <u>✓</u> |
| E. Is the time, date, and year correct (re-set if necessary)? | <u>✓</u> |
| F. Print test. (Attach test record.) | <u>✓</u> |
| G. Low AC. Use <0.03 AC in ACA mode. (Attach test record.)
Sim SN <u>DR5132</u> Lot No. <u>—</u> Rep AC <u>0.021 AC</u> | <u>✓</u> |
| H. Linearity Test. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN <u>DR3375</u> Lot No. <u>—</u> Rep AC <u>0.292 AC</u> | <u>✓</u> |
| I. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode.
(Attach test record.)
Sim SN <u>DR5186</u> Lot No. <u>—</u> Rep AC <u>0.10 AC + Acetone (0.05%)</u> | <u>✓</u> |
| J. RFI check. Run CMS mode. Key radio on first room air.
(Attach test record.) | <u>✓</u> |
| K. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) | |

Lot No. 709352 Cyl No. 45 Exp. Date 12/23/11

Test 1 0.077 Test 4 0.076 Test 7 0.077

Test 2 0.076 Test 5 0.076 Test 8 MAP
0.076 0.077

Test 3 0.076 Test 6 0.076 Test 9 0.076

Average 0.076 AC

L. Tank pressure, Level 3 Func D Sub G P: 980 psi Reg 1000 psi

M. Remarks/Maintenance record: Annual Inspection; Download, Clear Forms

13 May 2010
Date

Margy Pearson
Field Inspector's Signature

15 May 2010
Date

Deb Kashner
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004935
Location = HARV 8164.13.00 06/09
05/13/2010 12:49

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:50
02 Std. Sol.	0.020	12:52
03 Room Air	0.000	12:52
04 Std. Sol.	0.020	12:53
05 Room Air	0.000	12:53
06 Std. Sol.	0.020	12:54
07 Room Air	0.000	12:54


08 Sim Temp = 34.0°C

Simul Ser No = 0000000

Std Sol No = 00000

County = 52

Oper No. = 107501

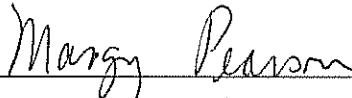

Operator Signature
MARGY PEARSON

Remarks:

Low AC 0.021 AC

Form 106-I8000

80-004935
P. 44/52 PA

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004935
Location = HARV 8164.13.00 06/09
05/13/2010 12:55

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	12:56
02 Std. Sol.	0.285	12:57
03 Room Air	0.000	12:57
04 Std. Sol.	0.287	12:58
05 Room Air	0.000	12:58
06 Std. Sol.	0.288	12:59
07 Room Air	0.000	13:00


08 Sim Temp = 34.0°C

Simul Ser No = 0000000

Std Sol No = 00000

County = 52

Oper No. = 107501

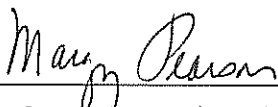

Operator Signature
MARGY PEARSON

Remarks:

Linearity 0.292 AC

Form 106-I8000

80-004935
P. 47/52 KR


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004935
Location = HARV 8164.13.00 06/09
05/13/2010 13:01

Test	AC	Time
01 Room Air	0.000	13:01
02 *Subject Test	INT*	13:02
03 Room Air	0.000	13:02

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = 0000000
Std Sol No = 00000
County = 52 Oper No. = 107501

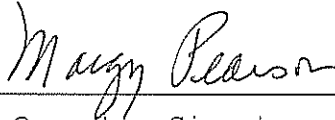

Operator Signature
MARGY PEARSON

Remarks:

Interferent - AC + Acetone

Form 106-I8000

80-004935
P. 48/52 *MA*

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004935
Location = HARV 8164.13.00 06/09
05/13/2010 13:05

Test	AC	Time
01 Diagnostic	OK	13:06
02 Room Air	RFI*	13:06
03 Room Air	0.000	13:07

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS821456
Simul Ser No = 0000000
Std Sol No = 00000
County = 52 Oper No. = 107501

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

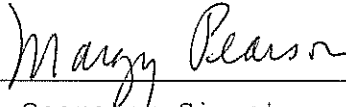

Operator Signature
MARGY PEARSON

Remarks:

RFI

Form 106-I8000

80-004935
P.49/52 FM

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004935
Location = HARV 8164.13.00 06/09
05/13/2010 13:09

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:10
02 Std. Gas	0.077	13:11
03 Room Air	0.000	13:11
04 Std. Gas	0.076	13:12
05 Room Air	0.000	13:12
06 Std. Gas	0.076	13:13
07 Room Air	0.000	13:13

Lot No = 709352
Cyl No = 45
Exp Date = 12/23/2011
County = 52 Oper No. = 107501

Margy Pearson

Operator Signature
MARGY PEARSON

Remarks:

Cal Check 0.080 AC

Form 106-I8000

80-004935
P. 50/52 FR


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004935
Location = HARV 8164.13.00 06/09
05/13/2010 13:14

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:15
02 Std. Gas	0.076	13:15
03 Room Air	0.000	13:16
04 Std. Gas	0.076	13:16
05 Room Air	0.000	13:16
06 Std. Gas	0.076	13:17
07 Room Air	0.000	13:17

Lot No = 709352
Cyl No = 45
Exp Date = 12/23/2011
County = 52 Oper No. = 107501

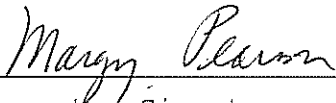

Operator Signature
MARGY PEARSON

Remarks:

Cal Check 0.080 AC

Form 106-I8000


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004935
Location = HARV 8164.13.00 06/09
05/13/2010 13:18

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:18
02 Std. Gas	0.077	13:18
03 Room Air	0.000	13:19
04 Std. Gas	0.077	13:19
05 Room Air	0.000	13:20
06 Std. Gas	0.076	13:20
07 Room Air	0.000	13:21

Lot No = 709352
Cyl No = 45
Exp Date = 12/23/2011
County = 52 Oper No. = 107501

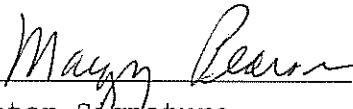

Operator Signature
MARGY PEARSON

Remarks:

Cal Check 0.080 AC

Form 106-I8000

80-004935
P. 52/52 YH