



ETHANOL GAS STANDARD CYLINDER REPORT
 OFFICE OF ATTORNEY GENERAL
 CRIME LABORATORY DIVISION
 Toxicology Section/Breath Alcohol Program
 SFN 59282 (10-2013)

Chemical Test Operator Name (Print) CHARLES E. EDER		
Location LINT	Intoxilyzer® Serial Number 80-004953	
Gas Lot Number 19615080AS	Gas Cylinder Number 005	Gas Expiration Date 9/5/2017

Check When Done:

- 1. Scan/Enter Gas Cylinder Information (Level 1, Function S).
- 2. Perform an ACA Test (Level 1, Function C).
 Write in the ACA Test Results (Reported to 3 Digits, Ex. 0.081).
 - A. 0. 079 AC
 - B. 0. 080 AC
 - C. 0. 079 AC
- 3. Complete the Top Portion of the Intoxilyzer® Record (SFN50496, Form 120-G) and Place it by the Intoxilyzer® for Use.
- 4. File Previous Intoxilyzer® Record (SFN50496, Form 120-G) at the Intoxilyzer® Location at the Agency.
- 5. Send the Following to the Crime Laboratory:
 - A. Ethanol Gas Standard Cylinder Report (SFN59282, Form 105-G).
 - B. ACA Test Record.

Chemical Test Operator Signature Charles E. Eder	Date 5/13/16 ✓
Reviewed By (Crime Laboratory Use Only) Deb Shanaver	Date 16 May 16


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004953
Location = LINT 8164.13.00 06/09
05/13/2016 09:44

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:46
02 Std. Gas	0.079	09:46
03 Room Air	0.000	09:47
04 Std. Gas	0.080	09:47
05 Room Air	0.000	09:48
06 Std. Gas	0.079	09:48
07 Room Air	0.000	09:49

Lot No = 19615080A5
Cyl No = 5
Exp Date = 09/05/2017
County = 15 Oper No. = 666666

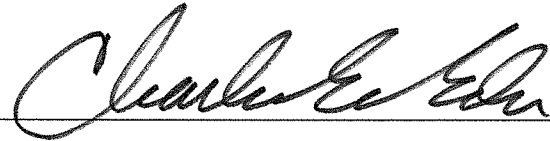

Operator Signature
CHARLES EDER

Remarks:

ACA
0.080AC