

INTOXILYZER 8000 INSPECTION

Serial No.: 80-003072

Location: Kidder Co SO

- Check When Done
- A. Is the warm-up time less than 20 minutes? ✓
 - B. Is a three-pronged grounded outlet used? ✓
 - C. Is the breath tube heated? ✓
 - D. Is the diagnostic check complete? ✓
 - E. Is the time, date, and year correct (re-set if necessary)? ✓
 - F. Print test. (Attach test record.) ✓
 - G. Low AC. Use <0.03 AC in ACA mode. (Attach test record.)
Sim SN DR5132 Lot No. 12010 Rep AC 0.020 ✓
 - H. Linearity Test. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN DR3803 Lot No. 12050 Rep AC 0.300 ✓
 - I. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode.
(Attach test record.)
Sim SN DR3846 Lot No. — Rep AC 0.100 ✓
 - J. RFI check. Run CMS mode. Key radio on first room air.
(Attach test record.) ✓
 - K. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.)

Lot No. 11611080A2 Cyl No. 90 Exp. Date 06/01/2013

Test 1 <u>0.079</u>	Test 4 <u>0.080</u>	Test 7 <u>0.080</u>
Test 2 <u>0.079</u>	Test 5 <u>0.080</u>	Test 8 <u>0.080</u>
Test 3 <u>0.080</u>	Test 6 <u>0.080</u>	Test 9 <u>0.080</u>

Average 0.080

L. Tank pressure, Level 3 Func D Sub G P: 872 psi Reg 800 psi

M. Remarks/Maintenance record: Annual inspection

15 May 2012
Date

Deb Kashur
Field Inspector's Signature

17 May 12
Date

Charles E. Eden
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-003072
Location = STEE 8164.13.00 06/09
05/15/2012 09:08

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-|=|
ABCDEFGHIJKLMNQPQRSTUVWXYZ!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-|=|
ABCDEFGHIJKLMNQPQRSTUVWXYZ!@#\$\$%^&*()_+?

Current Instrument Setup

Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 11611080A2
Standard Cyl #? 90
Standard Expiration? 06/01/2013
Oper No? 777777

Flow Cal. Date: 03/09/2011
Slope 682
Intercept -620512

IR Calibration Date: 03/09/2011
 3um 9um

	3um	9um
0th Coef(*100):	-24766	-17499
1st Coef(*100):	264559	138630
2nd Coef(*100):	1876	1114
H2O adj(mg/l*10k):	497	331

***** Printer Test End *****



Operator Signature
DEB KASHUR

Remarks:



Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-003072
Location = STEE 8164.13.00 06/09
05/15/2012 09:24

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:25
02 Std. Sol.	0.016	09:25
03 Room Air	0.000	09:26
04 Std. Sol.	0.016	09:27
05 Room Air	0.000	09:27
06 Std. Sol.	0.016	09:28
07 Room Air	0.000	09:28


08 Sim Temp = 34.0°C

Simul Ser No = 0

Std Sol No = 0

County = 22

Oper No. = 777777

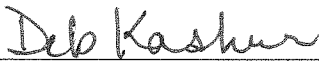

Operator Signature
DEB KASHUR

Remarks:

Low AC 0.020 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-003072
Location = STEE 8164.13.00 06/09
05/15/2012 09:36

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:36
02 Std. Sol.	0.285	09:38
03 Room Air	0.000	09:38
04 Std. Sol.	0.288	09:39
05 Room Air	0.000	09:39
06 Std. Sol.	0.291	09:40
07 Room Air	0.000	09:41

08 Sim Temp = 34.0°C

Simul Ser No = 0

Std Sol No = 0

County = 22

Oper No. = 777777

Deb Kashur

Operator Signature
DEB KASHUR

Remarks:

Linearity test 0.300 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-003072
Location = STEE 8164.13.00 06/09
05/15/2012 09:47

Test	AC	Time
01 Room Air	0.000	09:48
02 *Subject Test	INT*	09:48
03 Room Air	0.000	09:49

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS921456
Lot No = 11611080A2
Cyl No = 90
Expiration Date = 06/01/2013
County = 22 Oper No. = 777777

Deb Kasher

Operator Signature
DEB KASHUR

Remarks:

Interferent check 0.100 AC + acetone

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-003072
Location = STEE 8164.13.00 06/09
05/15/2012 09:44

Test	AC	Time
01 Diagnostic	OK	09:45
02 Room Air	RFI*	09:45
03 Room Air	0.000	09:46

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 02/01/1992
Sub Sex = Female Weight = 150
Test = DUI Cit = NA
Dr. Lic. = ND/DIS921456
Lot No = 11611080A2
Cyl No = 90
Expiration Date = 06/01/2013
County = 22 Oper No. = 777777

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

Deb Kasher

Operator Signature
DEB KASHUR

Remarks:

RFI test

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-003072
Location = STEE 8164.13.00 06/09
05/15/2012 09:09

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:10
02 Std. Gas	0.079	09:10
03 Room Air	0.000	09:11
04 Std. Gas	0.079	09:11
05 Room Air	0.000	09:12
06 Std. Gas	0.080	09:12
07 Room Air	0.000	09:13

Lot No = 11611080A2
Cyl No = 90
Exp Date = 06/01/2013
County = 22 Oper No. = 777777

Deb Kasher

Operator Signature
DEB KASHUR

Remarks:

Calibration check 0.080 AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-003072
Location = STEE 8164.13.00 06/09
05/15/2012 09:13

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:14
02 Std. Gas	0.080	09:14
03 Room Air	0.000	09:14
04 Std. Gas	0.080	09:15
05 Room Air	0.000	09:15
06 Std. Gas	0.080	09:16
07 Room Air	0.000	09:16

Lot No = 11611080A2
Cyl No = 90
Exp Date = 06/01/2013
County = 22 Oper No. = 777777

Deb Kasher

Operator Signature
DEB KASHUR

Remarks:

Calibration check 0.080 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-003072
Location = STEE 8164.13.00 06/09
05/15/2012 09:16

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:17
02 Std. Gas	0.080	09:17
03 Room Air	0.000	09:18
04 Std. Gas	0.080	09:18
05 Room Air	0.000	09:19
06 Std. Gas	0.080	09:19
07 Room Air	0.000	09:19

Lot No = 11611080A2
Cyl No = 90
Exp Date = 06/01/2013
County = 22 Oper No. = 777777

Deb Kasher

Operator Signature
DEB KASHUR

Remarks:

Calibration check 0.080 AC

Form 106-I8000