

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-004937

Location: Cass Co Jail

- | | Check When Done |
|--|-----------------|
| A. Is the breath tube heated? | <u>✓</u> |
| B. Is the diagnostic check complete? | <u>✓</u> |
| C. Is the time, date, and year correct (re-set if necessary)? | <u>✓</u> |
| D. Print test. (Attach test record.) | <u>✓</u> |
| E. Low AC. Use ≤ 0.03 AC in ACA mode. (Attach test record.)
Sim SN <u>DR5191</u> Lot No. <u>15160</u> Rep AC <u>0.020</u> | <u>✓</u> |
| F. High AC. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN <u>DR3847</u> Lot No. <u>16060</u> Rep AC <u>0.300</u> | <u>✓</u> |
| G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.
(Attach test record.)
Sim SN <u>DR5190</u> Lot No. <u>ICS4</u> Solution <u>0.10AC + 0.05% acetone</u> | <u>✓</u> |
| H. RFI check. Run CMS mode. Key radio during CMS test.
(Attach test record.) | <u>✓</u> |
| I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) <u>0.080</u> AC | |

Lot No. <u>29614080A2</u>	Cyl No. <u>36</u>	Exp. Date <u>12/05/2016</u>
Test 1 <u>0.080</u>	Test 4 <u>0.080</u>	Test 7 <u>0.080</u>
Test 2 <u>0.080</u>	Test 5 <u>0.081</u>	Test 8 <u>0.081</u>
Test 3 <u>0.080</u>	Test 6 <u>0.079</u>	Test 9 <u>0.081</u>
Average <u>0.080</u>		

J. Tank pressure, Level 3 Func D Sub G P: 735 psi Reg 750 psi

K. Remarks/Maintenance record: Annual inspection replaced breath hose simulator return o-rings

18 May 16 ✓
Date
5/19/2016 ✓
Date

Deb Shanaver
Field Inspector's Signature
Charles E. [Signature]
Reviewed by

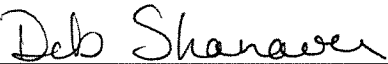
Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004937
Location = CASS 8164.13.00 06/09
05/18/2016 09:07

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:07
02 Std. Sol.	0.017	09:08
03 Room Air	0.000	09:08
04 Std. Sol.	0.017	09:09
05 Room Air	0.000	09:10
06 Std. Sol.	0.018	09:10
07 Room Air	0.000	09:11
08 Sim Temp = 34.0°C		

Simul Ser No = 15160
Std Sol No = DR5191
County = 09 Oper No. = 777777


Operator Signature
DEB KASHUR

Remarks:

low alcohol 0.020 AC

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004937
Location = CASS 8164.13.00 06/09
05/18/2016 09:14

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	09:15
02 Std. Sol.	0.290	09:16
03 Room Air	0.000	09:16
04 Std. Sol.	0.292	09:17
05 Room Air	0.000	09:18
06 Std. Sol.	0.292	09:18
07 Room Air	0.000	09:19

08 Sim Temp = 34.0°C

Simul Ser No = 16060

Std Sol No = DR3847

County = 09

Oper No. = 777777

Deb Sheuaver

Operator Signature
DEB KASHUR

Remarks:

High alcohol 0.300 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004937
Location = CASS 8164.13.00 06/09
05/18/2016 08:54

Test	AC	Time
01 Room Air	0.000	08:56
02 *Subject Test	INT*	08:56
03 Room Air	0.000	08:57

*Invalid Test
Interferent Detected

Sub Name = INTERFERENT, NA NA
Sub DOB = 08/23/1978
Sub Sex = Unknown Weight = NA
Test = DUI Cit = NA
Dr. Lic. = ND/NA
Lot No = 29614080A2
Cyl No = 36
Expiration Date = 12/05/2016
County = 09 Oper No. = 777777

Deb Shanaver

Operator Signature
DEB KASHUR

Remarks:

Interferent detected 0.10AC + 0.05% acetone

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004937
Location = CASS 8164.13.00 06/09
05/18/2016 08:59

Test	AC	Time
01 Diagnostic	OK	08:59
02 Room Air	RFI*	09:00
03 Room Air	0.000	09:00

*Invalid Test
 Inhibited - RFI

Sub Name = RFI, NA NA
Sub DOB = 08/23/1978
Sub Sex = Unknown Weight = NA
Test = DUI Cit = NA
Dr. Lic. = ND/NA
Lot No = 29614080A2
Cyl No = 36
Expiration Date = 12/05/2016
County = 09 Oper No. = 777777

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

Deb Shanaver

Operator Signature
DEB KASHUR

Remarks:

RFI test

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004937
Location = CASS 8164.13.00 06/09
05/18/2016 08:42

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	08:43
02 Std. Gas	0.080	08:43
03 Room Air	0.000	08:44
04 Std. Gas	0.080	08:44
05 Room Air	0.000	08:44
06 Std. Gas	0.080	08:45
07 Room Air	0.000	08:45

Lot No = 29614080A2

Cyl No = 36

Exp Date = 12/05/2016

County = 09

Oper No. = 777777

Deb Shanaver

Operator Signature
DEB KASHUR

Remarks:

Calibration check 0.080AC

Form 106-I8000


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004937
Location = CASS 8164.13.00 06/09
05/18/2016 08:46

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	08:46
02 Std. Gas	0.080	08:47
03 Room Air	0.000	08:47
04 Std. Gas	0.081	08:47
05 Room Air	0.000	08:48
06 Std. Gas	0.079	08:48
07 Room Air	0.000	08:49

Lot No = 29614080A2
Cyl No = 36
Exp Date = 12/05/2016
County = 09 Oper No. = 777777



Operator Signature
DEB KASHUR

Remarks:

Calibration check 0.080AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004937
Location = CASS 8164.13.00 06/09
05/18/2016 08:51

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	08:51
02 Std. Gas	0.080	08:52
03 Room Air	0.000	08:52
04 Std. Gas	0.081	08:53
05 Room Air	0.000	08:53
06 Std. Gas	0.081	08:53
07 Room Air	0.000	08:54

Lot No = 29614080A2
Cyl No = 36
Exp Date = 12/05/2016
County = 09 Oper No. = 777777

Debs Shanaver
Operator Signature
DEB KASHUR

Remarks:

Calibration check 0.080AC

Form 106-I8000