

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-004189

Location: Bott. Co SO

- Check When Done
- A. Is the breath tube heated? ✓
  - B. Is the diagnostic check complete? ✓
  - C. Is the time, date, and year correct (re-set if necessary)? ✓
  - D. Print test. (Attach test record.) ✓
  - E. Low AC. Use  $\leq 0.03$  AC in ACA mode. (Attach test record.)  
Sim SN DR5143 Lot No. 261501B Rep AC 0.020 ✓
  - F. High AC. Use  $\geq 0.25$  AC in ACA mode. (Attach test record.)  
Sim SN DR5132 Lot No. 14270 Rep AC 0.300 ✓
  - G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.  
(Attach test record.)  
Sim SN DR5132 Lot No. 1CS 2 Solution 0.10 AC + 0.05% acetone ✓
  - H. RFI check. Run CMS mode. Key radio during CMS test.  
(Attach test record.) ✓
  - I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.  
(Attach test record.) 0.080 AC

Lot No. 13314080A2 Cyl No. 52 Exp. Date 06/05/2016  
Test 1 0.082 Test 4 0.082 Test 7 0.082  
Test 2 0.082 Test 5 0.082 Test 8 0.082  
Test 3 0.082 Test 6 0.082 ~~0.083~~ DK Test 9 0.083 ~~0.082~~ DK  
Average 0.082

J. Tank pressure, Level 3 Func D Sub G P: 1607 psi Reg 1000 psi

K. Remarks/Maintenance record: Annual inspection

19 May 15  
Date

20 MAY 15 ✓  
Date

Deb Kasher  
Field Inspector's Signature

Charles E. Eden  
Reviewed by

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-004189  
Location = BOTT                    8164.13.00 06/09  
05/19/2015                            13:26

\*\*\*\*\* Printer Test \*\*\*\*\*

abcdefghijklmnopqrstuvwxyz1234567890-=\_|  
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&\*()\_+?

abcdefghijklmnopqrstuvwxyz1234567890-=\_|  
ABCDEFGHIJKLMNPOQRSTUVWXYZ!@#\$\$%^&\*()\_+?

Current Instrument Setup

Data Entry Mode:            Enabled  
Start Test Sequence:        DABACABA  
Display Prelim Rslt?        Yes  
Display Third Digit?        Yes  
Inhib Printer(Y/N)?        No  
Display Volume?            No  
Disable On Memfull?        Yes  
# of Print Copies?          1  
Select Std (D/W/I)?        Dry  
Standard Value?            0.080  
Standard Lot #?            13314080A2  
Standard Cyl #?            52  
Standard Expiration?        06/05/2016  
Oper No?                    777777

Flow Cal. Date:            02/02/2011  
Slope                      628  
Intercept                  -452608

IR Calibration Date:        02/02/2011  
                                  3um                    9um

-----  
0th Coef(\*100):            -22075                    -26801  
1st Coef(\*100):            257433                    132480  
2nd Coef(\*100):            2132                      1061  
H2O adj(mg/l\*10k):        437                        430

\*\*\*\*\* Printer Test End \*\*\*\*\*



Operator Signature  
DEB KASHUR

Remarks:



Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-004189  
Location = BOTT            8164.13.00 06/09  
05/19/2015                            13:38

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:39
02 Std. Sol.	0.018	13:40
03 Room Air	0.000	13:40
04 Std. Sol.	0.019	13:41
05 Room Air	0.000	13:41
06 Std. Sol.	0.019	13:42
07 Room Air	0.000	13:43

08 Sim Temp = 34.0°C

Simul Ser No = DR5143

Std Sol No = 201501B

County = 05

Oper No. = 777777

Deb Kasher

Operator Signature  
DEB KASHUR

Remarks:

Low AC 0.020 AC

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004189  
Location = BOTT      8164.13.00 06/09  
05/19/2015      13:43

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:44
02 Std. Sol.	0.298	13:44
03 Room Air	0.000	13:45
04 Std. Sol.	0.297	13:46
05 Room Air	0.000	13:46
06 Std. Sol.	0.297	13:47
07 Room Air	0.000	13:48

08 Sim Temp = 34.0°C

Simul Ser No = DR3455

Std Sol No = 14270

County = 05

Oper No. = 777777

Deb Kashur

Operator Signature

DEB KASHUR

Remarks:

High AC 0.300 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-004189  
Location = BOTT                    8164.13.00 06/09  
05/19/2015                            13:36

Test	AC	Time
01 Room Air	0.000	13:36
02 *Subject Test	INT*	13:37
03 Room Air	0.000	13:37

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 02/01/1992  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = NA  
Dr. Lic. = ND/DIS921456  
Lot No = 13314080A2  
Cyl No = 52  
Expiration Date = 06/05/2016  
County = 05                            Oper No. = 777777



Operator Signature  
DEB KASHUR

Remarks:

*Interferent check 0.10Ac + 0.05% acetone*

Form 106-I8000

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

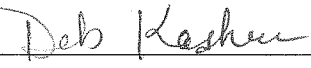
CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-004189  
Location = BOTT                    8164.13.00 06/09  
05/19/2015                            13:34

Test	AC	Time
01 Diagnostic	OK	13:35
02 Room Air	RFI*	13:35
03 Room Air	0.000	13:35

\*Invalid Test  
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 02/01/1992  
Sub Sex = Female                    Weight = 150  
Test = DUI                            Cit = NA  
Dr. Lic. = ND/DIS921456  
Lot No = 13314080A2  
Cyl No = 52  
Expiration Date = 06/05/2016  
County = 05                            Oper No. = 777777

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.

  
\_\_\_\_\_  
Operator Signature  
DEB KASHUR

Remarks:

*RFI check*

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer            Alcohol Analyzer  
North Dakota Model 8000            SN 80-004189  
Location = BOTT                    8164.13.00 06/09  
05/19/2015                            13:22

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:23
02 Std. Gas	0.082	13:23
03 Room Air	0.000	13:24
04 Std. Gas	0.082	13:24
05 Room Air	0.000	13:25
06 Std. Gas	0.082	13:25
07 Room Air	0.000	13:26

Lot No = 13314080A2  
Cyl No = 52  
Exp Date = 06/05/2016  
County = 05                            Oper No. = 777777

Deb Kasher  
Operator Signature  
DEB KASHUR

Remarks:

Calibration check 0.080 AC


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004189  
Location = BOTT      8164.13.00 06/09  
05/19/2015      13:15

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:16
02 Std. Gas	0.082	13:16
03 Room Air	0.000	13:17
04 Std. Gas	0.082	13:17
05 Room Air	0.000	13:17
06 Std. Gas	0.082	13:18
07 Room Air	0.000	13:18

Lot No = 13314080A2  
Cyl No = 52  
Exp Date = 06/05/2016  
County = 05      Oper No. = 777777

  
\_\_\_\_\_  
Operator Signature  
DEB KASHUR

Remarks:

Calibration check 0.080AC



Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-004189  
Location = BOTT      8164.13.00 06/09  
05/19/2015      13:19

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:19
02 Std. Gas	0.082	13:20
03 Room Air	0.000	13:20
04 Std. Gas	0.082	13:21
05 Room Air	0.000	13:21
06 Std. Gas	0.083	13:21
07 Room Air	0.000	13:22

Lot No = 13314080A2  
Cyl No = 52  
Exp Date = 06/05/2016  
County = 05      Oper No. = 777777

Deb Kashur  
Operator Signature  
DEB KASHUR

Remarks:

Calibration check 0.080AC