

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-005956

Location: TOXL

- Check When Done
- A. Is the breath tube heated? ✓
  - B. Is the diagnostic check complete? ✓
  - C. Is the time, date, and year correct (re-set if necessary)? ✓
  - D. Print test. (Attach test record.) ✓
  - E. Low AC. Use  $\leq 0.03$  AC in ACA mode. (Attach test record.)  
Sim SN MP3066 Lot No. 201507F Rep AC 0.020 ✓
  - F. High AC. Use  $\geq 0.25$  AC in ACA mode. (Attach test record.)  
Sim SN MP3069 Lot No. 201603C Rep AC 0.400 ✓
  - G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.  
(Attach test record.)  
Sim SN DR5135 Lot No. IC84 Solution 0.10AC+0.05% ACETONE ✓
  - H. RFI check. Run CMS mode. Key radio during CMS test.  
(Attach test record.) ✓
  - I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.  
(Attach test record.) 0.080 AC

Lot No. 13314080A3 Cyl No. 034 Exp. Date 06/05/2014

Test 1 <u>0.079</u>	Test 4 <u>0.079</u>	Test 7 <u>0.079</u>
Test 2 <u>0.079</u>	Test 5 <u>0.079</u>	Test 8 <u>0.079</u>
Test 3 <u>0.079</u>	Test 6 <u>0.078</u>	Test 9 <u>0.079</u>

Average 0.079

J. Tank pressure, Level 3 Func D Sub G P: 821 psi Reg 850 psi

K. Remarks/Maintenance record: ANNUAL INSPECTION

23 MAY 2016  
Date

Charles E. Ed  
Field Inspector's Signature

23 May 16  
Date

Deb Shanaver  
Reviewed by



Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-005956  
Location = TOXL      8164.13.00 06/09  
05/23/2016      13:31

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:32
02 Std. Sol.	0.018	13:32
03 Room Air	0.000	13:33
04 Std. Sol.	0.018	13:34
05 Room Air	0.000	13:34
06 Std. Sol.	0.018	13:35
07 Room Air	0.000	13:35

08 Sim Temp = 34.0°C

Simul Ser No = MP3066  
Std Sol No = 201507F  
County = 08      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks:      Low AC  
                 0.020 AC

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-005956  
Location = TOXL      8164.13.00 06/09  
05/23/2016      13:57

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:57
02 Std. Sol.	0.389	13:58
03 Room Air	0.000	13:59
04 Std. Sol.	0.392	14:00
05 Room Air	0.000	14:00
06 Std. Sol.	0.391	14:01
07 Room Air	0.000	14:01

08 Sim Temp = 34.0°C

Simul Ser No = MP3069  
Std Sol No = 201603C  
County = 08      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks:      HIGH AC  
                    0.400 AC


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-005956  
Location = TOXL      8164.13.00 06/09  
05/23/2016      14:06

Test	AC	Time
01 Room Air	0.000	14:07
02 *Subject Test	INT*	14:07
03 Room Air	0.000	14:08

\*Invalid Test  
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female      Weight = 150  
Test = DUI      Cit = INTERFERENT  
Dr. Lic. = ND/DIS821456  
Lot No = 13314080A3  
Cyl No = 34  
Expiration Date = 06/05/2016  
County = 08      Oper No. = 666666

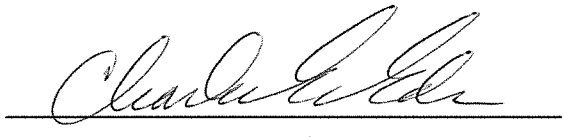

Operator Signature  
CHARLES EDER

Remarks: *Interferent Check*  
*0.10AC + 0.05% ACETONE*

Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-005956  
Location = TOXL      8164.13.00 06/09  
05/23/2016      14:09

Test	AC	Time
01 Diagnostic	OK	14:09
02 Room Air	0.000	14:10
03 *Subject Test	RFI*	14:10
04 Room Air	0.000	14:10

\*Invalid Test  
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT  
Sub DOB = 01/01/1982  
Sub Sex = Female      Weight = 150  
Test = DUI      Cit = RFI CHECK  
Dr. Lic. = ND/DIS821456  
Lot No = 13314080A3  
Cyl No = 34  
Expiration Date = 06/05/2016  
County = 08      Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature  
CHARLES EDER

Remarks:



Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-005956  
Location = TOXL      8164.13.00 06/09  
05/23/2016      14:15

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:16
02 Std. Gas	0.079	14:16
03 Room Air	0.000	14:16
04 Std. Gas	0.079	14:17
05 Room Air	0.000	14:17
06 Std. Gas	0.079	14:18
07 Room Air	0.000	14:18

Lot No = 13314080A3  
Cyl No = 34  
Exp Date = 06/05/2016  
County = 08      Oper No. = 666666



Operator Signature  
CHARLES EDER

Remarks:

*CALIBRATION CHECK*  
*0.080 AC*


Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-005956  
Location = TOXL      8164.13.00 06/09  
05/23/2016      14:20

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:20
02 Std. Gas	0.079	14:20
03 Room Air	0.000	14:21
04 Std. Gas	0.079	14:21
05 Room Air	0.000	14:22
06 Std. Gas	0.078	14:22
07 Room Air	0.000	14:23

Lot No = 13314080A3  
Cyl No = 34  
Exp Date = 06/05/2016  
County = 08      Oper No. = 666666

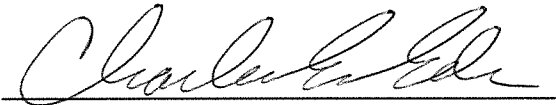

Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
0.080 AC




Intoxilyzer Test Record and Checklist  
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer      Alcohol Analyzer  
North Dakota Model 8000      SN 80-005956  
Location = TOXL      8164.13.00 06/09  
05/23/2016      14:23

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:24
02 Std. Gas	0.079	14:24
03 Room Air	0.000	14:25
04 Std. Gas	0.079	14:25
05 Room Air	0.000	14:25
06 Std. Gas	0.079	14:26
07 Room Air	0.000	14:26

Lot No = 13314080A3  
Cyl No = 34  
Exp Date = 06/05/2016  
County = 08      Oper No. = 666666

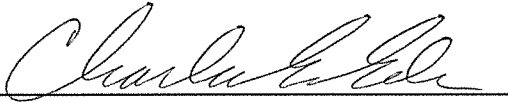  
Operator Signature  
CHARLES EDER

Remarks: CALIBRATION CHECK  
0.080 AC