

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-004186

Location: OLIV (OLIVER County S.C.)

- Check When Done
- A. Is the breath tube heated? ✓
 - B. Is the diagnostic check complete? ✓
 - C. Is the time, date, and year correct (re-set if necessary)? ✓
 - D. Print test. (Attach test record.) ✓
 - E. Low AC. Use ≤ 0.03 AC in ACA mode. (Attach test record.)
Sim SN DR5191 Lot No. 15160 Rep AC 0.200 ✓
 - F. High AC. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN DR3847 Lot No. 16060 Rep AC 0.300 ✓
 - G. Interferent check. Use 0.05% acetone in ethanol solution in ABA mode.
(Attach test record.)
Sim SN DR5190 Lot No. ICS4 Solution 0.10AC + 0.05% ACETONE ✓
 - H. RFI check. Run CMS mode. Key radio during CMS test.
(Attach test record.) ✓
 - I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.) 0.080 AC

Lot No. <u>19615080A5</u>	Cyl No. <u>033</u>	Exp. Date <u>09/05/2017</u>
Test 1 <u>0.078</u>	Test 4 <u>0.078</u>	Test 7 <u>0.079</u>
Test 2 <u>0.078</u>	Test 5 <u>0.078</u>	Test 8 <u>0.078</u>
Test 3 <u>0.078</u>	Test 6 <u>0.078</u>	Test 9 <u>0.078</u>
Average <u>0.078</u>		

J. Tank pressure, Level 3 Func D Sub G P: 662 psi Reg -650 psi

K. Remarks/Maintenance record: ANNUAL INSPECTION

5/24/2016
Date

Charles E. ...
Field Inspector's Signature

25 May 16 ✓
Date

Deb Shanaver
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004186
Location = OLIV 8164.13.00 06/09
05/24/2016 13:18

***** Printer Test *****

abcdefghijklmnopqrstuvwxyz1234567890-|=|
ABCDEFGHIJKLMNPNOPQRSTUVWXYZ!@#\$\$%^&*()_+?

abcdefghijklmnopqrstuvwxyz1234567890-|=|
ABCDEFGHIJKLMNPNOPQRSTUVWXYZ!@#\$\$%^&*()_+?

Current Instrument Setup

Data Entry Mode: Enabled
Start Test Sequence: DABACABA
Display Prelim Rslt? Yes
Display Third Digit? Yes
Inhib Printer(Y/N)? No
Display Volume? No
Disable On Memfull? Yes
of Print Copies? 1
Select Std (D/W/I)? Dry
Standard Value? 0.080
Standard Lot #? 19615080A5
Standard Cyl #? 33
Standard Expiration? 09/05/2017
Oper No? 666666

Flow Cal. Date: 04/22/2014
Slope 589
Intercept -71396

IR Calibration Date: 12/10/2015
 3um 9um

0th Coef(*100): -37549 -32658
1st Coef(*100): 272347 137615
2nd Coef(*100): 2049 1261
H2O adj(mg/l*10k): 399 507

***** Printer Test End *****



Operator Signature
CHARLES EDER

Remarks: PRINT TEST

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004186
Location = OLIV 8164.13.00 06/09
05/24/2016 13:40

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:41
02 Std. Sol.	0.017	13:42
03 Room Air	0.000	13:42
04 Std. Sol.	0.016	13:43
05 Room Air	0.000	13:43
06 Std. Sol.	0.016	13:44
07 Room Air	0.000	13:45

08 Sim Temp = 34.0°C

Simul Ser No = DR5191

Std Sol No = 15160

County = 33

Oper No. = 666666



Operator Signature

CHARLES EDER

Remarks:

Low AC

0.020 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004186
Location = OLIV 8164.13.00 06/09
05/24/2016 13:50

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:51
02 Std. Sol.	0.303	13:51
03 Room Air	0.000	13:52
04 Std. Sol.	0.305	13:53
05 Room Air	0.000	13:53
06 Std. Sol.	0.305	13:54
07 Room Air	0.000	13:54


08 Sim Temp = 34.0°C

Simul Ser No = DR3847

Std Sol No = 16060

County = 33

Oper No. = 666666

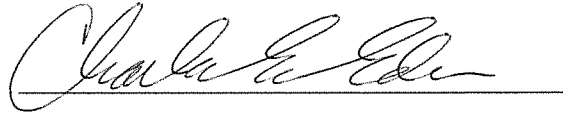

Operator Signature

CHARLES EDER

Remarks:

HIGH AC
0.300 AC

Form 106-18000


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004186
Location = OLIV 8164.13.00 06/09
05/24/2016 13:27

Test	AC	Time
01 Room Air	0.000	13:28
02 *Subject Test	INT*	13:28
03 Room Air	0.000	13:29

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = INTERFERENT
Dr. Lic. = ND/DIS821456
Lot No = 19615080A5
Cyl No = 33
Expiration Date = 09/05/2017
County = 33 Oper No. = 666666

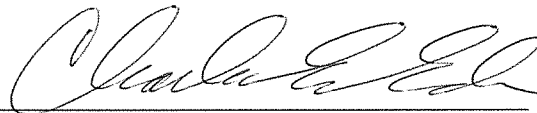

Operator Signature
CHARLES EDER

Remarks:

Interferent Check
0.10 Ac + 0.05% ACETONE

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004186
Location = OLIV 8164.13.00 06/09
05/24/2016 13:24

Test	AC	Time
01 Diagnostic	OK	13:25
02 Room Air	0.000	13:25
03 *Subject Test	RFI*	13:25
04 Room Air	0.000	13:26

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = RFI CHECK
Dr. Lic. = ND/DIS821456
Lot No = 19615080A5
Cyl No = 33
Expiration Date = 09/05/2017
County = 33 Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature
CHARLES EDER

Remarks:




Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004186
Location = OLIV 8164.13.00 06/09
05/24/2016 13:56

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	13:57
02 Std. Gas	0.078	13:57
03 Room Air	0.000	13:58
04 Std. Gas	0.078	13:58
05 Room Air	0.000	13:58
06 Std. Gas	0.078	13:59
07 Room Air	0.000	13:59

Lot No = 19615080A5
Cyl No = 33
Exp Date = 09/05/2017
County = 33 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
0.080 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004186
Location = OLIV 8164.13.00 06/09
05/24/2016 14:00

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:00
02 Std. Gas	0.078	14:01
03 Room Air	0.000	14:01
04 Std. Gas	0.078	14:02
05 Room Air	0.000	14:02
06 Std. Gas	0.078	14:02
07 Room Air	0.000	14:03

Lot No = 19615080A5
Cyl No = 33
Exp Date = 09/05/2017
County = 33 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
0.080 AC

Form 106-I8000


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004186
Location = OLIV 8164.13.00 06/09
05/24/2016 14:03

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	14:04
02 Std. Gas	0.079	14:04
03 Room Air	0.000	14:05
04 Std. Gas	0.078	14:05
05 Room Air	0.000	14:05
06 Std. Gas	0.078	14:06
07 Room Air	0.000	14:06

Lot No = 19615080A5
Cyl No = 33
Exp Date = 09/05/2017
County = 33 Oper No. = 666666

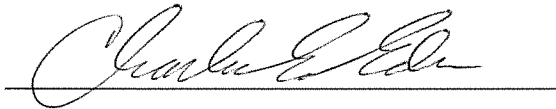

Operator Signature
CHARLES EDER

Remarks: CALIBRATION CHECK
0.080 AC