

INTOXILYZER® 8000 INSPECTION

Serial No.: 80-004950

Location: WILL

- Check When Done
- A. Is the breath tube heated? ✓
 - B. Is the diagnostic check complete? ✓
 - C. Is the time, date, and year correct (re-set if necessary)? ✓
 - D. Print test. (Attach test record.) ✓
 - E. Low AC. Use <0.03 AC in ACA mode. (Attach test record.)
Sim SN DR5115 Lot No. 12111 Rep AC 0.020 ✓
 - F. High AC. Use ≥ 0.25 AC in ACA mode. (Attach test record.)
Sim SN DR5143 Lot No. 13080 Rep AC 0.300 ✓
 - G. Interferent check. Use 0.05% acetone plus 0.10 AC ethanol in ABA mode.
(Attach test record.)
Sim SN DR3375 Lot No. IC5#2 Rep AC 0.100AC + 0.05% ACETONE ✓
 - H. RFI check. Run CMS mode. Key radio on first room air.
(Attach test record.) ✓
 - I. Calibration Check. Use Ethanol Breath Standard cylinder. Do three sets in ACA mode.
(Attach test record.)

Lot No. 16512080A2 Cyl No. 061 Exp. Date 08/01/14
Test 1 0.080 Test 4 0.080 Test 7 0.081
Test 2 0.079 Test 5 0.080 Test 8 0.081
Test 3 0.080 Test 6 0.080 Test 9 0.080
Average 0.080

J. Tank pressure, Level 3 Func D Sub G P: 75 psi Reg 75 psi

K. Remarks/Maintenance record: REPLACED BREATH HOSE AND SIMULATOR
RETURN O-RING, ANNUAL INSPECTION

28 MAY 14
Date
29 May 14 ✓
Date

Chanda E. Eden
Field Inspector's Signature
Deb Kashner
Reviewed by

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004950
Location = WILL 8164.13.00 06/09
05/28/2014 15:45

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:46
02 Std. Sol.	0.019	15:46
03 Room Air	0.000	15:47
04 Std. Sol.	0.018	15:48
05 Room Air	0.000	15:48
06 Std. Sol.	0.020	15:49
07 Room Air	0.000	15:49

08 Sim Temp = 34.0°C

Simul Ser No = DR5115

Std Sol No = 12111

County = 53

Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

Low AC
0.020 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004950
Location = WILL 8164.13.00 06/09
05/28/2014 15:50

WET CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:51
02 Std. Sol.	0.293	15:52
03 Room Air	0.000	15:52
04 Std. Sol.	0.298	15:53
05 Room Air	0.000	15:53
06 Std. Sol.	0.300	15:54
07 Room Air	0.000	15:55

08 Sim Temp = 34.0°C

Simul Ser No = DR5143
Std Sol No = 13080
County = 53 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

HIGH AC
0.300 AC

Form 106-I8000


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004950
Location = WILL 8164.13.00 06/09
05/28/2014 15:38

Test	AC	Time
01 Room Air	0.000	15:38
02 *Subject Test	INT*	15:39
03 Room Air	0.000	15:40

*Invalid Test
Interferent Detected

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = INTERFERENT
Dr. Lic. = ND/DIS821456
Lot No = 16512080A2
Cyl No = 61
Expiration Date = 08/01/2014
County = 53 Oper No. = 666666

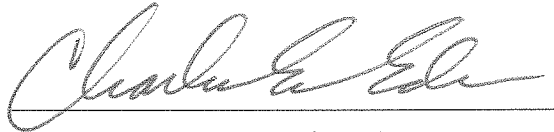

Operator Signature
CHARLES EDER

Remarks:

*INTERFERENT
CHECK
0.100 AC + 0.05% ACETONE*

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501


CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004950
Location = WILL 8164.13.00 06/09
05/28/2014 15:42

Test	AC	Time
01 Diagnostic	OK	15:43
02 Room Air	RFI*	15:43
03 Room Air	0.000	15:44

*Invalid Test
Inhibited - RFI

Sub Name = DISCOVER, THE SPIRIT
Sub DOB = 01/01/1982
Sub Sex = Female Weight = 150
Test = DUI Cit = RFI CHECK
Dr. Lic. = ND/DIS821456
Lot No = 16512080A2
Cyl No = 61
Expiration Date = 08/01/2014
County = 53 Oper No. = 666666

I followed the Approved Method and the instructions displayed by the Intoxilyzer in conducting this test.



Operator Signature
CHARLES EDER

Remarks:

RFI CHECK

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004950
Location = WILL 8164.13.00 06/09
05/28/2014 15:56

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	15:57
02 Std. Gas	0.080	15:57
03 Room Air	0.000	15:58
04 Std. Gas	0.079	15:58
05 Room Air	0.000	15:59
06 Std. Gas	0.080	15:59
07 Room Air	0.000	16:00

Lot No = 16512080A2
Cyl No = 61
Exp Date = 08/01/2014
County = 53 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks: CALIBRATION
CHECK
0.080 AC

Form 106-I8000

Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004950
Location = WILL 8164.13.00 06/09
05/28/2014 16:00

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	16:00
02 Std. Gas	0.080	16:01
03 Room Air	0.000	16:01
04 Std. Gas	0.080	16:02
05 Room Air	0.000	16:02
06 Std. Gas	0.080	16:02
07 Room Air	0.000	16:03

Lot No = 16512080A2
Cyl No = 61
Exp Date = 08/01/2014
County = 53 Oper No. = 666666



Operator Signature
CHARLES EDER

Remarks:

CALIBRATION
CHECK
0.080 AC

Form 106-I8000


Intoxilyzer Test Record and Checklist
NDOAG Crime Lab. Div., Bismarck, ND 58501

CMI, Inc. Intoxilyzer Alcohol Analyzer
North Dakota Model 8000 SN 80-004950
Location = WILL 8164.13.00 06/09
05/28/2014 16:03

DRY CAL CHECK

Test	AC	Time
01 Room Air	0.000	16:04
02 Std. Gas	0.081	16:04
03 Room Air	0.000	16:05
04 Std. Gas	0.081	16:05
05 Room Air	0.000	16:05
06 Std. Gas	0.080	16:06
07 Room Air	0.000	16:06

Lot No = 16512080A2
Cyl No = 61
Exp Date = 08/01/2014
County = 53 Oper No. = 666666

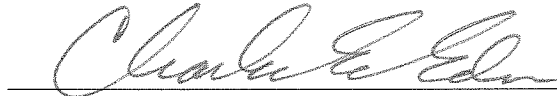

Operator Signature
CHARLES EDER

Remarks: CALIBRATION
CHECK
0.080 AC

Form 106-I8000